

Click-iT® EdU Flow Cytometry Assay Kits

100 tests

Catalog numbers C10419, C10420

Table 1 Contents and storage

Material	Amount	Concentration	Storage*	Stability	
EdU (Component A)	20 mg	NA		When stored as directed, this kit is stable for up to 1 year after receipt.	
Alexa Fluor [®] 488 azide (Cat. no. C10420) or Alexa Fluor [®] 647 azide (Cat. no. C10419) (Component B)	1 vial	NA	- - - • 2°C-8°C		
Dimethylsulfoxide (DMSO) (Component C)	8.5 mL	NA			
Click-iT [®] fixative (Component D)	10 mL	4% paraformaldehyde in PBS	DesiccateProtect from light		
Click-iT [®] saponin-based permeabilization and wash reagent (Component E)	100 mL	10X solution	Do not freeze		
CuSO ₄ (Component F)	1 mL	100 mM aqueous solution			
Click-iT [®] EdU buffer additive (Component G)	400 mg	NA			

^{*}These storage conditions are appropriate when storing the entire kit upon receipt. For optimal storage conditions for each component, see vial labels. NA = Not applicable.

Number of assays: Sufficient material is supplied for 100 reactions, based on the protocol below.


Approximate fluorescence excitation/emission maxima: Alexa Fluor® 647 azide: 650/670 nm; Alexa Fluor® 488 azide: 495/519 nm.

MAN0003748 | MP10419 Revision A.0

Measuring a cell's ability to proliferate is a fundamental method for assessing cell health, determining genotoxicity, and evaluating anti-cancer drugs. The most accurate method is by directly measuring DNA synthesis. Initially, this was performed by incorporation of radioactive nucleosides, i.e., ³H-thymidine. This method was replaced by antibody-based detection of the nucleoside analog bromo-deoxyuridine (BrdU). The Click-iT[®] EdU Flow Cytometry Assay Kits are novel alternatives to the BrdU assay. EdU (5-ethynyl-2'-deoxyuridine) is a nucleoside analog to thymidine and is incorporated into DNA during active DNA synthesis. Detection is based on a click reaction, ¹⁻⁴ a copper catalyzed covalent reaction between an azide and an alkyne. In this application, the alkyne is found in the ethynyl moiety of EdU, while the azide is coupled to Alexa Fluor® 647 dye or Alexa Fluor® 488 dye. Standard flow cytometry methods are used for determining the percentage of S-phase cells in the population (Figure 1).

The advantage of Click-iT[®] EdU labeling is that the small size of the dye azide allows for efficient detection of the incorporated EdU using mild conditions. Standard aldehyde-based fixation and detergent permeabilization are sufficient for the Click-iT® detection reagent to gain access to the DNA. This is in contrast to BrdU assays that require DNA denaturation (using acid, heat, or digestion with DNase) to expose the BrdU so that it may be detected with an anti-BrdU antibody. Sample processing for the BrdU assay can result in signal alteration of the cell cycle distribution as well as the destruction of antigen recognition sites when using the acid denaturation method. In contrast, the EdU cell proliferation kit is compatible with cell cycle dyes (Figure 2, page 3). The EdU assay can also be multiplexed with antibodies against surface and intracellular markers. However, some reagents or antibody conjugates may not be compatible with the Click-iT® EdU detection reaction and may need some additional steps to ensure compatibility (see Table 2, page 3 for details).

Figure 1 Fluorescence signal from Alexa Fluor® 488 and Alexa Fluor® 647 Click-iT® EdU Flow Cytometry Assay Kits. Jurkat (human T-cell leukemia) cells were treated with 10 µM EdU for 2 hours and detected according to the recommended staining protocol. The figures show a clear separation of proliferating cells which have incorporated EdU and nonproliferating cells which have not. Panel A shows data from cells labeled with Alexa Fluor® 488 azide analyzed on an Attune® Acoustic Focusing Cytometer using 488 nm excitation and a 530/30 nm bandpass emission filter; Panel B shows data from cells labeled with Alexa Fluor® 647 azide analyzed on a flow cytometer using 633 nm excitation and a 660/20 nm bandpass emission filter.

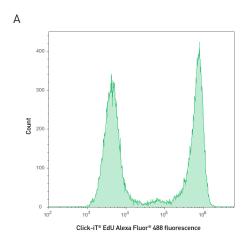

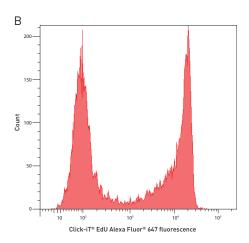


Figure 2 Dual parameter plot of Click-iT® EdU Alexa Fluor® 488 and FxCycle™ Violet. Jurkat (human T-cell leukemia) cells were treated with 10 μ M EdU for 2 hours and detected according to the recommended staining protocol. Data were collected and analyzed using an Attune® Acoustic Cytometer using 488 nm excitation and a 530/30 bandpass for detection of the EdU Alexa Fluor® 488 azide and 405 nm excitation and a 450/40 bandpass for detection of the FxCycle™ Violet fluorescence. This figure combines DNA content with EdU; cells that are positive for both labels are in S-phase of the cell cycle.

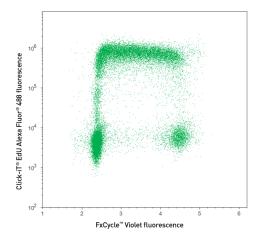


Table 2 Click-iT® EdU detection reagent compatibility

Fluorescent molecule	Compatibility*
Qdot [®] nanocrystals	Use Qdot [®] nanocrystals after the Click-iT [®] detection reaction.
Fluorescent proteins (GFP)	Use anti-GFP antibodies before the Click-iT® detection reaction or use organic dye-based reagents for protein expression detection.
Organic dyes such as Alexa Fluor® dyes, fluorescein (FITC)	Compatible
PerCP, allophycocyanin (APC) and APC-based tandems (i.e., Alexa Fluor® 680-APC)	Compatible
R-phycoerythrin (R-PE) and R-PE based tandems (i.e., Alexa Fluor® 610-RPE)	Use R-PE and R-PE-based tandems after the Click-iT [®] detection reaction.
TC-FIAsH™/TC-ReAsH™ reagents	Detect the tetracysteine (TC) tag with TC-FlAsH [™] or TC-ReAsH [™] reagents before the Click-iT [®] detection reaction.

^{*}Compatibility indicates whether fluorescent molecule itself or the detection methods involve components that are unstable in the presence of copper catalyst used for the Click-iT $^{\otimes}$ EdU detection reaction. Not all anti-GFP antibodies recognize the same antigen site. Rabbit and chicken anti-GFP antibodies perform well, but the mouse monoclonal antibodies tested do not generate an acceptable amount of fluorescence and are not recommended for this application.

Materials required but not provided

- 1% bovine serum albumin (BSA) in phosphate buffered saline (PBS), pH 7.1–pH 7.4
- Buffered saline solution, such as PBS, D-PBS, or TBS
- Deionized water or 18 MΩ purified water
- 12×75 -mm tubes, or other flow cytometry tubes

Cautions

- DMSO (Component C), provided as a solvent in this kit, is known to facilitate the entry of organic molecules into tissues. Handle reagents containing DMSO using equipment and practices appropriate for the hazards posed by such materials. Dispose of the reagents in compliance with all pertaining local regulations.
- Click-iT® fixative (Component D) contains paraformaldehyde, which is harmful. Use with appropriate precautions.
- Click-iT[®] saponin-based permeabilization and wash reagent (Component E) contains sodium azide, which yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.

Preparing reagents

- 1.1 Allow vials to warm to room temperature before opening.
- 1.2 To prepare a 10 mM solution of EdU, add 8 mL of DMSO (Component C) or aqueous solution (PBS) to Component A and mix well. After use, store any remaining stock solution at <-20°C. When stored as directed, the stock solution is stable for up to 1 year.
- 1.3 To prepare a working solution of Alexa Fluor® 647 azide (Cat. no. C10419) or Alexa Fluor® 488 azide (Cat. no. C10420), add 260 µL of DMSO to Component B and mix well. After use, store any remaining working solution at ≤−20°C. When stored as directed, this working solution is stable for up to 1 year.
- 1.4 To prepare 1 mL of 1X Click-iT® saponin-based permeabilization and wash reagent, add 100 mL of Component E to 900 mL of 1% BSA in PBS. Smaller amounts can be prepared by diluting a volume of Component E 1:10 with 1% BSA in PBS. After use, store any remaining solutions at 2°C-8°C. When stored as directed, the 1X solution is stable for 6 months and the 10X solution is stable for 12 months after receipt.

Note: Component E contains sodium azide (see Cautions, page 4).

1.5 To make a 10X stock solution of the Click-iT® EdU buffer additive (Component G), add 2 mL of deionized water to the vial and mix until the Click-iT® EdU buffer additive is fully dissolved. After use, store any remaining stock solution at ≤-20°C. When stored as directed, the stock solution is stable for up to 1 year.

The following protocol was developed with Jurkat cells, a human T cell line, and using an EdU concentration of 10 µM, and can be adapted for any cell type. Growth medium, cell density, cell type variations, and other factors may influence labeling. In initial experiments, we recommend testing a range of EdU concentrations to determine the optimal concentration for your cell type and experimental conditions. If currently using a BrdU based assay for cell proliferation, a similar concentration to BrdU is a good starting concentration for EdU. If using whole blood as the sample, we recommend heparin as the anticoagulant for collection.

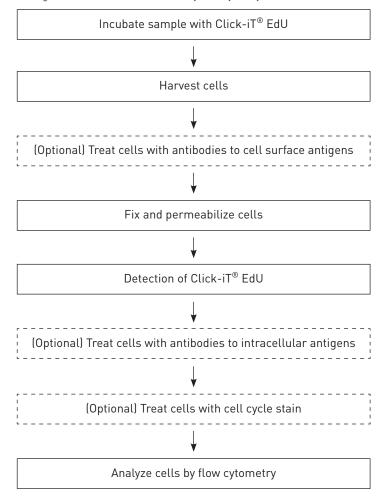


Figure 3 Workflow diagram for the Click-iT® EdU Flow Cytometry Assay Kits.

Labeling cells with EdU

- 2.1 Suspend the cells in an appropriate tissue culture medium to obtain optimal conditions for cell growth. Disturbing the cells by temperature changes or washing prior to incubation with EdU slows the growth of the cells during incorporation.
- 2.2 Add EdU to the culture medium at the desired final concentration and mix well. We recommend a starting concentration of 10 µM for 1-2 hours. For longer incubations, use lower concentrations. For shorter incubations, higher concentrations may be required. For a negative staining control, include cells from the same population that have not been treated with EdU.

- 2.3 Incubate under conditions optimal for cell type for the desired length of time. Altering the amount of time the cells are exposed to EdU or subjecting the cells to pulse labeling with EdU allows the evaluation of various DNA synthesis and proliferation parameters. Effective time intervals for pulse labeling and the length of each pulse depend on the cell growth rate.
- 2.4 Harvest cells and proceed immediately to step 3.1 if performing antibody surface labeling; otherwise continue to step 4.1.

Staining cell-surface antigens with antibodies (optional)

- 3.1 Wash cells once with 3 mL of 1% BSA in PBS, pellet cells by centrifugation, and remove supernatant.
- 3.2 Dislodge the pellet and resuspend cells at 1×10^7 cells/mL in 1% BSA in PBS.
- 3.3 Add 100 µL of cell suspension or whole blood sample to flow tubes.
- 3.4 Add surface antibodies and mix well (Table 2, page 3).

Note: Do not use PE, PE-tandem, or Qdot® antibody conjugates before performing the click reaction; wait until step 6.1 for labeling with these fluorophores.

- 3.5 Incubate for the recommended time and temperature, protected from light.
- 3.6 Proceed to step 4.1 for cell fixation.

Fixation and permeabilization

The Click-iT® saponin-based permeabilization and wash reagent can be used with whole blood or cell suspensions containing red blood cells, as well as with cell suspensions containing more than one cell type. This permeabilization and wash reagent maintains the morphological light scatter characteristics of leukocytes while lysing red blood cells.

- 4.1 Wash the cells once with 3 mL of 1% BSA in PBS, pellet the cells, and remove the supernatant.
- 4.2 Dislodge the pellet, add 100 μL of Click-iT[®] fixative (Component D), and mix well.
- 4.3 Incubate the cells for 15 minutes at room temperature, protected from light.
- 4.4 Wash the cells with 3 mL of 1% BSA in PBS, pellet the cells, and remove the supernatant. Repeat the wash step if red blood cells or hemoglobin are present in the sample. Remove all residual red blood cell debris and hemoglobin before proceeding.
- 4.5 Dislodge the cell pellet and resuspend the cells in 100 µL of 1X Click-iT® saponin-based permeabilization and wash reagent (prepared in step 1.4), and mix well. Incubate the cells for 15 minutes or proceed directly to step 5.1 for click labeling.

Click-iT® reaction

- 5.1 Prepare 1X Click-iT® EdU buffer additive by diluting the 10X stock solution (prepared in step 1.5) 1:10 in deionized water.
- 5.2 Prepare the Click-iT[®] reaction cocktail according to Table 3 (page 7).

Table 3 Click-iT® reaction cocktails

Reaction	Number of reactions							
components	1	2	5	10	15	30	50	100
PBS, D-PBS, or TBS	438 µL	875 µL	2.19 mL	4.38 mL	6.57 mL	13.2 mL	21.9 mL	43.8 mL
CuSO ₄ (Component F)	10 μL	20 µL	50 μL	100 μL	150 µL	300 µL	500 μL	1 mL
Fluorescent dye azide (prepared in step 1.3)	2.5 µL	5 μL	12.5 µL	25 µL	37.5 μL	75 µL	125 µL	250 μL
Reaction Buffer Additive (prepared in step 5.1)	50 μL	100 μL	250 μL	500 μL	750 μL	1.5 mL	2.5 mL	5 mL
Total reaction volume	500 μL	1 mL	2.5 mL	5 mL	7.5 mL	15 mL	25 mL	50 mL

Note: Use the Click-iT[®] reaction cocktail within 15 minutes of preparation.

- 5.3 Add 0.5 mL of Click-iT[®] reaction cocktail to each tube and mix well.
- 5.4 Incubate the reaction mixture for 30 minutes at room temperature, protected from light.
- 5.5 Wash the cells once with 3 mL of 1X Click-iT[®] saponin-based permeabilization and wash reagent (prepared in step 1.4), pellet the cells, and remove the supernatant. Dislodge the cell pellet and resuspend the cells in 100 µL of 1X Click-iT® saponin-based permeabilization and wash reagent, if proceeding with intracellular antibody labeleing in step 6.1. Otherwise, add 500 µL of 1X Click-iT® saponin-based permeabilization and wash reagent and proceed with step 7.1 for staining the cells for DNA content, or with step 8.1 for analyzing the cells on a flow cytometer.

Staining intracellular or surface antigens (optional)

- 6.1 Add antibodies against intracellular antigens or against surface antigens that use RPE, PE-tandem, or Qdot® antibody conjugates. Mix well.
- 6.2 Incubate the tubes for the time and temperature required for antibody staining, protected from light.
- 6.3 Wash each tube with 3 mL of 1X Click-iT® saponin-based permeabilization and wash reagent (prepared in step 1.4), pellet the cells, and remove the supernatant. Dislodge the cell pellet and resuspend the cells in 500 µL of 1X Click-iT® saponin-based permeabilization and wash reagent, and proceed with step 7.1 for staining the cells for DNA content, or with step 8.1 for analyzing the cells on a flow cytometer.

7.1 If necessary, add Ribonuclease A to each tube and mix (Table 4).

Table 4 Click-iT® EdU compatibility with DNA content stains

	Click-iT® EdU St		
DNA Content Stain	Alexa Fluor® 647	Alexa Fluor® 488	RNase required?
FxCycle [™] Violet	Yes	Yes	No
Propidium iodide (PI)	Yes	No*	Yes
SYTOX [®] AADvanced [™]	Yes	Yes	Yes
FxCycle [™] Far Red	No	Yes	Yes

^{*} You can use propidium iodide (PI) for cell cycle analysis with Alexa Fluor® 488 dye; however, this may reduce the Click-iT® Alexa Fluor® 488 EdU signal. Other DNA stains are recommended

Note: The complete excitation and emission spectra for the DNA content stains are available at www.lifetechnologies.com/spectraviewer.

7.2 Add the appropriate DNA stain to each tube, mix well, and incubate as recommended for each DNA stain.

Analysis by flow cytometry


If measuring total DNA content on a traditional flow cytometer using hydrodynamic focusing, use a low flow rate during acquisition. If using the Attune® Acoustic Focusing Cytometer, all collection rates may be used without loss of signal integrity if the event rate is kept below 10,000 events per second. However, for each sample within an experiment, the same collection rate and cell concentration should be used. The fluorescent signal generated by DNA content stains is best detected with linear amplification. The fluorescent signal generated by Click-iT® EdU labeling is best detected with logarithmic amplification.

- 8.1 Analyze the cells using a flow cytometer.
 - For the detection of EdU with Alexa Fluor® 647 azide use 633/635 nm excitation with a red emission filter (660/20 nm or similar).
 - For the detection of EdU with Alexa Fluor® 488 azide, use 488 nm excitation with a green emission filter (530/30 nm or similar).

Follow these guidelines to perform dual labeling of cultured cells by combining EdU with BrdU labeling.

- Use EdU for the first pulse and BrdU for the second pulse.
- Removal of EdU from the cell culture media is not required when BrdU is added as the second label.
- Addition of BrdU to culture media containing EdU results in preferential incorporation of BrdU into the DNA with the exclusion of EdU, while simultaneous addition of EdU with equimolar or half equimolar BrdU to the media results in only BrdU incorporation. This simplifies the dual labeling protocol by eliminating the wash step normally required to remove the first label from the culture media prior to addition of the second label.
- Process the cells after dual pulse labeling using a proven BrdU protocol.
- After the DNA denaturation step in the BrdU protocol, click label the cells first for the detection of EdU, and then follow with an antibody labeling protocol for the detection of BrdU.
- Select a BrdU antibody which does not have cross-reactivity to EdU, such as clone MoBu-1 (Cat. nos. B35129, B35139, B35140, B35141).

Figure 4 Dual pulse labeling with EdU and BrdU. Jurkat cells were first treated with 20 mM EdU for 1 hour. Without washing or removal of the EdU, BrdU was added at a 10 µM concentration, and the cells were incubated for 1 hour. The cells were harvested, washed, fixed with 70% ice-cold ethanol and stored at 4°C for 96 hours, followed by an HCl-based denaturation procedure. EdU was detected with Alexa Fluor® 488 azide using the Click-iT® EdÚ Flow Cytometry Kit (Cat. no. C10420). BrdU was then detected with anti-BrdU, Alexa Fluor® 647 conjugate (Cat. no. A21305). SYTOX® Blue nucleic acid stain (Cat. no. S11348) with RNase was used to detect DNA content. The labeled cells were analyzed by flow cytometry using 488 nm excitation with a 530/30 nm bandpass to detect EdU, 633 nm excitation with a 660/20 nm bandpass to detect BrdU, and 405 nm excitation with a 450/50 nm bandpass to detect DNA content. Four populations of cells are distinguished in the EdU vs BrdU plot: cells which are positive for both (Q2, upper right), cells which are negative for both (Q3, lower left), EdU positive but BrdU negative (Q4, lower right), and cells which are positive for BrdU but negative for EdU (Q1, upper left).

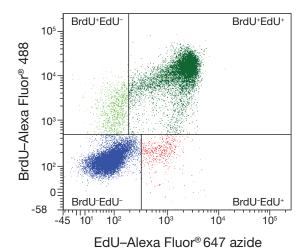

References

1. Chembiochem 4, 1147 (2003); 2. J Am Chem Soc 125, 3192 (2003); 3. Angew Chem Int Ed Engl 41, 2596 (2002); 4. Angew Chem Int Ed Engl 40, 2004 (2001); **5.** BioTechniques 44, 927 (2008); **6.** Curr Protoc Cytom 55, 7.38.1 (2011).

Product List Current prices may be obtained from our website or from our Customer Service Department.

Cat no.	Product name	Unit size
C10419	Click-iT [®] EdU Alexa Fluor [®] 647 Flow Cytometry Assay Kit *100 assays*	1 kit
C10420	Click-iT® EdU Alexa Fluor® 488 Flow Cytometry Assay Kit *100 assays*	1 kit
Related prod	ducts	
C10418	Click-iT [®] EdU Pacific Blue [™] Flow Cytometry Assay Kit *50 assays*	1 kit
C10424	Click-iT® EdU Alexa Fluor® 647 Flow Cytometry Assay Kit *50 assays*	1 kit
C10425	Click-iT® EdU Alexa Fluor® 488 Flow Cytometry Assay Kit *50 assays*	1 kit
C10632	Click-iT® Plus EdU Alexa Fluor® 488 Flow Cytometry Assay Kit *50 assays*	1 kit
C10633	Click-iT® Plus EdU Alexa Fluor® 488 Flow Cytometry Assay Kit *100 assays*	1 kit
C10634	Click-iT® Plus EdU Alexa Fluor® 647 Flow Cytometry Assay Kit *50 assays*	1 kit
C10635	Click-iT® Plus EdU Alexa Fluor® 647 Flow Cytometry Assay Kit *100 assays*	1 kit
A10044	EdU (5-ethynyl-2´-deoxyuridine)	50 mg
B35129	BrdU mouse monoclonal antibody (Clone MoBU-1) - Pacific Blue [™] *for flow cytometry* *100 tests*	1 each
B35139	BrdU mouse monoclonal antibody (Clone MoBU-1) Alexa Fluor® 488 *for flow cytometry* *100 tests*	1 each
B35140	BrdU mouse monoclonal antibody (Clone MoBU-1) Alexa Fluor® 647 *for flow cytometry* *100 tests*	1 each
B35141	BrdU mouse monoclonal antibody (Clone MoBU-1) unconjugated *for flow cytometry* *100 tests*	1 each
F10347	FxCycle [™] Violet Stain *for flow cytometry* *500 assays*	1 kit
F10348	FxCycle [™] Far Red Stain *for flow cytometry* *500 assays*	1 kit
F10797	FxCycle [™] PI/RNase Staining Solution *200 assays*	100 mL
H3570	Hoechst® 33342, trihydrochloride, trihydrate *10 mg/mL solution in water*	10 mL
P3566	Propidium Iodide – 1.0mg/ml solution in water	10 mL
S10349	SYTOX® AADvanced™ dead cell stain *for 488 excitation* *for flow cytometry* *100 tests*	1 kit
V35003	Vybrant® DyeCycle [™] Violet stain *5 mM in water* *200 assays*	200 μL
12091-021	RNase A (20 mg/mL)	10 mL
14190-144	Dulbecco's Phosphate Buffered Saline 1X without Calcium Chloride without Magnesium Chloride	500 mL
14190-250	Dulbecco's Phosphate Buffered Saline 1X without Calcium Chloride without Magnesium Chloride	x 500 mL

Purchaser notification

Corporate headquarters

5791 Van Allen Way Carlsbad, CA 92008

USA

Phone: +1 760 603 7200 Fax: +1 760 602 6500

Email: techsupport@lifetech.com

European headquarters

Inchinnan Business Park 3 Fountain Drive Paisley PA4 9RF UK Phone: +44 141 814 6100

Toll-Free Phone: 0800 269 210
Toll-Free Tech: 0800 838 380
Fax: +44 141 814 6260
Teh Fax: +44 141 814 6117
Email: euroinfo@invitrogen.com
Email Tech: eurotech@invitrogen.com

Japanese headquarters

LOOP-X Bldg. 6F 3-9-15, Kaigan Minato-ku, Tokyo 108-0022 Japan Phone: +81 3 5730 6509 Fax: +81 3 5730 6519

Email: jpinfo@invitrogen.com

Additional international offices are listed at www.lifetechnologies.com

These high-quality reagents and materials must be used by, or directly under the supervision of, a technically qualified individual experienced in handling potentially hazardous chemicals. Read the Safety Data Sheet provided for each product; other regulatory considerations may apply.

Obtaining support

For the latest services and support information for all locations, go to www.lifetechnologies.com.

At the website, you can:

- Access worldwide telephone and fax numbers to contact Technical Support and Sales facilities
- Search through frequently asked questions (FAQs)
- Submit a question directly to Technical Support (techsupport@lifetech.com)
- Search for user documents, SDSs, vector maps and sequences, application notes, formulations, handbooks, certificates of analysis, citations, and other product support documents
- · Obtain information about customer training
- · Download software updates and patches

SDS

Safety Data Sheets (SDSs) are available at www.lifetechnologies.com/sds.

Certificate of Analysis

The Certificate of Ánalysis provides detailed quality control and product qualification information for each product. Certificates of Analysis are available on our website. Go to www.lifetechnologies.com/support and search for the Certificate of Analysis by product lot number, which is printed on the product packaging (tube, pouch, or box).

Limited product warranty

Life Technologies Corporation and/or its affiliate(s) warrant their products as set forth in the Life Technologies' General Terms and Conditions of Sale found on Life Technologies' website at www.lifetechnologies.com/termsandconditions. If you have any questions, please contact Life Technologies at www.lifetechnologies.com/support.

Disclaimer

LIFE TECHNOLOGIES CORPORATION AND/OR ITS AFFILIATE(S) DISCLAIM ALL WARRANTIES WITH RESPECT TO THIS DOCUMENT, EXPRESSED OR IMPLIED, INCLUDING BUT NOT LIMITED TO THOSE OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE, OR NON-INFRINGEMENT. TO THE EXTENT ALLOWED BY LAW, IN NO EVENT SHALL LIFE TECHNOLOGIES AND/OR ITS AFFILIATE(S) BE LIABLE, WHETHER IN CONTRACT, TORT, WARRANTY, OR UNDER ANY STATUTE OR ON ANY OTHER BASIS FOR SPECIAL, INCIDENTAL, INDIRECT, PUNITIVE, MULTIPLE OR CONSEQUENTIAL DAMAGES IN CONNECTION WITH OR ARISING FROM THIS DOCUMENT, INCLUDING BUT NOT LIMITED TO THE USE THEREOF.

Important licensing information

These products may be covered by one or more Limited Use Label Licenses. By use of these products, you accept the terms and conditions of all applicable Limited Use Label Licenses.

Hoechst is a registered trademark of Hoechst GmbH. All other trademarks are the property of Thermo Fisher Scientific and its subsidiaries.

Life Technologies is a Thermo Fisher Scientific brand. © 2014 Thermo Fisher Scientific Inc. All rights reserved.

