
thermoscientific



Higher success and simpler setup with a trusted PCR enzyme

Thermo Scientific™ Phusion™ Plus DNA Polymerase is the latest addition to our family of Phusion™ products, which are *Pyrococcus* polymerase–like high-fidelity enzymes fused to a DNA-binding domain. This unique protein fusion technology enables Phusion Plus DNA Polymerase to generate PCR sequences with high accuracy, sensitivity, and inhibitor tolerance. In addition, its specialized reaction buffer allows PCR amplification without the need to calculate primer annealing temperatures, helping you save time and avoid mistakes in PCR runs.

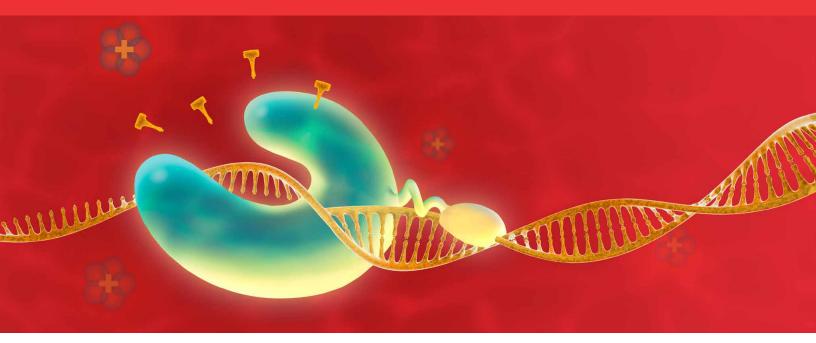

Highlights

- High fidelity—provides >100x higher sequence accuracy than Tag polymerase
- Convenient setup—allows for a universal annealing temperature of 60°C
- Enhanced specificity—reduces nonspecific amplification and primer degradation via hot-start modification
- Benchtop stability—enables lab automation setup since assembled reactions are stable at room temperature for 24 hours
- High inhibitor tolerance—works with DNA of suboptimal purity



thermo scientific

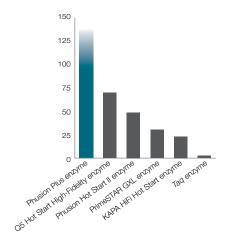


Figure 1. Fidelity relative to *Taq* DNA polymerase. Fidelity values were determined by next-generation sequencing with molecular barcodes. Phusion Plus DNA Polymerase generates PCR amplicons with very few errors, due to its extremely high fidelity >100x that of *Taq* enzyme.

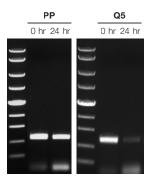


Figure 4. Reaction stability on the benchtop. Assembled PCR reactions were loaded immediately (0 hr) onto a thermal cycler or incubated at room temperature for 24 hours before thermal cycling. With a stringent hot start, of Phusion Plus DNA Polymerase (PP) produced PCR amplicons with high specificity and yields even after 24 hours, in contrast to NEB Q5™ Hot Start High-Fidelity DNA Polymerase (Q5).

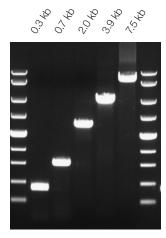


Figure 2. Specific and efficient amplification using a universal annealing temperature. Targets of different lengths were amplified efficiently from human genomic DNA (gDNA) using Phusion Plus DNA Polymerase and an annealing temperature of 60°C. Calculated annealing temperatures of the primers range from 61°C to 69°C.

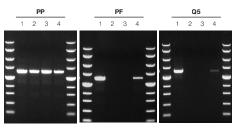


Figure 5. Enzyme tolerance to common PCR inhibitors. A 2 kb target was amplified from human gDNA, with the reactions containing: (1) no inhibitor; (2) humic acid (0.5 μg/mL); (3) hemin (2.5 μM); (4) xylan (250 μg/mL). Phusion Plus DNA Polymerase (PP) amplified the targets more efficiently than Phusion™ Hot Start Flex DNA Polymerase (PF) and Q5 Hot Start High-Fidelity DNA Polymerase (Q5) in the presence of the inhibitors.


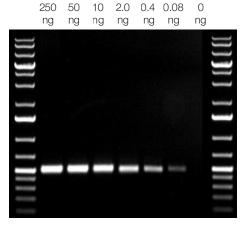


Figure 3. Highly sensitive and specific PCR. A 0.5 kb target can be detected from as little as 0.08 ng of human gDNA using Phusion Plus DNA Polymerase.







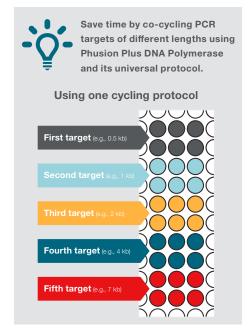

Visit our molecular biology resource library for scientific webinars, videos, articles, and more at

thermofisher.com/mbresources

Ordering information

Description	Size*	Cat. No.
Phusion Plus DNA Polymerase	100 reactions	F630S
Phusion Plus PCR Master Mix	100 reactions	F631S
Phusion Plus Green PCR Master Mix	100 reactions	F632S

^{*} Additional product sizes available.



