

Corning® Cryogenic Storage Solutions

A new and improved way to freeze your cells

CORNING




A New Standard in Cell Cryopreservation

At Corning, we continuously look for ways to help our customers improve or streamline steps in the cell culture workflow. One such area is cell cryopreservation. While current methods exist, they require chemicals and maintenance.

Now there is a new and improved way to freeze cells for cryogenic storage – we call it Corning® CoolCell®.

Corning CoolCell is an alcohol-free cell freezing container, which controls the rate of freezing to $-1^{\circ}\text{C}/\text{minute}$ when placed in a -80°C freezer. CoolCell has been performance tested with a variety of cell types including stem cells, primary cells, PBMC cell lines, insect cells, and yeast. The CoolCell technology utilizes a thermo-conductive alloy core and highly insulative outer material to control the rate of heat removal and provide reproducible cell cryopreservation. CoolCell units are easy to use and deliver comparable results.

Corning CoolCell Reproducibility

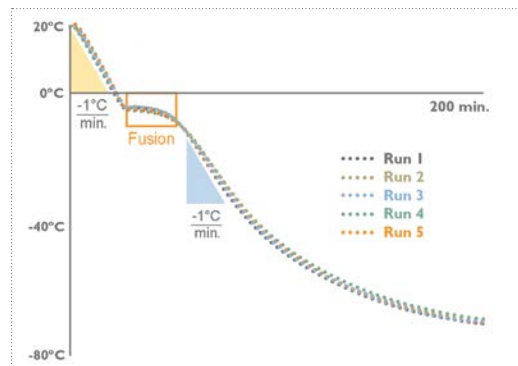


Performance test: A temperature probe was placed into a 2.0 mL cryogenic vial containing 1.0 mL of water and the tube was inserted into a room temperature Corning CoolCell. The CoolCell was placed into a -80°C freezer and the temperature rate and profile was recorded over a 3-hour period. The test was repeated 5 consecutive times.

Conclusion: Corning CoolCell generated identical fusion time and cooling profiles over five consecutive freeze cycles.


Alcohol-free with No Ongoing Costs or Maintenance

Isopropanol (IPA) containers used for cryogenic freezing require costly alcohol replacements every 5 uses, can be cumbersome to handle, and may have inconsistent freezing rates. Corning CoolCell is different, because it's a reusable, alcohol-free way to uniformly freeze your cells at a lower cost of use. With CoolCell, you can depend on high reproducibility and high cell viability, to ensure you preserve the most cells possible for your research.

Corning CoolCell Performance Versus IPA Container

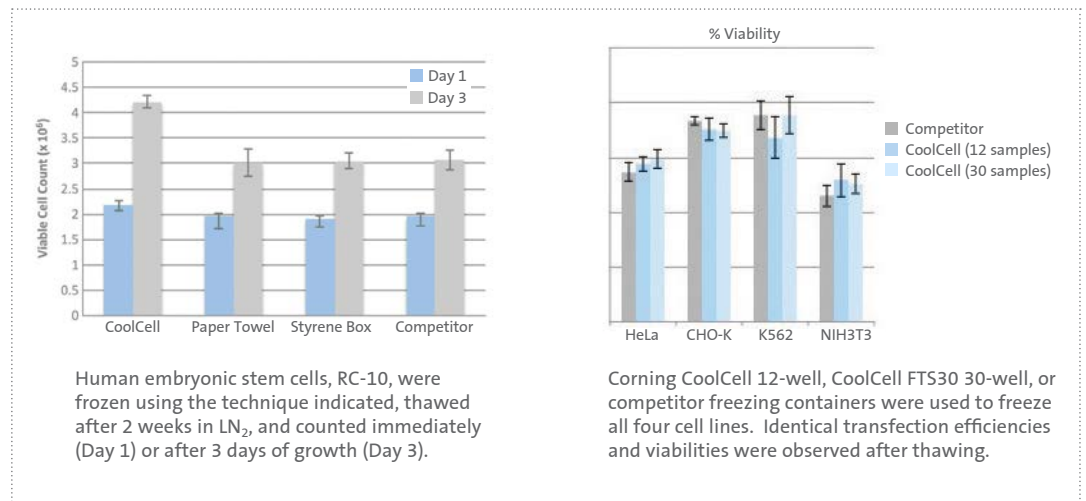

How Corning CoolCell LX Works

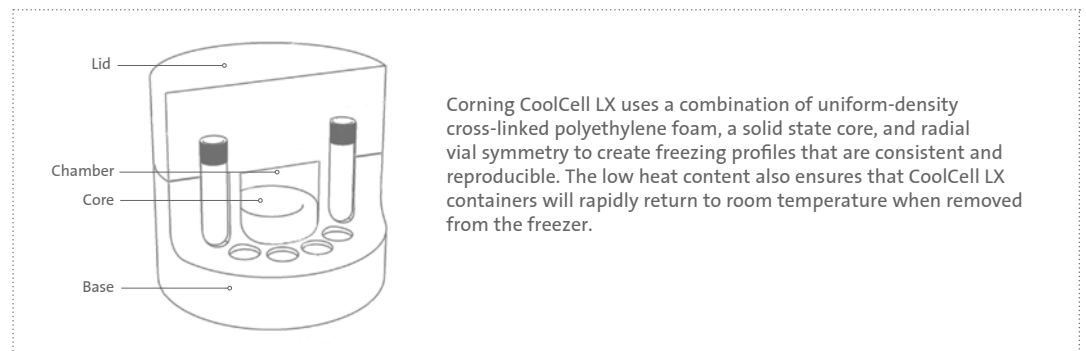



Corning 2D bar coded cryogenic tube, 2 mL, self-standing (Cat. Nos. 8670, 8671)



DMSO media (Cat. No. 25-950-CQC)



Cryopreservation bags

Corning® CoolCell® Features

Unique features of controlled-rate freezing with Corning CoolCell include:

- ▶ Ease of use
- ▶ Alcohol and fluid-free freezing
- ▶ Lower cost of use than alcohol-based devices
- ▶ High cell recovery and cell viability
- ▶ Reproducibility
- ▶ Simple, consistent way to standardize controlled-rate freezing

Keep Your Samples Safe

Use Corning cryogenic vials and DMSO media with Corning CoolCell to further protect your valuable cell lines, biological, and aqueous solutions in ultra-low temperature storage. Choose from external or internal cryogenic vial thread caps or assorted color cap options to suit your needs. For added convenience, Corning also offers reusable cryogenic racks.

New Improve Your Sample Management

Manage and manipulate multiple storage tubes more efficiently with Corning 1D/2D bar coded cryogenic vials. Our cryogenic vials have a permanent 2D bar code on the bottom and a standard linear 1D bar code on the side of the vial. Corning 1D/2D bar coded cryogenic vials are temperature-resistant polypropylene vials that can withstand temperatures down to -196°C and are compatible with most scanning and capper/decapper systems.

Further Protect Your Valuable Cells with a Complete Cryopreservation Solution

Combine **Corning CoolCell** with **Corning cryogenic vials, grippers, and DMSO** to further protect your valuable cell lines, biological, and aqueous solutions in ultra-low temperature storage.

Also Available for Cell Freezing

Corning Cryopreservation Bags

Corning also offers cryogenic storage containers designed for the storage, preservation, and transfer of cells. Features include: a unique bag film material that remains flexible at low temperatures and proprietary port designs that allow for increased flexibility. Learn more at www.corning.com/lifesciences.



Ordering Information

Corning® CoolCell® Containers

Cat. No.	Description	Capacity (Vials)	Exposed Vial Tops	Qty/Pk	Qty/Cs
432000	CoolCell, purple	12	No	1	1
432001	CoolCell LX, purple	12	Yes	1	1
432002	CoolCell LX, green	12	Yes	1	1
432003	CoolCell LX, orange	12	Yes	1	1
432004	CoolCell LX, pink	12	Yes	1	1
432138	CoolCell LX, 4 colors (purple, green, orange, pink)	12	Yes	1	4
432005	CoolCell 5 mL LX, purple	12	Yes	1	1
432006	CoolCell FTS30, purple	30	Yes	1	1
432007	CoolCell FTS30, orange	30	Yes	1	1
432008	CoolCell FTS30, green	30	Yes	1	1
432009	CoolCell FTS30, pink	30	Yes	1	1
432010	CoolCell SV2	12	Yes	1	1
432011	CoolCell SV10	6	Yes	1	1

Corning CoolCell Accessories

Cat. No.	Description	Capacity (Vials)	Exposed Vial Tops	Qty/Pk	Qty/Cs
432076	CoolCell filler vials, 2.0 mL	—	—	6	6
432077	CoolCell filler vials, 5.0 mL	—	—	6	6
432078	CoolCell FTS30 vial module	30	—	10	10
432136	Cryogenic vial grippers, multi-colored	—	—	5	5

Corning Cryogenic Vials and Accessories

External Thread Cryogenic Vials

Cat. No.	Capacity (mL)	Style	Self-standing	Qty/Pk	Qty/Cs
8671	2.0	1D and 2D bar coded, round bottom	Yes	50	500
8676	2.0	1D bar coded, round bottom	Yes	50	500
430658	1.2	Conical bottom	Yes	50	500
430659	2.0	Round bottom	Yes	50	500
430661	2.0	Round bottom	No	50	500
430662	4.0	Round bottom	Yes	50	500
430663	5.0	Round bottom	Yes	50	500

Internal Thread Orange Cap Cryogenic Vials

8670	2.0	1D and 2D bar coded, round bottom	Yes	50	500
8672	2.0	1D bar coded, round bottom	Yes	50	500
430487	1.2	Conical bottom	Yes	50	500
430488	2.0	Round bottom	Yes	50	500
430489	2.0	Round bottom	No	50	500
430490	4.0	Round bottom	No	50	500
430491	4.0	Round bottom	Yes	50	500
430656	5.0	Round bottom	Yes	50	500

Warning! Do not use cryogenic vials for storage in the liquid phase of liquid nitrogen. Only store vials in the vapor phase above the liquified gas. Always use appropriate safety equipment when removing vials from cryogenic storage.

Cryogenic Storage Boxes

Cat. No.	Description	Qty/Pk	Qty/Cs
8673	Cryogenic storage box, polycarbonate, holds 81 vials, designed to fit Corning 2D bar coded cryogenic vials	5	10
8674	Cryogenic storage box, polycarbonate, holds 100 vials, designed to fit Corning 2D bar coded cryogenic vials	5	10

Cryopreservation Bags

Cat. No.	Size (mL)	Fill Volume (mL)	Qty/Cs
91-200-88	50	10 - 20	1
91-200-89	250	30 - 70	1
91-200-90	500	55 - 100	1
91-200-91	750	80 - 190	1

For more specific information on claims, visit www.corning.com/certificates.

Warranty/Disclaimer: Unless otherwise specified, all products are for research use only.* Not intended for use in diagnostic or therapeutic procedures. Not for use in humans. Corning Life Sciences makes no claims regarding the performance of these products for clinical or diagnostic applications. *Cryogenic vials are considered US Class I medical devices.

For additional product or technical information, visit www.corning.com/lifesciences or call 800.492.1110. Outside the United States, call +1.978.442.2200 or contact your local Corning sales office.

CORNING

Corning Incorporated
Life Sciences

836 North St.
Building 300, Suite 3401
Tewksbury, MA 01876
t 800.492.1110
t 978.442.2200
f 978.442.2476

www.corning.com/lifesciences

ASIA/PACIFIC
Australia/New Zealand
t 61 427286832

Chinese Mainland
t 86 21 3338 4338
f 86 21 3338 4300

India
t 91 124 4604000
f 91 124 4604099

Japan
t 81 3-3586 1996
f 81 3-3586 1291

Korea
t 82 2-796-9500
f 82 2-796-9300

Singapore
t 65 6572-9740
f 65 6735-2913

Taiwan
t 886 2-2716-0338
f 886 2-2516-7500

EUROPE
CSEurope@corning.com

France
t 0800 916 882
f 0800 918 636

Germany
t 0800 101 1153
f 0800 101 2427

The Netherlands
t 020 655 79 28
f 020 659 76 73

United Kingdom
t 0800 376 8660
f 0800 279 1117

All Other European Countries
t +31 (0) 206 59 60 51
f +31 (0) 206 59 76 73

LATIN AMERICA
grupoLA@corning.com

Brazil
t 55 (11) 3089-7400

Mexico
t (52-81) 8158-8400