TaqMan[®] Gene Expression Cells-to-CT[™] Kit

Two-Step RT-PCR Directly from Cultured-Cell Lysates

Part Numbers 4399002, AM1728, AM1729



TaqMan® Gene Expression Cells-to-CT™ Kit

(P/N 4399002, AM1728, AM1729)

Protocol

I.	Introduction
	A. Product Description and Background
	B. Materials Provided with the Kit and Storage Conditions
	C. Materials Not Provided with the Kit
	D. Related Products Available from Applied Biosystems
II.	TaqMan® Gene Expression Cells-to-CT™ Procedure
	A. Suggested Controls and Replicates
	B. Cell Lysis
	C. Reverse Transcription (RT)
	D. Real-Time PCR
III.	Troubleshooting
	A. No PCR Product or Unexpected PCR Products
	B. RT-PCR Products in the Negative Control Reactions
IV.	Appendix
	A. Pilot Experiment
	B. Safety Information
	C. Quality Control

P/N 4385117 Revision D

For research use only. Not for use in diagnostic procedures. By use of this product, you accept the terms and conditions of all applicable Limited Label Licenses. For statement(s) and/or disclaimer(s) applicable to this product, see below.

Revision Date: August 23, 2010

Information in this document is subject to change without notice. Applied Biosystems assumes no responsibility for any errors that may appear in this document.

Applied Biosystems disclaims all warranties with respect to this document, expressed or implied, including but not limited to those of merchantability or fitness for a particular purpose. In no event shall Applied Biosystems be liable, whether in contract, tort, warranty, or under any statute or on any other basis for special, incidental, indirect, punitive, multiple or consequential damages in connection with or arising from this document, including but not limited to the use thereof.

Literature Citation When describing a procedure using this product for publication, please refer to it as the TaqMan $^{\circ}$ Gene Expression Cells-to-CT $^{\infty}$ Kit.

Warranty and Liability Applied Biosystems is committed to delivering superior product quality and performance, supported by industry-leading global service and technical support teams. Warranty information for the accompanying consumable product is available at www.ambion.com/info/warranty in "Limited Warranty for Consumables," which is subject to the exclusions, conditions, exceptions, and limitations set forth under the caption "EXCLUSIONS, CONDITIONS, EXCEPTIONS, AND LIMITATIONS" in the full warranty statement. Please contact Applied Biosystems if you have any questions about our warranties or would like information about post-warranty support.

Notice to Purchaser: Limited License A license to perform the patented 5' Nuclease Process for research is obtained by the purchase of (i) both Authorized 5' Nuclease Core Kit and Licensed Probe, (ii) a Licensed 5' Nuclease Kit, or (iii) license rights from Applied Biosystems.

This product is an Authorized 5' Nuclease Core Kit. Use of this product is covered by US patent claims and patent claims outside the US. The purchase of this product includes a limited, non-transferable immunity from suit under the foregoing patent claims for using only this amount of product for the purchaser's own internal research. Separate purchase of a Licensed Probe would convey rights under the applicable claims of US Patents, and corresponding claims outside the United States. No right under any other patent claim and no right to perform commercial services of any kind, including without limitation reporting the results of purchaser's activities for a fee or other commercial consideration, is conveyed expressly, by implication, or by estoppel. This product is for research use only. Diagnostic uses under Roche patents require a separate license from Roche. Further information on purchasing licenses may be obtained from the Director of Licensing, Applied Biosystems, 850 Lincoln Centre Drive, Foster City, California 94404, USA.

Trademarks Applied Biosystems, AB (Design), Ambion, GeneAmp, MicroAmp and RNase*Zap* are registered trademarks, and Cells-to-CT, DNA*Zap*, MultiScribe, ROX, StepOne, StepOnePlus, Veriti, FAM and XenoRNA are trademarks of Applera Corporation or its subsidiaries in the US and/or certain other countries. TaqMan is a registered trademark of Roche Molecular Systems, Inc. All other trademarks are the sole property of their respective owners.

I. Introduction

A. Product Description and Background

Reverse transcription (RT)-PCR with real-time detection of amplification products is a robust, simple, and quantitative way to measure mRNA levels in biological samples. Traditionally, the first step in gene expression experiments has been to recover pure RNA from experimental samples. Even using the quickest and simplest techniques, however, RNA isolation is fairly time-consuming, often requiring 30 minutes or more of hands-on sample manipulation. Furthermore, with small samples, it can lead to loss of RNA. Cells-to-CT™ technology enables reverse transcription of lysates from 10–10⁵ cultured cells without isolating or purifying RNA. Real-time PCR analysis is carried out directly afterwards. Eliminating the RNA isolation step substantially expedites and simplifies gene expression analysis of cultured cells.


Cells-to-CT lysates exhibit sensitivity and specificity similar to that from purified RNA in real-time RT-PCR using TaqMan® Gene Expression Assays. The lysis procedure simultaneously prepares cell lysates for RT-PCR and removes genomic DNA in under ten minutes. The lysis step is simple to automate with robotic platforms for high-throughput processing of 96- or 384-well plates because it takes place entirely at room temperature. The procedure is also economical; there are only a few pipetting steps and, with cells grown in 96- or 384-well plates, no sample transfers.

Procedure overview

The TaqMan $^{\circ}$ Gene Expression Cells-to-CT $^{\text{TM}}$ Kit procedure is shown in Figure $\underline{1}$ and described below.

- First, 10–10⁵ cultured cells are washed with phosphate-buffered saline (PBS), mixed with Lysis Solution, and incubated at room temperature for 5 min. Cells are lysed during this incubation and RNA is released into the Lysis Solution which contains reagents to inactivate endogenous RNases. If DNase I is added to the Lysis Solution (optional), genomic DNA is also degraded at this step.
- Next, Stop Solution is mixed into the lysate to inactivate the lysis reagents so that they will not inhibit the RT or PCR.
- Cell lysates are reverse transcribed to synthesize cDNA using a convenient 20X RT Enzyme Mix and 2X RT Buffer.
- Finally, the cDNA is amplified by real-time PCR using the included TaqMan Gene Expression Master Mix and the TaqMan Gene Expression Assay for your target-of-interest (assays sold separately).

Figure 1. TaqMan® Gene Expression Cells-to-CT™ Procedure Overview

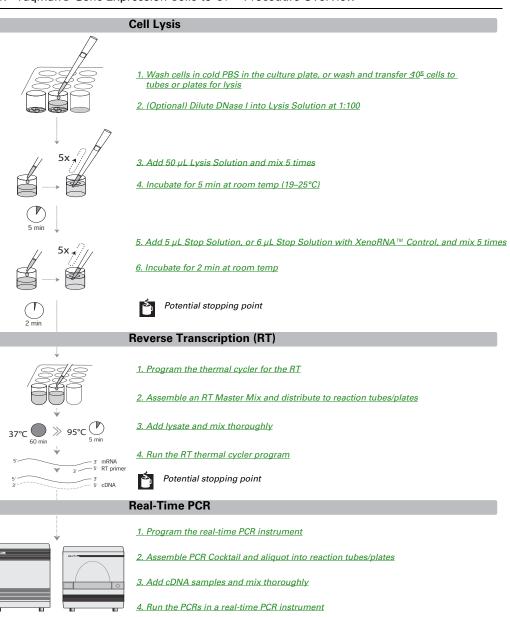

TaqMan® Cells-to-CT™ Control Kit

Applications

The TaqMan Cells-to-CT Control Kit (P/N 4386995, available separately) is designed for use with the TaqMan Gene Expression, TaqMan PreAmp, and TaqMan Fast Cells-to-CT Kits. It includes XenoRNA™ Control, an RNA transcript that has no homology to any known sequences, and a XenoRNA TaqMan Gene Expression Assay to detect the XenoRNA Control target. It also includes a TaqMan Assay for the highly expressed endogenous control gene, β-actin. These assays are designed to be run in parallel (not in multiplex PCR) with assay(s) for the gene(s)-of-interest, to serve as positive controls for RT and real-time PCR.

The TaqMan® Gene Expression Cells-to-CT™ Kit can be used in any application where real-time RT-PCR is used to analyze mRNA from cultured cells. It was developed using a process that included extensive testing for sensitivity and specificity with a broad selection of TaqMan Gene Expression Assays. The kit is well suited for large experiments, such as real-time PCR analysis of gene expression in differentially treated cell cultures or RNAi screening experiments using siRNA to modulate gene expression. Other applications include screening a library of compounds for their effects on mRNA expression, following the regulation of mRNA as cells are treated with increasing concentrations of a particular chemical, or evaluating the expression of mRNA in time course experiments.

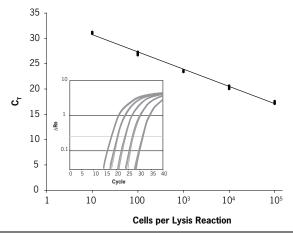


Figure 2. Real-time RT-PCR Using the TaqMan® Gene Expression Cells-to- $C\tau^{TM}$ Kit

A dilution series of $10{\text -}10^5$ HeLa cells was processed in triplicate with the TaqMan* Gene Expression Cells-to-CTT* Kit. The endogenous control gene β -actin was amplified from the cDNA in triplicate 20 μ L reactions. The inset shows amplification plots, and the main figure shows the threshold cycle (C_T) compared to the input number of cells. Amplification was linear over a cell input range of $10{\text -}10^5$ cells per lysis.



Cell type compatibility

Table 1 lists the cell types that have been shown to be compatible with the Cells-to-CT technology.

Table 1. Cell Types Compatible with Cells-to-CT™ Technology

Cell Line	Growth	Source Species	Source Tissue
A549	adherent	H. sapiens	Lung Carcinoma
ВЈ	adherent	H. sapiens	Foreskin fibroblast
CHO-K1	adherent	C. griseus (hamster)	Ovary
COS-7	adherent	C. aethiops (monkey)	Kidney
DU-145	adherent	H. sapiens	Prostate Carcinoma
HEK-293	adherent	H. sapiens	Kidney
HeLa	adherent	H. sapiens	Cervical Adenocarcinoma
HepG2	adherent	H. sapiens	Liver Carcinoma
Huh-7	adherent	H. sapiens	Liver carcinoma
Jurkat	suspension	H. sapiens	Acute T-Cell Leukemia
K-562	suspension	H. sapiens	Chronic Myelogenous Leukemia
ME-180	adherent	H. sapiens	Cervical Epidermoid Carcinoma
NCI-H460	adherent	H. sapiens	Large Cell Lung Cancer
Neuro 2A	adherent	M. musculus (mouse)	Brain blastoma
NIH/3T3	adherent	M. musculus (mouse)	Embryonic Fibroblast
PC-12	adherent	R. norvegicus (rat)	Adrenal Pheochromocytoma
Primary Hepatocytes	adherent	H. sapiens	Liver
PT-K75	adherent	S. scrofa (pig)	Nasal Turbinate Mucosa
Raji	suspension	H. sapiens	B Lymphocyte
SK-N-AS	adherent	H. sapiens	Brain Neuroblast
SK-N-SH	adherent	H. sapiens	Brain Fibroblast
U-87 MG	adherent	H. sapiens	Brain Glioblastoma
U-2 OS	adherent	H. sapiens	Bone osteosarcoma

B. Materials Provided with the Kit and Storage Conditions

The TaqMan® Gene Expression Cells-to-CT™ Kit is available in sizes for lysis of 40, 100, or 400 samples. For each lysis reaction, the kit includes RT and PCR reagents for the following numbers of reactions:

- 40 reaction size (P/N 4399002): For each lysis reaction, RT reagents for one reverse transcription reaction (RT) and one minus-RT control are included, and Gene Expression Master Mix for either 40, 50 μL real-time PCR assays or 100, 20 μL reactions.
- 100 reaction size (P/N AM1728): For each lysis reaction, RT reagents for one reverse transcription reaction (RT) and one minus-RT control are included, and Gene Expression Master Mix for either 200, 50 μL real-time PCR assays or 500, 20 μL reactions.

400 reaction size (P/N AM1729): For each lysis reaction, RT reagents for one reverse transcription reaction (RT) and one minus-RT control are included, and Gene Expression Master Mix for either 800, 50 μL real-time PCR assays or 2000, 20 μL reactions.

Note that TaqMan Gene Expression Assays, which contain the primers and probe needed for reverse transcription and real-time PCR are not included with the TaqMan Gene Expression Cells-to-CT Kit.

P/N 4399002 40 rxn	P/N AM1728 100 rxn	P/N AM1729 400 rxn	Component	Storage
200 μL	500 μL	2 x 1 mL	Stop Solution	–20°C
22 µL	55 μL	220 μL	DNase I	-20°C
110 µL	275 μL	1.1 mL	20X RT Enzyme Mix	–20°C
2.2 mL	5.5 mL	22 mL	Lysis Solution	4°C
2.2 mL	5.5 mL	22 mL	2X RT Buffer	4°C
1 mL	5 mL	4 x 5 mL	TaqMan® Gene Expression Master Mix	4°C



The product is shipped on dry ice, but once received, the kit components should be stored at the temperatures indicated in the table above. Storage at -20°C should be in a non-frost-free freezer.

C. Materials Not Provided with the Kit

Lab equipment and supplies

- General laboratory equipment including vortex mixer, microcentrifuge, and pipettors
- Nuclease-free pipette tips, nuclease-free microcentrifuge tubes and/or U-bottom 96-well plates (for cells *not* grown in 96- or 384-well culture plates)
- Real-time PCR tubes or 96-well plates appropriate for your instrument
- Thermal cycling instrument: The TaqMan Cells-to-CT Kits were developed using Applied Biosytems thermal cyclers (for the RT) and real-time PCR instruments, however the technology is expected to be compatible with any thermal cycler for the RT and any real-time PCR system compatible with the fluorophores used in TaqMan Gene Expression Assays. TaqMan Gene Expression Cells-to-CT reactions can be run on any currently available Applied Biosystems instrument platforms, including the following:

Reverse transcription

- GeneAmp® PCR System 9700
- Veriti[™] 96-Well Thermal Cycler

TaqMan® Gene Expression Cells-to-CT™Kit

Real-Time PCR

- Applied Biosystems 7900 HT Fast Real-Time PCR System (in "standard" mode)
- Applied Biosystems 7500 Real-Time PCR System
- StepOne[™] or StepOnePlus[™] Real-Time PCR Systems

Reagents

- RT-PCR grade water, for example Ambion P/N AM9935
- Phosphate-buffered saline (PBS), for example diluted from Ambion 10X PBS, P/N AM9624 and AM9625
- TaqMan* Gene Expression Assay for real-time PCR of the target(s) of interest (Applied Biosystems P/N 4331182, 4351372).
- (optional) TaqMan[®] Cells-to-CT[™] Control Kit (P/N 4386995)

D. Related Products Available from Applied Biosystems

TaqMan [®] Cells-to-CT [™] Control Kit P/N 4386995	Designed for use with TaqMan Cells-to-Ct Kits, the Control Kit includes XenoRNA $^{\text{\tiny M}}$ Control, an RNA transcript that has no homology to any known sequence, and a corresponding XenoRNA TaqMan Gene Expression Assay. It also includes a TaqMan Assay for the highly expressed endogenous control gene β -Actin. These assays are designed to be run in parallel with assay(s) for the gene(s)-of-interest to serve as positive controls for reverse transcription and real-time PCR.
TaqMan [®] Gene Expression Assays P/N 4331182 P/N 4351372	Pre-designed and pre-optimized TaqMan* probe and primer sets for quantitative gene expression analysis using real-time PCR. Each TaqMan Assay is a ready-to-use 20X mixture of PCR primers and TaqMan probe designed for amplification using universal cycling conditions so that any combination of assays can be run using the same thermal cycling conditions. There are over 700,000 inventoried assays for human, mouse, rat, <i>Arabidopsis</i> , <i>Drosophila</i> , <i>C. elegans</i> , <i>C. familiaris</i> (dog), and <i>Rhesus macaque</i> genes, and Custom TaqMan Gene Expression Assays are also available.
Applied Biosystems Real-Time PCR Instruments See web or print catalog for P/Ns	Since pioneering real-time PCR, Applied Biosystems has continued to develop the technology to provide more powerful solutions for labs of all sizes. The Applied Biosystems family of real-time platforms, which includes the Applied Biosystems 7300, 7500, 7500 Fast, 7900HT Fast, StepOne™, and StepOnePlus™ Real-Time PCR Systems, provide cutting-edge tools while making real-time PCR more accessible than ever. These systems are easy to use with next generation software and the flexibility to run the real-time chemistry of your choice.
RNase <i>Zap</i> [®] Solution P/N AM9780, AM9782, AM9784	RNaseZap* RNase Decontamination Solution is simply sprayed, poured, or wiped onto surfaces to instantly inactivate RNases. Rinsing twice with distilled water will eliminate all traces of RNase and RNaseZap Solution.
RNase-free Tubes & Tips See web or print catalog for P/Ns	Ambion® RNase-free tubes and tips are available in most commonly used sizes and styles. They are guaranteed RNase- and DNase-free. See our latest catalog or our website (www.ambion.com/prod/tubes) for specific information.

RT-PCR Grade Water P/N AM9935	Ambion* RT-PCR Grade Water is certified free of nucleases, and free of nucleic acid contamination that may cause false-positive signals in RT-PCR. The RT-PCR Grade Water is ready to use and requires no preparation, mixing, or autoclaving. Like all Ambion water products, RT-PCR Grade Water is deionized, autoclaved, and sterile filtered. Each lot is subjected to 2 rounds of rigorous quality control testing before being certified nuclease-free.
Plastic Consumables for PCR See web or print catalog for P/Ns	Applied Biosystems MicroAmp® disposables are optimized to provide unmatched temperature accuracy and uniformity for fast, efficient PCR amplification. The plastics are compatible with various Applied Biosystems platforms, including real-time PCR instruments. Choose your format based on throughput requirements or your thermal cycler block.
DNA <i>Zap</i> ™ Solution P/N AM9890	DNAZap™ PCR DNA degradation solution consists of two solutions that are innocuous separately but, when mixed, are able to degrade high levels of contaminating DNA and RNA from surfaces instantly.

II. TaqMan® Gene Expression Cells-to-C[™] Procedure

A. Suggested Controls and Replicates

Biological and technical replicates

The definitions of biological and technical replicates can have different interpretations; here we define biological replicates as separate lysis reactions and technical replicates as separate PCRs. Our minimum and recommended number of replicates are the following:

- Biological replicates: 2–4
- Technical (PCR) replicates: 2-4

Negative controls for RT-PCR

Minus-RT control

Minus-RT controls contain all the RT reaction components except the 20X RT Enzyme Mix (substitute water); they demonstrate that the template for the PCR was cDNA, and not genomic DNA.



Minus-RT controls may not be necessary if the cDNA will be amplified using a TaqMan[®] Gene Expression Assay with the suffix "_m*"; these assays are designed to span an exon-exon junction and will not amplify genomic DNA.

No-template control

No-template controls contain all the PCR components except the cell lysate (substitute water). If the no-template control yields a fluorescent signal, it indicates that the RT or PCR reagents were contaminated with DNA, for example; PCR product from previous reactions.

B. Cell Lysis



NOTE

These instructions describe using fresh cultured cells. If your experiment requires frozen cultured cells, we recommend washing the cells in cold PBS as described in step 1 (below) before freezing. When you are ready to start the procedure, allow the cells to thaw on ice, and begin the procedure at step 2.

Before you start

- Thaw Stop Solution, invert or flick the tube several times to mix thoroughly (do *not* vortex), and place on ice.
- Chill 1X PBS to 4°C.



MPORTANT

The maximum number of cells that can be used in this procedure varies somewhat according to cell type, but is generally 10⁵ cells. Instructions for a pilot experiment to determine the maximum number of cells per reaction is outlined in section <u>IV.A.</u> starting on page 17. Using too many cells per lysis reaction may result in incomplete lysis and/or inhibition of RT-PCR.

Adherent cells grown in 96- or 384-well plates

It is important to start with cells that have been cultured until they are fully adhered to the plate, otherwise they will be lost during washing.

- a. Estimate (or count) the number of cells per well. Make sure that each well contains $10-10^5$ cells.
- b. Aspirate and discard the culture medium from the wells.
- c. Add 50 µL of cold (4°C) 1X PBS to each well.
- d. Aspirate the PBS from the well. Remove as much PBS as possible without disturbing the cells. *Proceed to step 2.*

Cells grown in other vessels (adherent and suspension cells)

- a. Adherent cells only (for suspension cells, start at step b): Detach cells using the subculturing method routinely employed in your laboratory for the cell type. If trypsin is used, inactivate it before proceeding.
- b. Count, then gently pellet the cells, aspirate and discard the growth medium, and place cells on ice.
- c. Wash cells in cold (4°C) PBS by resuspending them in ~0.5 mL PBS per 10⁶ cells (or ~50 μL PBS for ≤0⁵ cells). Gently pellet the cells, then aspirate and discard as much of the PBS as possible without disturbing the cell pellet. Place the cells on ice.
- d. Resuspend cells in fresh, cold 1X PBS so that 5 μ L will contain the desired number of cells for a single lysis reaction (10–10⁵ cells/lysis).



Alternatively, cells can be split into aliquots for individual lysis reactions before washing in PBS.

- e. Split the cell suspension into individual lysis reactions: distribute $5~\mu L$ of the suspension to wells of a U-bottom multiwell plate or microcentrifuge tubes. Place the cells on ice.
- 2. (Optional) Dilute DNase I into Lysis Solution at 1:100

To remove genomic DNA during cell lysis, determine the volume of Lysis Solution needed for the experiment, and dilute DNase I into the Lysis Solution at 1:100 for use in the next step. Include ~10% overage; examples are shown in Table 2 below:

Table 2. Preparation of Lysis Solution Containing DNase I

Component	each rxn	96 rxns	384 rxns
Lysis Solution	49.5 µL	5.23 mL	20.91 mL
DNase I	0.5 μL	52.8 μL	211 μL

- 3. Add 50 µL Lysis Solution and mix 5 times
- a. Add 50 µL Lysis Solution (with or without DNase I) to each sample.
- b. Mix the lysis reaction by pipetting up and down 5 times. To avoid bubble formation, mix with the pipettor set at 35 μ L and expel the solution without emptying the pipette tip completely.



NOTE

Lysis Solution and Cells-to-C⊤ lysates may appear cloudy at room temperature—this is expected.

- 4. Incubate for 5 min at room temp (19-25°C)
- Add 5 µL Stop Solution, or 6 µL Stop Solution with XenoRNA™ Control, and mix 5 times



NOTE

To dispense Stop Solution using a multichannel pipettor, pipet Stop Solution from a set of strip tubes or adjacent wells in a 96-well plate rather than from a reagent reservoir. If desired, unused Stop Solution can be stored frozen in these tubes when properly sealed.

Incubate the lysis reactions for 5 min at room temperature (19–25°C).

- a. (Optional) To include an endogenous control using the TaqMan Cells-to-CT Control Kit, prepare Stop Solution containing XenoRNA™ Control for the samples in the experiment at hand. Add 1 µL XenoRNA Control per 5 µL Stop Solution.
- b. Pipet 5 μL of Stop Solution alone or 6 μL of Stop Solution with XenoRNA Control directly into each lysis reaction.
 - Do this by touching the surface of the lysate with the opening of the pipet tip to ensure that all of the Stop Solution is added to the lysate.
- c. Mix the lysis reaction by pipetting up and down 5 times. To avoid bubble formation, mix with the pipettor set at 35 μ L and expel the solution without emptying the pipette tip completely.



IMPORTANT

It is very important to thoroughly mix the Stop Solution into the lysate.

TagMan® Gene Expression Cells-to-CT™ Procedure

6. Incubate for 2 min at room temp

Incubate for 2 min at room temp (19-25°C).

Do not allow Cells-to-CT lysates to remain at room temp for longer than 20 min after adding the Stop Solution.



STOPPING POINT

Lysates can be stored on ice for \$\alpha\$ hr, or at -20°C or -80°C for \$\alpha\$ months.

C. Reverse Transcription (RT)

1. Program the thermal cycler for the RT

Program the thermal cycler as shown in Table 3.

Table 3. Thermal Cycler Settings for RT

	Stage	Reps	Temp	Time
Reverse transcription (hold)	1	1	37°C	60 min
RT inactivation (hold)	2	1	95°C	5 min
Hold	3	1	4°C	indefinite

2. Assemble an *RT Master Mix* and distribute to reaction tubes/plates



If a 50 μ L RT reaction will not provide enough material for all of the planned PCR assays, the RT reaction can be scaled up proportionally. a. Calculate the number of RT reactions in the experiment. Using Table 4 below, assemble an *RT Master Mix* for all the reactions plus ~10% overage in a nuclease-free microcentrifuge tube on ice.

Table $\underline{4}$ shows the reaction setup using the recommended amount of lysate in each RT reaction, 10 μ L. If desired, up to 45% of the RT reaction volume (22.5 μ L) can be Cells-to-CT lysate; adjust the volume of Nuclease-free Water accordingly.

Table 4. RT Master Mix: 50 µL Final Reaction Volume

Component	Each rxn	96 rxns	384 rxns
2X RT Buffer	25 µL	2.64 mL	10.56 mL
20X RT Enzyme Mix*	2.5 µL	264 µL	1.056 mL
Nuclease-free Water	12.5 µL	1.32 mL	5.28 mL
final volume RT master mix	40 µL	4.22 mL	16.9 mL

^{*} For the minus-RT control, use Nuclease-free Water in place of 20X RT Enzyme Mix.

- Mix gently, but thoroughly, centrifuge briefly, then place the RT Master Mix on ice.
- Distribute RT Master Mix to nuclease-free PCR tubes or wells of a multiwell plate.
- 3. Add lysate and mix thoroughly
- a. Add sample lysate to each aliquot of RT Master Mix for a final 50 μL reaction volume.
- b. Once assembled, mix reactions gently, then centrifuge briefly to collect the contents at the bottom of the reaction vessel.





Assembled RT reactions can be stored at 4°C for up to 4 hr.

4. Run the RT thermal cycler program

Using a thermal cycler (or real-time PCR instrument), incubate at 37°C for 60 min, then at 95°C for 5 min to inactivate the RT enzyme.



STOPPING POINT

Completed RT reactions may be stored at -20°C.

D. Real-Time PCR

1. Program the real-time PCR instrument

Program the real-time PCR instrument as in Table 5.



IAADODTA NIT

On Applied Biosystems real-time PCR instruments capable of Fast mode thermal cycling, select Standard mode.

- TaqMan® Gene Expression Master Mix contains ROX™ passive reference dye.
- Specify the fluorescent dye(s) used in the TaqMan Gene Expression Assays for the experiment.

The ACTB and XenoRNA™ Gene Expression Assays in the TaqMan Cells-to-CT™ Control Kit are labeled with FAM™ dye and a nonfluorescent quencher.

Table 5. Real-Time PCR Cycling Conditions

	Stage	Reps	Temp	Time
UDG Incubation (hold)	1	1	50°C	2 min
Enzyme Activation (hold)	2	1	95°C	10 min
PCR (cycle)	3	40	95°C	15 sec
			60°C	1 min

2. Assemble *PCR Cocktail* and aliquot into reaction tubes/plates

a. Calculate the number of PCR assays in the experiment. Using Table 6 below, assemble a *PCR Cocktail* for all the reactions plus ~10% overage in a nuclease-free microcentrifuge tube at room temp. Table 6 shows the reaction setup using the recommended amount of cDNA in each PCR, 20%. If desired, up to 45% of the PCR volume can be cDNA; adjust the quantity of Nuclease-free Water accordingly.

TaqMan® Gene Expression Cells-to-CT™ Procedure



Before use, mix the TaqMan Gene Expression Master Mix by swirling the bottle. Mix TaqMan Assays by vortexing briefly or flicking the tube a few times, and then centrifuging.

Table 6. PCR Cocktail Setup

	20 μL PCRs	50 μL PCRs
Component	Each rxn	Each rxn
TaqMan® Gene Expression Master Mix (2X)	10 μL	25 µL
TaqMan® Gene Expression Assay* (20X)	1 μL	2.5 µL
Nuclease-free water	5 μL	12.5 µL
final volume PCR cocktail	16 µL	40 µL

^{*} Not included in the TaqMan Gene Expression Cells-to-CT Kit.



Be sure to use reaction containers and lids that are compatible with your real-time PCR instrument.

- 3. Add cDNA samples and mix thoroughly
- b. Distribute the PCR Cocktail into individual PCR tubes or wells of a real-time PCR plate at room temp.
- a. Add a portion of the RT reaction (cDNA sample) to each aliquot of PCR Cocktail as indicated in Table Z.

Table 7. PCR Setup

Component	20 μL PCRs	50 μL PCRs
PCR Cocktail	16 µL	40 μL
RT Reaction (cDNA)	4 μL	10 μL

- b. Cover the plate or close the tubes, and mix gently. Then centrifuge briefly to remove bubbles and collect the contents at the bottom of the wells/tubes.
- 4. Run the PCRs in a real-time PCR instrument

Place the reactions in a real-time PCR instrument and start the run using the settings programmed in step $\underline{D.1}$ on page 12.

Refer to your real-time PCR instrument guide for information on evaluating the data.

III. Troubleshooting

A. No PCR Product or Unexpected PCR Products

Problems with adding or mixing the Stop Solution

Components in the Lysis Solution may inhibit RT-PCR if they are not fully inactivated by the Stop Solution.

- Be sure to add the Stop Solution directly to the lysate, in other words, touch the lysate with the opening of the pipet tip when adding the Stop Solution to make sure that the entire 5 or 6 μ L of Stop Solution is added to each sample.
- Also, we recommend mixing by pipetting up and down five times.

RNA was degraded before starting the procedure

To avoid RNA degradation, keep cells in PBS on ice before starting the cell lysis procedure. Take cells off ice just prior to adding Lysis Solution.

RNase in the sample was not completely inactivated

Too many cells were used in the lysis reaction

If too many cells per sample are used in the procedure, the RNase in the sample may not be totally inactivated and/or cellular components or debris could inhibit reverse transcription or PCR.

- Generally ≤105 cells can be used successfully in the Cells-to-CT[™] procedure, but if RT or PCR fails, try using fewer cells (e.g., 5–10-fold fewer cells).
- Also, consider performing a pilot experiment to determine the optimal number of cells for your cell type (section <u>IV.A</u> on page 17).

Too much PBS was left on the cells, diluting the Lysis Solution If >5 μ L of PBS remains in samples when the Lysis Solution is added, the Lysis Solution may be too dilute to fully inactivate cellular RNases. To avoid this, remove as much PBS as possible before adding Lysis Solution to the cells, or if you split the cells after the PBS wash, resuspend cells in μ L PBS for each sample of $10-10^5$ cells.

Lysates sat too long before going into RT

Do not allow lysates to sit longer than 20 min at room temperature once the Stop Solution has been added: either freeze the lysates at -20° C or -80° C, or start the RT reactions. Alternatively, lysates can be safely stored on ice for up to 2 hr after lysis.

The sample does not contain the target RNA

Negative results are often difficult to confirm as valid. Consider running the following experiments before concluding that the sample does not contain the RNA of interest:

Verify that the TaqMan® Gene Expression Cells-to-CT™ procedure
is working by including XenoRNA™ Control (from the TaqMan
Cells-to-CT Control Kit, sold separately) in the sample in
step II.B.5.a on page 10. Then use the XenoRNA TaqMan Gene
Expression Assay to amplify a XenoRNA target following the

- instructions in section <u>II.D</u> on page 12. If product is generated in the XenoRNA amplification, but no product is seen in the PCR for the gene of interest, then it is possible that the RNA of interest is not expressed in the cells and/or is undetectable with this procedure.
- For experiments with samples consisting of <100 cells per lysis, it may be desirable to verify that each sample did, in fact, contain cells. Check that samples contained cells with intact RNA by real-time RT-PCR with a TaqMan Gene Expression Assay for a highly-expressed endogenous control such as β -actin using the ACTB assay included in the TaqMan Cells-to-CT Control Kit (sold separately).
- Check that the PCR for your target works with your PCR primers, reagents, and equipment by using cDNA generated from purified RNA from the same source (or a similar one) in PCR. If the amplification does not give good results using cDNA from purified RNA, it will not work with Cells-to-CT lysate.

B. RT-PCR Products in the Negative Control Reactions

PCR products in the no-template PCR control

PCR products in the no-template PCR control indicate that DNA contamination of the sample occurred. Contamination of PCR reagents, pipettors, and benchtops with DNA should be considered.

- Careful laboratory practices are essential to avoid contaminating reactions with PCR products. Keep concentrated DNA solutions (PCR products, plasmid prep, etc.) away from the area where PCRs are assembled. Clean the lab bench and the pipettors routinely with the Ambion DNAZap™ Solution (P/N AM9890) or another DNA decontamination product. Use barrier tips to pipette PCR reagents, and store completed PCRs in a different location from the PCR reagents. Unfortunately, the only way to remedy contaminated reagents is to replace them.
- It is always a good idea to routinely include a no-template negative control reaction with experimental PCRs. If no-template controls routinely yield PCR products, more stringent steps may be taken to control contamination.

PCR products in the minus-RT control

If PCR products are seen in the minus-RT control reaction, but not in the no-template control, it indicates that genomic DNA remains in the samples and that genomic DNA was amplified in real-time PCR. Consider the following suggestions:

 Be sure to mix DNase I into the Lysis Solution thoroughly in step <u>II.B.2</u> on page 10 for efficient removal of genomic DNA from samples.

TaqMan® Gene Expression Cells-to-CT™Kit

- Consider using a TaqMan Gene Expression Assay designed to span an exon-exon boundary; such assays are designated with _m* at the end of the part number.
- Use fewer cells per lysis reaction.
- Lyse cells (step <u>II.B.3</u> on page 10) using Lysis Solution that is at room temperature, and make sure that lysis reactions occur at room temperature (19–25°C).
 - If PCR products are still seen in minus-RT control reactions, try the following two suggestions in the order shown:
 - Increase the lysis reaction incubation time in step $\underline{\text{II.B.4}}$ on page 10 to 8 minutes.
 - –Use Lysis Solution that has been warmed up to 25°C for cell lysis (step II.B.3 on page 10).

IV. Appendix

A. Pilot Experiment

The purpose of this pilot experiment is to identify the maximum number of cells to use in TaqMan® Gene Expression Cells-to- $CT^{\text{\tiny TM}}$ reactions. Using too many cells can result in inefficient cell lysis and RT-PCR inhibition, and the maximum number of cells varies somewhat according to the cell type.

In this simple experiment, cells are serially diluted and lysed following the normal protocol. If you have the TaqMan Cells-to-CT Control Kit (P/N 4386995), we recommend adding XenoRNA™ Control to the Stop Solution used to prepare Cells-to-CT lysates. The lysates are then subjected to real-time RT-PCR for an endogenous control gene, such as β-actin, (an ACTB TaqMan Gene Expression Assay is included in the TaqMan Cells-to-CT Control Kit), and the C_T values are plotted against the log of the number of cells in the lysis reaction. The resulting line will be linear for cell numbers that are compatible with the procedure and will deviate from linearity at concentrations that result in incomplete lysis or RT-PCR inhibition. The cDNA can be amplified in parallel using the XenoRNA TaqMan Gene Expression Assay (from the Cells-to-CT Control Kit). In contrast to the series of PCRs for an endogenous control gene, the number of cells in the lysis reaction should have no effect on the C_T value seen in XenoRNA amplification reactions because each reaction will contain the same amount of Xeno-RNA target.

- Harvest, count, and wash cells
- Follow the instructions for harvesting cells in step $\underline{\text{II.B.1}}$ on page 9, but follow the instructions in step $\underline{2}$ (next) for resuspending them in PBS.
- 2. Dilute cells to 2 x 10⁵ cells/µL

Prepare a cell suspension containing 2 x 10^5 cells/ μ L. Keep the cell suspension on ice.

- 3. Make 5 serial dilutions of the cells in 5-fold increments
- a. Prepare 5 tubes containing 45 μL of cold 1X PBS in ice.
- b. Transfer 5 μ L of the 2 x 10⁵ cells/ μ L to the first tube (1:5 dilution) and mix gently but thoroughly. Continue making the serial dilutions by transferring 5 μ L of each solution to the subsequent tube to finish with 5 suspensions containing 2 x 10⁴, 2000, 200, 20, and 2 cells per μ L.
- 4. Transfer cells to reaction containers with 3 replicates

Transfer $5 \,\mu L$ of each cell suspension to individual reaction tubes or wells of a multiwell plate. Include 3 biological replicates of each cell concentration.

The final cell counts will be 10^5 , 10^4 , 1000, 100, and 10 cells.

TagMan® Gene Expression Cells-to-CT™Kit

- 5. Lyse cells in 50 μL of Lysis Solution and incubate at room temp for 5 min
- a. Add 50 μL of Cells-to-CT Lysis Solution to each of the samples prepared in step $\underline{4}$.
- b. Mix thoroughly by pipetting up and down 5 times.
- c. Incubate at room temperature for 5 min.
- Add 5–6 µL Stop Solution, mix, and incubate at room temp for 2 min
- a. (Optional) Add XenoRNA Control to the Stop Solution. Mix 16.5 μ L XenoRNA Control with 82.5 μ L Stop Solution for use in the next step.
- b. Add 5 μL Stop Solution alone or 6 μL Stop Solution with XenoRNA Control to each cell lysate, and mix thoroughly by pipetting up and down 5 times.
- c. Incubate at room temp for 2 min.

7. Perform RT-PCR

Follow the protocols described in sections <u>II.C</u>–<u>D</u> to reverse transcribe and perform PCR using all samples.

- Amplify with any Applied Biosystems Endogenous Control TaqMan Gene Expression Assay, for example with the ACTB assay included in the Cells-to-CT Control Kit.
- If XenoRNA Control was added to samples, it can be amplified in a separate PCR using the XenoRNA TaqMan Gene Expression Assay included in the Cells-to-CT Control Kit.
- You can also evaluate real-time PCR of the cell titration using the TaqMan Gene Expression Assay for the target-of-interest to help determine the minimum number of cells required for its detection.

8. Evaluate results

Endogenous control

Create a plot of C_T versus the log of the number of cells in the lysis. The C_T values should decrease in a linear fashion as the number of cells increase, for cell numbers that are compatible with the procedure. When the number of cells per lysis reaction exceeds the capacity of the system, resulting in incomplete lysis or inhibition of RT-PCR, the data will not be linear. In future experiments, do not exceed the number of cells per lysis reaction that provided results within the linear range in the pilot experiment.

XenoRNA™ Control

The C_T values from the XenoRNA Control should be consistent (±1 C_T) regardless of the number of cells in the lysis reaction, indicating that no RT-PCR inhibitors are present in the Cells-to-CT lysate. If C_T values begin to increase at higher numbers of cells per lysis reaction, it indicates that inhibitors were introduced into RT-PCR with this number of cells. For future experiments, use only the number of cells per lysis reaction that did not show an increase in C_T value.

Target-of-interest

The pilot experiment can provide useful information about the number of cells required to detect the target-of-interest. Examine the results carefully and choose cell numbers that will provide sufficient signal for the experiment.

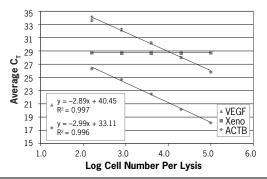


Figure 3. Example Pilot Experiment

HeLa cells were diluted and lysed following the instructions for the pilot experiment. Lysates were reverse transcribed and amplified using the indicated Taq-Man Gene Expression Assay in eight replicate reactions. The resulting $C_{\rm T}$ values were plotted against the log of the number of cells used in the lysis reaction (as recommended). This experiment shows that even at 10^5 cells per lysis reaction, RT-PCR is not compromised by reaction inhibitors carried over from the cell lysate. It also shows that the lowest number of cells/lysis reaction tested, 160 cells, is approaching the limit for detection of vascular endothelial growth factor (VEGF).

B. Safety Information

The MSDS for any chemical supplied by Applied Biosystems or Ambion is available to you free 24 hours a day.



IMPORTANT

For the MSDSs of chemicals not distributed by Applied Biosystems or Ambion, contact the chemical manufacturer.

To obtain Material Safety Data Sheets

- Material Safety Data Sheets (MSDSs) can be printed or downloaded from product-specific links on our website at the following address: www.ambion.com/techlib/msds
- Alternatively, e-mail your request to:
 MSDS_Inquiry_CCRM@appliedbiosystems.com. Specify the catalog
 or part number(s) of the product(s), and we will e-mail the associated
 MSDSs unless you specify a preference for fax delivery.

TaqMan® Gene Expression Cells-to-CT™Kit

 For customers without access to the internet or fax, our technical service department can fulfill MSDS requests placed by telephone or postal mail. (Requests for postal delivery require 1–2 weeks for processing.)

Chemical safety guidelines

To minimize the hazards of chemicals:

- Read and understand the Material Safety Data Sheets (MSDS) provided by the chemical manufacturer before you store, handle, or work with any chemicals or hazardous materials.
- Minimize contact with chemicals. Wear appropriate personal protective equipment when handling chemicals (for example, safety glasses, gloves, or protective clothing). For additional safety guidelines, consult the MSDS.
- Minimize the inhalation of chemicals. Do not leave chemical containers open. Use only with adequate ventilation (for example, fume hood). For additional safety guidelines, consult the MSDS.
- Check regularly for chemical leaks or spills. If a leak or spill occurs, follow the manufacturer's cleanup procedures as recommended on the MSDS.
- Comply with all local, state/provincial, or national laws and regulations related to chemical storage, handling, and disposal.

C. Quality Control

Functional testing

Lysis components are functionally tested by producing cDNA from a serial dilution of $10-10^5$ cultured cells, detecting PPIA using real-time PCR, and comparing amplification signal from Cells-to-CT lysates to those obtained from purified RNA from the same source.

Nuclease testing

Relevant kit components are tested in the following nuclease assays:

RNase activity

Meets or exceeds specification when a sample is incubated with 25 ng labeled RNA and analyzed by PAGE.

Nonspecific endonuclease activity

Meets or exceeds specification when a sample is incubated with 300 ng supercoiled plasmid DNA and analyzed by agarose gel electrophoresis.

Exonuclease activity

Meets or exceeds specification when a sample is incubated with 40 ng labeled *Sau3A* fragments of pUC19 and analyzed by PAGE.

Protease testing

Meets or exceeds specification when a sample is incubated with 1 μg protease substrate and analyzed by fluorescence.