

# Monkey TNF $\alpha$ Total ELISA Kit

Enzyme-linked Immunosorbent Assay for quantitative detection of monkey TNF  $\alpha$  total

**Catalog Numbers** BMS654 or BMS654TEN

**Pub. No.** MAN0016921 **Rev.** A.0 (30)

**WARNING!** Read the Safety Data Sheets (SDSs) and follow the handling instructions. Wear appropriate protective eyewear, clothing, and gloves. Safety Data Sheets (SDSs) are available from [thermofisher.com/support](http://thermofisher.com/support).

## Product description

The Monkey TNF $\alpha$  Total ELISA Kit is an enzyme-linked immunosorbent assay for the quantitative detection of monkey TNF  $\alpha$ .

## Summary

TNF  $\alpha$  is a multifunctional cytokine involved in many different pathways, in homeostasis and pathophysiology of mammals. It can show opposing biological effects suggesting complex regulatory mechanisms.

TNF  $\alpha$ , also known as cachectin, was first detected as a cytotoxic factor inducing lysis of certain tumor cells. The TNF  $\alpha$  gene is member 2 of the TNF-superfamily (consisting of at least 20 distinct members).

TNF  $\alpha$  release is mainly triggered by viral infections and endotoxins, lipopolysaccharides or other bacterial components, by tissue injury, DNA-damage and by IL-1, PDGF and TNF  $\alpha$  itself. It is primarily expressed in macrophages, but also in monocytes, neutrophils, NK-cells, mast-cells, endothelial cells and activated lymphocytes. TNF  $\alpha$  expression in endothelial cells and fibroblasts can be induced by IL-17. The expression of other cytokines, chemokines, reactive oxygen intermediates, nitric oxide and prostaglandins is stimulated by TNF  $\alpha$ .

The initially membrane bound TNF  $\alpha$  is enzymatically cleaved by TACE (=ADAM17). The soluble monomers aggregate to homotrimers and are secreted into blood and other biological fluids. The membrane bound and the soluble form are biologically active and bind to the TNF-receptors TNFR1 (=TNFRSF1A, p55-60) and TNFR2 (=TNFRSF1B, TNFBR2, p75-80).

Upon ligand binding, receptors form trimers leading to conformational changes, protein dissociation (SODD = silencer of death domains, BAG4, Bcl2-associated athanogene 4) and association (TRADD = TNF-R1 associated death domain protein) and yielding the following biological activities:

- transcription of anti-apoptotic factors and proteins involved in cell proliferation and inflammation via binding of TRAF2 (TNF-R associated factor 2) and RIPK1 (TNF-R interacting serine-threonine kinase 1) and activation of the transcription factor NF  $\kappa$  B.
- cell proliferation, differentiation but also apoptosis via TRAF2 binding, kinases activation, activation of c-Jun and ATF2 (JNK-MAPK-pathway).
- apoptosis via the binding of FADD (Fas associating protein with death domain) to TRADD and activation of caspases (including caspase 8 = FLICE).
- necrosis, a caspase independent cell death, mediated by NADPH oxidases, which form a complex with TRADD and RIPK1, leading to the generation of oxygen species.

TNF-R2 contains no DD (death domain), but exhibits its function via direct TRAF binding.

Thus the multiple biological functions of TNF  $\alpha$  comprise cellular proliferation and differentiation, tumorigenesis, apoptotic or necrotic

cell death (including certain tumor cell lines), immunoregulatory activities, lipid metabolism, coagulation and endothelial function. It promotes local or systemic inflammation (TNF  $\alpha$  is a potent pyrogen) and stimulates the acute phase response. Very high expressions of TNF  $\alpha$  after infection can lead to septic shock (TNF  $\alpha$  is highly cytotoxic.), whereas sustained low levels induce cachexia and inflammation.

Dysregulation of TNF  $\alpha$  is involved in many diseases.

For literature update, refer to our website.

## Principles of the test

An anti-monkey TNF  $\alpha$  coating antibody is adsorbed onto microwells.

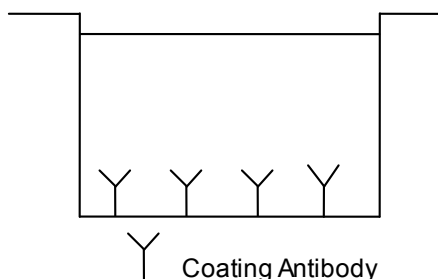


Fig. 1 Coated microwell

Monkey TNF  $\alpha$  present in the sample or standard binds to antibodies adsorbed to the microwells. A biotin-conjugated anti-monkey TNF  $\alpha$  antibody is added and binds to monkey TNF  $\alpha$  captured by the first antibody.

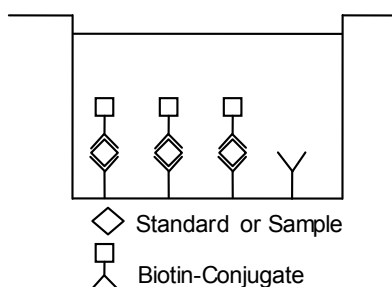


Fig. 2 First incubation

Following incubation unbound biotin-conjugated anti-monkey TNF  $\alpha$  antibody is removed during a wash step. Streptavidin-HRP is added and binds to the biotin-conjugated anti-monkey TNF  $\alpha$  antibody.

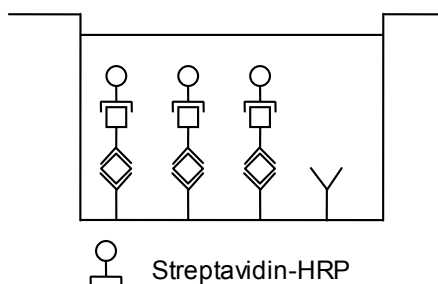


Fig. 3 Second incubation

Following incubation unbound Streptavidin-HRP is removed during a wash step, and substrate solution reactive with HRP is added to the wells.

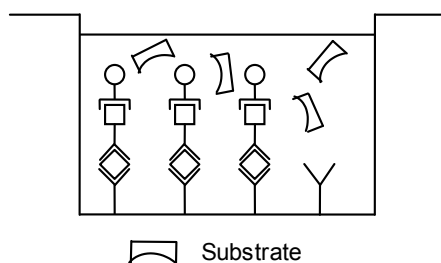


Fig. 4 Third incubation

A colored product is formed in proportion to the amount of monkey TNF  $\alpha$  present in the sample or standard. The reaction is terminated by addition of acid and absorbance is measured at 450 nm. A standard curve is prepared from 7 monkey TNF  $\alpha$  standard dilutions and monkey TNF  $\alpha$  sample concentration determined.

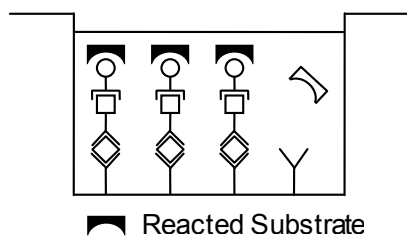


Fig. 5 Stop reaction

## Reagents provided

### Reagents for monkey TNF $\alpha$ total ELISA BMS654 (96 tests)

- 1 aluminum pouch with a Microwell Plate (12 strips of 8 wells each) coated with monoclonal antibody to monkey TNF  $\alpha$
- 1 vial (70  $\mu$ L) Biotin-Conjugate anti-monkey TNF  $\alpha$  polyclonal antibody
- 1 vial (150  $\mu$ L) Streptavidin-HRP
- 2 vials monkey TNF  $\alpha$  Standard lyophilized, 3600 U/mL upon reconstitution
- 1 vial (12 mL) Sample Diluent
- 1 vial (5 mL) Assay Buffer Concentrate 20x (PBS with 1% Tween<sup>™</sup> 20, 10% BSA)
- 1 bottle (50 mL) Wash Buffer Concentrate 20x (PBS with 1% Tween<sup>™</sup> 20)
- 1 vial (15 mL) Substrate Solution (tetramethyl-benzidine)
- 1 vial (15 mL) Stop Solution (1M Phosphoric acid)
- 4 Adhesive Films

### Reagents for monkey TNF $\alpha$ total ELISA BMS654TEN (10x96 tests)

- 10 aluminum pouches with a Microwell Plate (12 strips of 8 wells each) coated with monoclonal antibody to monkey TNF  $\alpha$
- 10 vials (70  $\mu$ L) Biotin-Conjugate anti-monkey TNF  $\alpha$  polyclonal antibody
- 10 vials (150  $\mu$ L) Streptavidin-HRP
- 10 vials monkey TNF  $\alpha$  Standard lyophilized, 3600 U/mL upon reconstitution
- 7 vials (12 mL) Sample Diluent
- 2 vials (5 mL) Assay Buffer Concentrate 20x (PBS with 1% Tween<sup>™</sup> 20, 10% BSA)
- 6 bottles (50 mL) Wash Buffer Concentrate 20x (PBS with 1% Tween<sup>™</sup> 20)
- 10 vials (15 mL) Substrate Solution (tetramethyl-benzidine)
- 1 vial (100 mL) Stop Solution (1M Phosphoric acid)

## 20 Adhesive Films

### Storage instructions – ELISA kit

Store kit reagents between 2°C and 8°C. Immediately after use remaining reagents should be returned to cold storage (2°C to 8°C). Expiry of the kit and reagents is stated on labels.

Expiry of the kit components can only be guaranteed if the components are stored properly, and if, in case of repeated use of one component, this reagent is not contaminated by the first handling.

### Sample collection and storage instructions

Cell culture supernatant and serum (baboon, cynomolgus) were tested with this assay. Other biological samples might be suitable for use in the assay. Remove serum or plasma from the clot or cells as soon as possible after clotting and separation.

Samples containing a visible precipitate must be clarified prior to use in the assay. Do not use grossly hemolyzed or lipemic specimens.

Samples should be aliquoted and must be stored frozen at -20°C to avoid loss of bioactive monkey TNF  $\alpha$ . If samples are to be run within 24 hours, they may be stored at 2°C to 8°C.

Avoid repeated freeze-thaw cycles. Prior to assay, the frozen sample should be brought to room temperature slowly and mixed gently. Do not thaw samples in a 37°C water bath. Do not vortex or sharply agitate samples.

### Materials required but not provided

- 5 mL and 10 mL graduated pipettes
- 5  $\mu$ L to 1000  $\mu$ L adjustable single channel micropipettes with disposable tips
- 50  $\mu$ L to 300  $\mu$ L adjustable multichannel micropipette with disposable tips
- Multichannel micropipette reservoir
- Beakers, flasks, cylinders necessary for preparation of reagents
- Device for delivery of wash solution (multichannel wash bottle or automatic wash system)
- Microwell strip reader capable of reading at 450 nm (620 nm as optional reference wave length)
- Glass-distilled or deionized water
- Statistical calculator with program to perform regression analysis

### Precautions for use

- All chemicals should be considered as potentially hazardous. We therefore recommend that this product is handled only by those persons who have been trained in laboratory techniques and that it is used in accordance with the principles of good laboratory practice. Wear suitable protective clothing such as laboratory overalls, safety glasses, and gloves. Care should be taken to avoid contact with skin or eyes. In the case of contact with skin or eyes wash immediately with water. See material safety data sheet(s) and/or safety statement(s) for specific advice.
- Reagents are intended for research use only and are not for use in diagnostic or therapeutic procedures.
- Do not mix or substitute reagents with those from other lots or other sources.
- Do not use kit reagents beyond expiration date on label.
- Do not expose kit reagents to strong light during storage or incubation.
- Do not pipet by mouth.
- Do not eat or smoke in areas where kit reagents or samples are handled.
- Avoid contact of skin or mucous membranes with kit reagents or samples.
- Rubber or disposable latex gloves should be worn while handling kit reagents or samples.
- Avoid contact of substrate solution with oxidizing agents and metal.
- Avoid splashing or generation of aerosols.

- In order to avoid microbial contamination or cross-contamination of reagents or specimens which may invalidate the test use disposable pipette tips and/or pipettes.
- Use clean, dedicated reagent trays for dispensing the conjugate and substrate reagent.
- Exposure to acid inactivates the conjugate.
- Glass-distilled water or deionized water must be used for reagent preparation.
- Substrate solution must be at room temperature prior to use.
- Decontaminate and dispose samples and all potentially contaminated materials as if they could contain infectious agents. The preferred method of decontamination is autoclaving for a minimum of 1 hour at 121.5°C.
- Liquid wastes not containing acid and neutralized waste may be mixed with sodium hypochlorite in volumes such that the final mixture contains 1.0% sodium hypochlorite. Allow 30 minutes for effective decontamination. Liquid waste containing acid must be neutralized prior to the addition of sodium hypochlorite.

## Preparation of reagents

1. Buffer concentrates should be brought to room temperature and should be diluted before starting the test procedure.
2. If crystals have formed in the Buffer Concentrates, warm them gently until they have completely dissolved.

## Wash buffer (1x)

1. Pour entire contents (50 mL) of the Wash Buffer Concentrate (20x) into a clean 1000 mL graduated cylinder. Bring to final volume of 1000 mL with glass-distilled or deionized water.
2. Mix gently to avoid foaming.
3. Transfer to a clean wash bottle and store at 2°C to 25°C. Please note that Wash Buffer (1x) is stable for 30 days.
4. Wash Buffer (1x) may also be prepared as needed according to the following table:

Number of Strips	Wash Buffer Concentrate (20x) (mL)	Distilled Water (mL)
1 - 6	25	475
1 - 12	50	950

## Assay buffer (1x)

1. Pour the entire contents (5 mL) of the Assay Buffer Concentrate (20x) into a clean 100 mL graduated cylinder. Bring to final volume of 100 mL with distilled water. Mix gently to avoid foaming.
2. Store at 2°C to 8°C. Please note that the Assay Buffer (1x) is stable for 30 days.
3. Assay Buffer (1x) may also be prepared as needed according to the following table:

Number of Strips	Assay Buffer Concentrate (20x) (mL)	Distilled Water (mL)
1 - 6	2.5	47.5
1 - 12	5.0	95.0

## Biotin-Conjugate

**Note:** The Biotin-Conjugate should be used within 30 minutes after dilution.

Make a 1:100 dilution of the concentrated Biotin-Conjugate solution with Assay Buffer (1x) in a clean plastic tube as needed according to the following table:

Number of Strips	Biotin-Conjugate (mL)	Assay Buffer (1x) (mL)
1 - 6	0.03	2.97
1 - 12	0.06	5.94

## Streptavidin-HRP

**Note:** The Streptavidin-HRP should be used within 30 minutes after dilution.

Make a 1:100 dilution of the concentrated Streptavidin-HRP solution with Assay Buffer (1x) in a clean plastic tube as needed according to the following table:

Number of Strips	Streptavidin-HRP (mL)	Assay Buffer (1x) (mL)
1 - 6	0.06	5.94
1 - 12	0.12	11.88


## Monkey TNF $\alpha$ total standard

1. Reconstitute monkey TNF  $\alpha$  total standard by addition of distilled water.
2. Reconstitution volume is stated on the label of the standard vial. Swirl or mix gently to insure complete and homogeneous solubilization (concentration of reconstituted standard = 3600 U/mL). Allow the standard to reconstitute for 10-30 minutes. Mix well prior to making dilutions.

The standard has to be used immediately after reconstitution and cannot be stored.

## External standard dilution

1. Label 7 tubes, one for each standard point: S1, S2, S3, S4, S5, S6, S7.
  2. Then prepare 1:2 serial dilutions for the standard curve as follows: Pipette 225  $\mu$ L of Sample Diluent into each tube.
  3. Pipette 225  $\mu$ L of reconstituted standard (concentration = 3600 U/mL) into the first tube, labeled S1, and mix (concentration of standard 1 = 1800 U/mL).
  4. Pipette 225  $\mu$ L of this dilution into the second tube, labeled S2, and mix thoroughly before the next transfer.
  5. Repeat serial dilutions 5 more times thus creating the points of the standard curve (see Figure 6).
- Sample Diluent serves as blank.

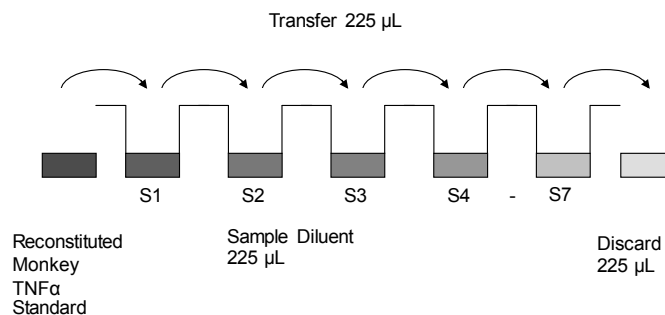

## Test protocol

1. Determine the number of microwell strips required to test the desired number of samples plus appropriate number of wells needed for running blanks and standards. Each sample, standard, blank and optional control sample should be assayed in duplicate. Remove extra microwell strips from holder and store in foil bag with the desiccant provided at 2°C to 8°C sealed tightly.
2. Wash the microwell strips twice with approximately 400  $\mu$ L Wash Buffer per well with thorough aspiration of microwell contents between washes. Allow the Wash Buffer to sit in the wells for about 10 – 15 seconds before aspiration. Take care not to scratch the surface of the microwells.

After the last wash step, empty wells and tap microwell strips on absorbent pad or paper towel to remove excess Wash Buffer. Use the microwell strips immediately after washing. Alternatively microwell strips can be placed upside down on a wet absorbent paper for not longer than 15 minutes. Do not allow wells to dry.

- Standard dilution on the microwell plate (Alternatively the standard dilution can be prepared in tubes - see "External standard dilution" on page 3):  
Add 100 µL of Sample Diluent in duplicate to all standard wells.  
Pipette 100 µL of prepared standard in duplicate into well A1 and A2 (see Table 1). Mix the contents of wells A1 and A2 by repeated aspiration and ejection (concentration of standard 1, S1 = 1800.0 U/mL), and transfer 100 µL to wells B1 and B2, respectively (see Figure 7). Take care not to scratch the inner surface of the microwells. Continue this procedure 5 times, creating two rows of monkey TNF α standard dilutions ranging from 1800.0 to 28.1 U/mL. Discard 100 µL of the contents from the last microwells (G1, G2) used.

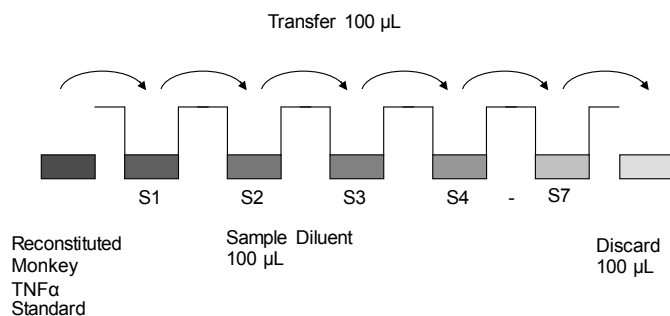


Fig. 7 Dilute standards - microwell plate

**Table 1** Example of the arrangement of blanks, standards and samples in the microwell strips.

	1	2	3	4
A	Standard 1 1800.0 U/mL	Standard 1 1800.0 U/mL	Sample 1	Sample 1
B	Standard 2 900.0 U/mL	Standard 2 900.0 U/mL	Sample 2	Sample 2
C	Standard 3 450.0 U/mL	Standard 3 450.0 U/mL	Sample 3	Sample 3
D	Standard 4 225.0 U/mL	Standard 4 225.0 U/mL	Sample 4	Sample 4
E	Standard 5 112.5 U/mL	Standard 5 112.5 U/mL	Sample 5	Sample 5
F	Standard 6 56.3 U/mL	Standard 6 56.3 U/mL	Sample 6	Sample 6
G	Standard 7 28.1 U/mL	Standard 7 28.1 U/mL	Sample 7	Sample 7
H	Blank	Blank	Sample 8	Sample 8

In case of an external standard dilution (see "External standard dilution" on page 3), pipette 100 µL of these standard dilutions (S1 - S7) in the standard wells according to Table 1.

- Add 100 µL of Sample Diluent in duplicate to the blank wells.
- Add 50 µL of Sample Diluent to the sample wells.
- Add 50 µL of each sample in duplicate to the sample wells.
- Prepare Biotin-Conjugate (see "Biotin-Conjugate" on page 3).
- Add 50 µL of Biotin-Conjugate to all wells.
- Cover with an adhesive film and incubate at room temperature (18°C to 25°C) for 2 hours on a microplate shaker.
- Prepare Streptavidin-HRP (refer to "Streptavidin-HRP" on page 3).
- Remove adhesive film and empty wells. Wash microwell strips 6 times according to point 2. of the test protocol. Proceed immediately to the next step.
- Add 100 µL of diluted Streptavidin-HRP to all wells, including the blank wells.
- Cover with an adhesive film and incubate at room temperature (18°C to 25°C) for 1 hour on a microplate shaker.
- Remove adhesive film and empty wells. Wash microwell strips 6 times according to point 2. of the test protocol. Proceed immediately to the next step.
- Pipette 100 µL of TMB Substrate Solution to all wells.

- Incubate the microwell strips at room temperature (18°C to 25°C) for about 10 minutes. Avoid direct exposure to intense light.

The color development on the plate should be monitored and the substrate reaction stopped before positive wells are no longer properly recordable. Determination of the ideal time period for color development has to be done individually for each assay.

It is recommended to add the stop solution when the highest standard has developed a dark blue color. Alternatively the color development can be monitored by the ELISA reader at 620 nm. The substrate reaction should be stopped as soon as Standard 1 has reached an OD of 0.9 – 0.95.

- Stop the enzyme reaction by quickly pipetting 100 µL of Stop Solution into each well. It is important that the Stop Solution is spread quickly and uniformly throughout the microwells to completely inactivate the enzyme. Results must be read immediately after the Stop Solution is added or within one hour if the microwell strips are stored at 2°C to 8°C in the dark.
- Read absorbance of each microwell on a spectro-photometer using 450 nm as the primary wave length (optionally 620 nm as the reference wave length; 610 nm to 650 nm is acceptable). Blank the plate reader according to the manufacturer's instructions by using the blank wells. Determine the absorbance of both the samples and the standards.

**Note:** In case of incubation without shaking the obtained O.D. values may be lower than indicated below. Nevertheless the results are still valid.

## Calculation of results

- Calculate the average absorbance values for each set of duplicate standards and samples. Duplicates should be within 20% of the mean value.
- Create a standard curve by plotting the mean absorbance for each standard concentration on the ordinate against the monkey TNF α concentration on the abscissa. Draw a best fit curve through the points of the graph (a 5-parameter curve fit is recommended).
- To determine the concentration of circulating monkey TNF α for each sample, first find the mean absorbance value on the ordinate and extend a horizontal line to the standard curve. At the point of intersection, extend a vertical line to the abscissa and read the corresponding monkey TNF α concentration.
- If instructions in this protocol have been followed, samples have been diluted 1:2 (50 µL sample + 50 µL Sample Diluent) and the concentration read from the standard curve must be multiplied by the dilution factor (x 2).
- Calculation of samples with a concentration exceeding standard 1 may result in incorrect, low monkey TNF α levels. Such samples require further external predilution according to expected monkey TNF α values with Sample Diluent in order to precisely quantitate the actual monkey TNF α level.
- It is suggested that each testing facility establishes a control sample of known monkey TNF α concentration and runs this additional control with each assay. If the values obtained are not within the expected range of the control, the assay results may be invalid.

- A representative standard curve is shown in Figure 8. This curve cannot be used to derive test results. Each laboratory must prepare a standard curve for each group of microwell strips assayed.

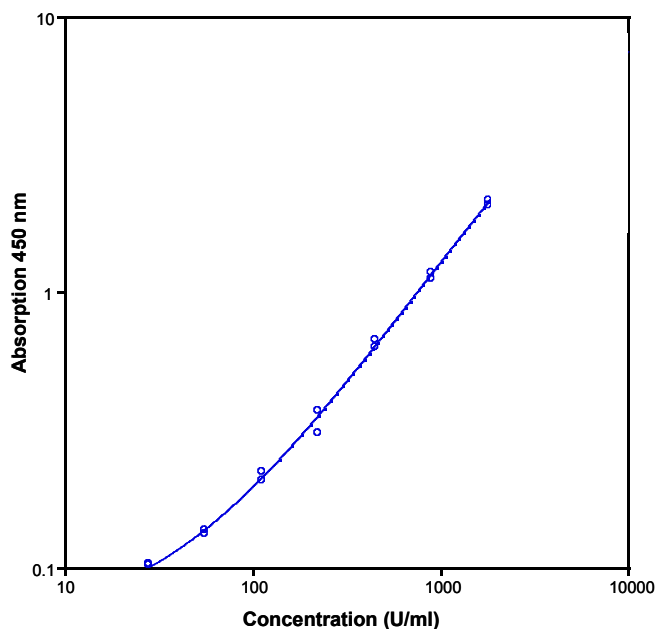


Fig. 8 Representative standard curve for monkey TNF  $\alpha$  total ELISA. Monkey TNF  $\alpha$  was diluted in serial 2-fold steps in Sample Diluent. Do not use this standard curve to derive test results. A standard curve must be run for each group of microwell strips assayed.

**Table 2** Typical data using the monkey TNF  $\alpha$  ELISA

Measuring wavelength: 450 nm

Reference wavelength: 620 nm

Standard	Monkey TNF $\alpha$ Concentration (U/mL)	O.D. at 450 nm	Mean O.D. at 450 nm	C.V. (%)
1	1800.0	2.055 2.145	2.100	2.1
2	900.0	1.180 1.120	1.150	2.6
3	450.0	0.667 0.631	0.649	2.8
4	225.0	0.371 0.306	0.339	9.6
5	112.5	0.223 0.206	0.215	4.0
6	56.3	0.137 0.133	0.135	1.5
7	28.1	0.102 0.103	0.103	0.5
Blank	0.0	0.044 0.046	0.045	2.2

The OD values of the standard curve may vary according to the conditions of assay performance (e.g. operator, pipetting technique, washing technique or temperature effects). Furthermore shelf life of the kit may affect enzymatic activity and thus color intensity. Values measured are still valid.

## Limitations

- Since exact conditions may vary from assay to assay, a standard curve must be established for every run.
- Bacterial or fungal contamination of either screen samples or reagents or cross-contamination between reagents may cause erroneous results.

- Disposable pipette tips, flasks or glassware are preferred, reusable glassware must be washed and thoroughly rinsed of all detergents before use.
- Improper or insufficient washing at any stage of the procedure will result in either false positive or false negative results. Empty wells completely before dispensing fresh wash solution, fill with Wash Buffer as indicated for each wash cycle and do not allow wells to sit uncovered or dry for extended periods.
- The use of mouse IgG antibodies or mouse IgG antibody coupled molecules in the course of animal experiments has significantly increased the number of animals with anti-mouse IgG antibodies (equivalent to human anti-mouse IgG antibodies - HAMA). These anti-mouse IgG antibodies may interfere with assays utilizing murine monoclonal antibodies leading to both false positive and false negative results. Serum samples containing antibodies to murine immunoglobulins can still be analyzed in such assays when murine immunoglobulins (serum, ascitic fluid, or monoclonal antibodies of irrelevant specificity) are added to the sample.

## Performance characteristics

### Sensitivity

The limit of detection of monkey TNF  $\alpha$  defined as the analyte concentration resulting in an absorbance significantly higher than that of the dilution medium (mean plus 2 standard deviations) was determined to be 6.0 U/mL (mean of 6 independent assays).

### Reproducibility

#### Intra-assay

Reproducibility within the assay was evaluated in 3 independent experiments. Each assay was carried out with 6 replicates of 8 samples containing different concentrations of monkey TNF  $\alpha$ . Two standard curves were run on each plate. Data below show the mean monkey TNF  $\alpha$  concentration and the coefficient of variation for each sample. The calculated overall intra-assay coefficient of variation was 8%.

#### Inter-assay

Assay to assay reproducibility within one laboratory was evaluated in 3 independent experiments. Each assay was carried out with 6 replicates of 8 samples containing different concentrations of monkey TNF  $\alpha$ . Two standard curves were run on each plate. Data below show the mean monkey TNF  $\alpha$  concentration and the coefficient of variation calculated on 18 determinations of each sample. The calculated overall inter-assay coefficient of variation was 8%.

### Spike recovery

The spike recovery was evaluated by spiking 4 levels of monkey TNF  $\alpha$  into samples. Recoveries were determined in 3 independent experiments with 6 replicates each.

Unspiked sample was used as blank in these experiments.

The overall mean recovery was 90%. Recoveries were shown to depend on the sample used.

### Dilution parallelism

4 samples with different levels of monkey TNF  $\alpha$  were analyzed at serial 2-fold dilutions with 4 replicates each.

The overall mean recovery was 113%. Recoveries were shown to depend on the sample used.

### Sample stability

#### Freeze-Thaw stability

Aliquots of samples (spiked or unspiked) were stored at -20°C and thawed 5 times, and the monkey TNF  $\alpha$  levels determined. A significant decrease of monkey TNF  $\alpha$  immunoreactivity was detected. Therefore samples should be stored in aliquots at -20°C and thawed only once.

## Storage stability

Aliquots of samples (spiked or unspiked) were stored at -20°C, 2°C to 8°C, room temperature, and at 37°C, and the monkey TNF  $\alpha$  level determined after 25 hours. There was no significant loss of monkey TNF  $\alpha$  immunoreactivity detected during storage at -20°C and 2°C to 8°C. A significant loss of monkey TNF  $\alpha$  immunoreactivity was detected during storage at room temperature and 37°C after 24 hours.

## Reagent preparation summary

### Wash buffer (1x)

Add Wash Buffer Concentrate 20x (50 mL) to 950 mL distilled water.

Number of Strips	Wash Buffer Concentrate (mL)	Distilled Water (mL)
1 - 6	25	475
1 - 12	50	950

### Assay buffer (1x)

Add Assay Buffer Concentrate 20x (5 mL) to 95 mL distilled water.

Number of Strips	Assay Buffer Concentrate (mL)	Distilled Water (mL)
1 - 6	2.5	47.5
1 - 12	5.0	95.0

### Biotin-Conjugate

Make a 1:100 dilution of Biotin-Conjugate in Assay Buffer (1x):

Number of Strips	Biotin-Conjugate (mL)	Assay Buffer (1x) (mL)
1 - 6	0.03	2.97
1 - 12	0.06	5.94

### Streptavidin-HRP

Make a 1:100 dilution of Streptavidin-HRP in Assay Buffer (1x):

Number of Strips	Streptavidin-HRP (mL)	Assay Buffer (1x) (mL)
1 - 6	0.06	5.94
1 - 12	0.12	11.88

### Monkey TNF $\alpha$ standard

Reconstitute lyophilized monkey TNF  $\alpha$  standard with distilled water. (Reconstitution volume is stated on the label of the standard vial.)

## Test protocol summary

1. Determine the number of microwell strips required.
2. Wash microwell strips twice with Wash Buffer.
3. Standard dilution on the microwell plate: Add 100  $\mu$ L Sample Diluent, in duplicate, to all standard wells. Pipette 100  $\mu$ L prepared standard into the first wells and create standard dilutions by transferring 100  $\mu$ L from well to well. Discard 100  $\mu$ L from the last wells.

Alternatively external standard dilution in tubes (see "External standard dilution" on page 3): Pipette 100  $\mu$ L of these standard dilutions in the microwell strips.

4. Add 100  $\mu$ L Sample Diluent, in duplicate, to the blank wells.
5. Add 50  $\mu$ L Sample Diluent to sample wells.
6. Add 50  $\mu$ L sample in duplicate, to designated sample wells.
7. Prepare Biotin-Conjugate.
8. Add 50  $\mu$ L Biotin-Conjugate to all wells.
9. Cover microwell strips and incubate 2 hours at room temperature (18°C to 25°C).
10. Prepare Streptavidin-HRP.
11. Empty and wash microwell strips 6 times with Wash Buffer.
12. Add 100  $\mu$ L diluted Streptavidin-HRP to all wells.
13. Cover microwell strips and incubate 1 hour at room temperature (18°C to 25°C).
14. Empty and wash microwell strips 6 times with Wash Buffer.
15. Add 100  $\mu$ L of TMB Substrate Solution to all wells.
16. Incubate the microwell strips for about 10 minutes at room temperature (18°C to 25°C).
17. Add 100  $\mu$ L Stop Solution to all wells.
18. Blank microwell reader and measure color intensity at 450 nm.

**Note:** If instructions in this protocol have been followed, samples have been diluted 1:2 (50  $\mu$ L sample + 50  $\mu$ L Sample Diluent) and the concentration read from the standard curve must be multiplied by the dilution factor ( $\times 2$ ).

## Customer and technical support

Visit [thermofisher.com/support](http://thermofisher.com/support) for the latest service and support information.

- Worldwide contact telephone numbers
- Product support information
  - Product FAQs
  - Software, patches, and updates
  - Training for many applications and instruments
- Order and web support
- Product documentation
  - User guides, manuals, and protocols
  - Certificates of Analysis
  - Safety Data Sheets (SDSs; also known as MSDSs)

**Note:** For SDSs for reagents and chemicals from other manufacturers, contact the manufacturer.

## Limited product warranty

Life Technologies Corporation and/or its affiliate(s) warrant their products as set forth in the Life Technologies' General Terms and Conditions of Sale at [www.thermofisher.com/us/en/home/global/terms-and-conditions.html](http://www.thermofisher.com/us/en/home/global/terms-and-conditions.html). If you have any questions, please contact Life Technologies at [www.thermofisher.com/support](http://www.thermofisher.com/support).



Bender MedSystems GmbH | Campus Vienna Biocenter 2 | 1030 Vienna, Austria  
For descriptions of symbols on product labels or product documents, go to [thermofisher.com/symbols-definition](http://thermofisher.com/symbols-definition).

The information in this guide is subject to change without notice.

**DISCLAIMER:** TO THE EXTENT ALLOWED BY LAW, THERMO FISHER SCIENTIFIC INC. AND/OR ITS AFFILIATE(S) WILL NOT BE LIABLE FOR SPECIAL, INCIDENTAL, INDIRECT, PUNITIVE, MULTIPLE, OR CONSEQUENTIAL DAMAGES IN CONNECTION WITH OR ARISING FROM THIS DOCUMENT, INCLUDING YOUR USE OF IT.

**Important Licensing Information:** These products may be covered by one or more Limited Use Label Licenses. By use of these products, you accept the terms and conditions of all applicable Limited Use Label Licenses.

©2019 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific and its subsidiaries unless otherwise specified. All other trademarks are the property of their respective owners.

[thermofisher.com/support](http://thermofisher.com/support) | [thermofisher.com/askaquestion](http://thermofisher.com/askaquestion)

[thermofisher.com](http://thermofisher.com)