


Catalog # CSC1753

Lot #: 102001
Expiration: 31.07.2011

- 1 **Coating Antibody:** **Anti-Swine TNF- α (0.125mg / 0.125 ml)**
 Part Number: 5S.175.09
 Lot Number: **10E8/1**
 Form: Liquid, 1 vial, contains 0.1% sodium azide
 Storage: Store at 2-8°C until expiration date.
 Recommended Dilution: Dilute to 1 μ g/mL with Coating Buffer B (Cat. # CB01100, or see Recommended Buffers). For example, to make 10 mL (enough to coat 1 plate), add 10 μ L coating antibody to 9.990 mL Coating Buffer B.
- 2 **Detection Antibody:** **Anti-Swine TNF- α Biotin (25 μ g/0.125 mL)**
 Part Number: 5S.175.03
 Lot Number: **10E8/1**
 Form: Liquid, 1 vial, contains 0.1% sodium azide
 Storage: Store at 2-8°C until expiration date.
 Recommended Dilution: Dilute to 0.4 μ g/mL with Assay Buffer (Cat. # DS98200, or see Recommended Buffers). For example, to make enough for 1 plate, add 11 μ L detection antibody to 5.489 mL Assay Buffer.
- 3 **Standard:** **Recombinant Swine TNF- α**
 Part Number: 5S.175.10 (additional vials of standard may be purchased using this part number)
 Lot Number: **8G3/1**
 Form: Lyophilized, 3 vials
 Storage: Store at 2-8°C.
 Reconstitution: Reconstitute with Assay Buffer (Cat. # DS98200 or see Recommended Buffers) to yield a stock of 20,000 pg/mL. After 10 minutes of rehydration, use the standard stock immediately or aliquot in polypropylene tubes and freeze at -80°C. Do not store at room temperature or at 4°C and do not subject to more than one freeze-thaw cycle.
 Standard Curve: Dilute standard stock to 4,000 pg/mL (120 μ L stock plus 480 μ L Assay Buffer) with Assay Buffer (Cat. # DS98200 or see Recommended Buffers). Add 300 μ L Assay Buffer to 6 tubes and label as 2,000, 1,000, 500, 250, 125 and 62.5 pg/mL. Make serial dilutions starting with 4,000 pg/mL by transferring 300 μ L of each standard to next tube and vortexing each tube. Assay Buffer should be used as the zero standard.
- 4 **Streptavidin-HRP:** **0.25 mL**
 Part Number: 4100024
 Lot Number: **10D6/1**
 Form: Liquid, 1vial, contains 0.05% thymol
 Storage: Store at 2-8°C until expiration date.
 Recommended Dilution: Dilute to 1/7500. For example, to make enough for 1 plate, add 1.3 μ L of streptavidin-HRP to 9.9987 mL of Assay Buffer (Cat. # DS98200 or see Recommended Buffers).

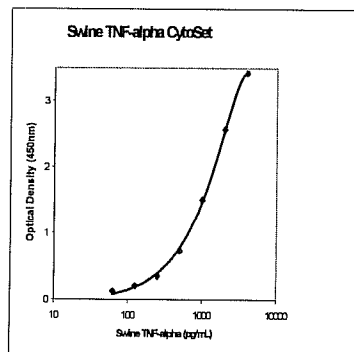

Representative standard curve was generated by following the recommended assay procedure, which includes the use of the **BioSource CytoSet™ Buffer Set (Cat. # CNB0011)**

Intended Use and Materials Provided

The CytoSet™ for Swine TNF- α contains components required to construct an enzyme-linked immunoassay for the specific and quantitative measurement of TNF- α . Sufficient quantities of all reagents are provided to yield 10 plates of 96 wells if the recommended assay procedure and recommended storage and handling of materials are followed as specified on this insert. The materials provided are **FOR RESEARCH USE ONLY**.

Recommended Buffers and Solutions

The BioSource CytoSet™ Buffer Set (Cat. # CNB0011) containing Coating Buffers A and B, Assay Buffer, Substrate Solution (TMB), Stop Solution, and Wash Buffer is recommended.

1. **Coating Buffer A:** *Coating Buffer A (Cat. # CB07100) from BioSource is recommended. Alternate buffer choice listed below.*
8.0 g NaCl, 1.13 g Na₂HPO₄, 0.2 g KH₂PO₄, 0.2 g KCl, 0.1% ProClin™; q.s. to 1.0 L with distilled H₂O, pH to 7.4.
2. **Coating Buffer B:** *Coating Buffer B (Cat. # CB01100) from BioSource is recommended. Alternate buffer choice listed below.*
4.3 g NaHCO₃, 5.3 g Na₂CO₃, 0.1% ProClin™; q.s. to 1.0 L with distilled H₂O, pH to 9.4.
3. **Assay Buffer:** *Assay Buffer (Cat. # DS98200) from BioSource is recommended. Alternate buffer choice listed below.*
8.0 g NaCl, 1.13 g Na₂HPO₄, 0.2 g KH₂PO₄, 0.2 g KCl, 5.0 g bovine serum albumin (fraction V), 1 mL Tween 20 and 0.5% ProClin™ as a preservative; q.s. to 1.0 L with distilled H₂O, pH to 7.4.
4. **Wash Buffer:** *Wash Buffer 25x (Cat. # WB01) from BioSource is recommended. Alternate buffer choice listed below.*
0.2 g KH₂PO₄, 1.9 g, K₂HPO₄ · 3H₂O, 0.4 g EDTA, 0.5 mL Tween 20; q.s. to 1.0 L with distilled H₂O, pH to 7.4.
5. **Substrate Solution:** *TMB (Cat. # SB01) from BioSource is recommended. Alternate solution choice listed below.*
Tetramethylbenzidine (TMB) and Hydrogen Peroxide.
6. **Stop Solution:** *Stop Solution (Cat.# SS01100) from BioSource is recommended. Alternate solution choice listed below.*
1.8 N H₂SO₄.

Assay Optimization

CytoSets™ from BioSource are designed to be very flexible for your experiments. Consequently, the assay procedure contains only recommendations. The assay procedure has been optimized for use with tissue culture samples. However, serum and plasma samples may be used but may require that certain assay parameters be modified. Investigators are advised to determine optimal buffer formulations, concentrations and incubation times for individual applications.

Recommended Assay Procedure

1. Prepare coating solution by diluting the coating antibody. See “coating antibody” section for the recommended coating antibody dilution.
2. Coat plates with 100 μ L per well of the coating solution. Cover plates and incubate overnight (12-18 hr.) at 4°C.
3. Aspirate wells and wash 1 time with > 400 μ L of Wash Buffer per well. Following wash, invert and tap on absorbent paper to remove excess liquid.
4. Block plate with 300 μ L per well of Assay Buffer for 1 hour at room temperature.
5. Aspirate, invert, and tap on absorbent paper to remove excess liquid.
6. Prepare standards and sample dilutions in Assay Buffer (or in a diluent that most closely matches the matrix of your sample). For recommended dilutions and storage of the standard, see “standard” section.
7. Pipette 100 μ L of standards (in duplicate) and samples into designated wells.
8. Immediately following step 7, add 50 μ L of the working detection antibody into each well. For recommended dilutions, see “detection antibody” section. **Incubate for 2 hours at room temperature.**
9. Aspirate and wash 5 times using the method in step 3.
10. Add 100 μ L of the working streptavidin-HRP solution into each well. For recommended dilutions, see “streptavidin-HRP conjugate” section. **Incubate for 30 minutes at room temperature.**
11. Aspirate and wash 5 times using the method in step 3.
12. Add 100 μ L of the TMB substrate to each well. Incubate plate for 30 minutes at room temperature.
13. Add 100 μ L of Stop Solution to each well.
14. Measure absorbance at 450 nm (reference absorbance: 650 nm) within 30 minutes of adding Stop Solution. Calculate results using a log-log or 4-parameter curve fit.

Additional Materials Required

- 96 well NUNC MaxiSorp microplates; NUNC Cat. # 434797.
- Pipettes, shaker and timer.
- Microplate reader with a detector that can measure absorbance at 450 nm.
- 1 L graduated cylinder; plate washer or wash bottle.
- Polypropylene tubes for standards and sample dilutions, if needed.

Invitrogen Corporation., USA • 542 Flynn Road, Camarillo, CA 93012 • (800) 242-0607 • FAX (805) 987-3385
Diasource ImmunoAssays, S.A. • Rue de l'Industrie 8, B-1400 Nivelles, Belgium • +32 67 88 99 99 • FAX +32 67 88 99 96
This product is for research use only. Not for use in diagnostic procedures.