

ProcartaPlex™ Human Neuroscience Panel 18-Plex USER GUIDE

Catalog Number EPX180-15837-901

Publication Number MAN0019827

Revision A.0 (30)



Bender MedSystems GmbH | Campus Vienna Biocenter 2 | 1030 Vienna, Austria

For descriptions of symbols on product labels or product documents, go to [thermofisher.com/symbols-definition](https://www.thermofisher.com/symbols-definition).

The information in this guide is subject to change without notice.

DISCLAIMER: TO THE EXTENT ALLOWED BY LAW, THERMO FISHER SCIENTIFIC INC. AND/OR ITS AFFILIATE(S) WILL NOT BE LIABLE FOR SPECIAL, INCIDENTAL, INDIRECT, PUNITIVE, MULTIPLE, OR CONSEQUENTIAL DAMAGES IN CONNECTION WITH OR ARISING FROM THIS DOCUMENT, INCLUDING YOUR USE OF IT.

Revision history: Pub. No. MAN0019827

Revision	Date	Description
A.0 (30)	26 January 2021	New manual

End User License Agreement for Luminex Assay Products: By opening the packaging containing this Assay Product (which contains fluorescently labeled microsphere beads authorized by Luminex Corporation) or using this Assay Product in any manner, you are consenting and agreeing to be bound by the End User Terms and Conditions and the End User License Agreement available at <https://www.luminexcorp.com/end-user-terms-and-conditions/#sec3>. If you do not agree to all of the terms and conditions, you must promptly return this Assay Product for a full refund prior to using it in any manner.

TRADEMARKS: xMAP, Luminex 100/200, MAGPIX, INTELLIFLEX and FLEXMAP 3D are trademarks of Luminex Corp. and Bio-Plex is a trademark of Bio-Rad Laboratories Inc. All trademarks are the property of Thermo Fisher Scientific and its subsidiaries unless otherwise specified.

©2021 Thermo Fisher Scientific Inc. All rights reserved.

Contents

■	CHAPTER 1	ProcartaPlex™ Human Neuroscience Panel 18-Plex	4
		Introduction	4
		Contents and storage	5
		Required materials not supplied	5
		Precautions and technical hints	6
		Workflow	7
■	CHAPTER 2	Methods	8
		Sample preparation	8
		Plasma sample preparation	8
		Serum sample preparation	8
		CSF sample preparation	8
		Preparation of reagents	8
		Prepare 1X Wash Buffer	9
		Prepare Standard Mixes	9
		Prepare 4-fold serial dilution	9
		Assay protocol	11
		Instrument settings	13
		Analyze results	13
■	APPENDIX A	Recommended plate layout	14



ProcartaPlex™ Human Neuroscience Panel 18-Plex



WARNING! Read the Safety Data Sheets (SDSs) and follow the handling instructions. Wear appropriate protective eyewear, clothing, and gloves. Safety Data Sheets (SDSs) are available from [thermofisher.com/support](https://www.thermofisher.com/support).

Introduction

The ProcartaPlex™ Human Neuroscience Panel 18-Plex has been optimized for detection of multiple analytes from serum, plasma, and cerebrospinal fluid (CSF).

The panel is provided in a ready-to-use format with individual vials of 1X capture and detection reagents that require less pipetting and experimental setup. These reagents are not combinable with simplexes or other panels.

ProcartaPlex™ preconfigured panels are extensively tested for analyte combinability, interference and cross-reactivity to provide the highest level of validation and precision. All ProcartaPlex™ panels are supplied with the necessary reagents to perform the assay.

Analytes		
Amyloid beta 1-42	BDNF	CNTF
FGF-21	GDNF	GFAP
Kallikrein-6 (KLK6)	MIF	NCAM-1
Neurogranin (NRGN)	NF-H	NGF beta
S100B	Tau (total)	Tau (pT181)
TDP-43	UCHL1	YKL-40 (CHI3L1)

For detailed product information, visit [thermofisher.com/procartaplex](https://www.thermofisher.com/procartaplex)

Contents and storage

Upon receipt, store the kit at 2°C to 8°C. When stored as indicated, all reagents are stable until the expiration date.

Contents	Amount
Standard Mix 2 (lyophilized)	2 each
Standard Mix 4 (lyophilized)	2 each
Standard Mix 5 (lyophilized)	2 each
Biotinylated Detection Antibody Mix (1X)	1 x 3.5 mL
Capture Bead Mix (1X)	1 x 5 mL
Streptavidin-PE (SA-PE) (1X)	1 x 5 mL
Wash Buffer (10X)	1 x 25 mL
Reading Buffer (1X)	1 x 40 mL
Universal Assay Buffer (1X)	1 x 10 mL
8-Tube Strip	2 each
Flat Bottom 96-well Plate, black	1 each
Microplate Lid	1 each
Plate Seals	8 each

Retain the lot-specific Certificate of Analysis that contains the product expiration date. The Certificate of Analysis also contains important information such as bead number, analyte names and highest standard concentration required for the assay setup on the xMAP instrument.



CAUTION! This kit contains materials with small quantities of sodium azide. Sodium azide reacts with lead and copper plumbing to form explosive metal azides. Upon disposal, flush drains with a large volume of water to prevent azide accumulation. Avoid ingestion and contact with eyes, skin and mucous membranes. In case of contact, rinse affected area with plenty of water. Observe all federal, state and local regulations for disposal.

Required materials not supplied

- xMAP™ instrument
- Hand-Held Magnetic Plate Washer (Cat. No. EPX-55555-000)
- Deionized water
- Vortex mixer (e.g., Cat. No. 88882010)
- Microcentrifuge


- Adjustable single and multichannel pipettes with disposable tips and low volume reservoirs (e.g., Cat. No. 95128093)
- Beakers, flasks, and cylinders necessary for preparation of reagents
- Orbital microplate shaker with at least 1.5 mm or 0.059 inch orbit diameter capable of maintaining a speed of 600 ± 50 rpm (e.g., Cat. No. 88882006)

Note: The use of rockers or large orbit shakers may cause adverse results.

Precautions and technical hints

1. Thoroughly read this User Guide and Certificate of Analysis prior to using the kit.
2. All chemicals should be considered potentially hazardous.
3. To avoid cross-contamination, do not invert the assay plate during the assay or allow contents from one well to mix with another well.
4. Use a multichannel pipette and reagent reservoirs whenever possible to achieve optimal assay precision.
5. Ensure that the xMAP™ instrument has been properly calibrated and set up prior to preparing and running the assay.

Workflow

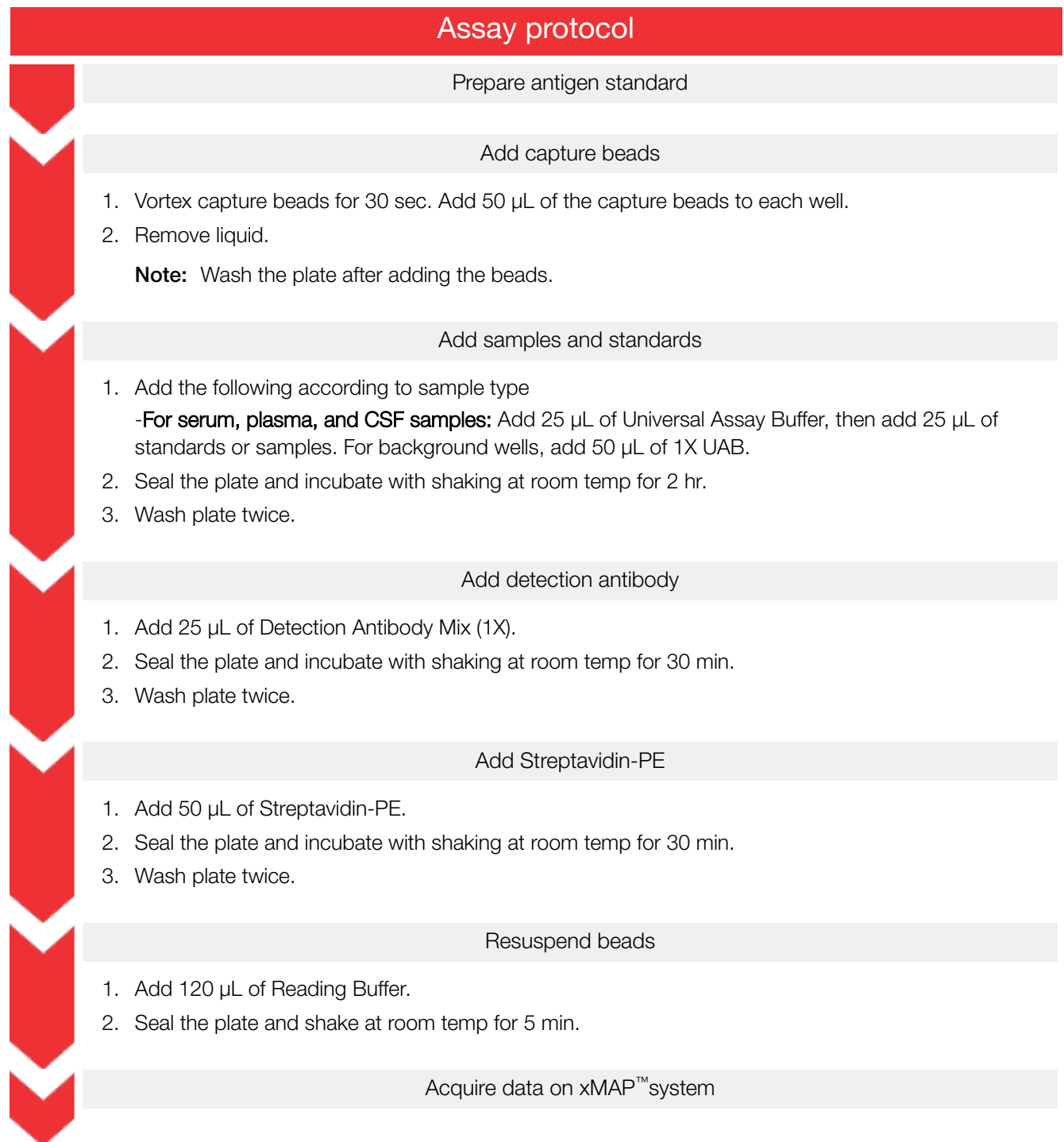

Sample preparation

Thaw frozen serum and plasma samples on ice and mix well by vortexing. Centrifuge at $10,000 \times g$ for 5–10 minutes to pellet out particulates. Avoid multiple freeze/thaw cycles. If samples are high in lipid content, centrifuge at $10,000 \times g$ for 10 minutes and transfer contents to a new tube.

Plasma sample preparation

1. Collect samples in sodium citrate or EDTA tubes. If using heparin as an anticoagulant, no more than 10 IU of heparin per mL of blood collected should be used to prevent assay interference that can result in a false positive signal.
2. Centrifuge samples at $1,000 \times g$ at 4°C for 10 minutes within 30 minutes of collection.
3. Collect the plasma fraction. Use immediately or store aliquots at -80°C .

Serum sample preparation

1. Allow blood to clot for 20–30 minutes at $20\text{--}25^{\circ}\text{C}$.
2. Centrifuge at $1,000 \times g$ for 10 minutes at $20\text{--}25^{\circ}\text{C}$.
3. Collect the serum fraction. Alternatively, a serum separator tube can be used following the manufacturer's instructions.
4. Use immediately or store aliquots at -80°C . Avoid multiple freeze/thaw cycles.

CSF sample preparation

1. Centrifuge samples at 1,400 rpm for 10 minutes at 4°C to remove particulates.
2. Use immediately or store aliquots at -80°C . Avoid multiple freeze/thaw cycles.

Preparation of reagents

Before starting with the assay protocol, define the plate map. Mark the standard, sample and background wells using the plate map found in Appendix A, "Recommended plate layout" to determine the number of wells used.

Prepare 1X Wash Buffer

Bring the Wash Buffer Concentrate (10X) to room temperature and vortex for 15 seconds. Mix 20 mL of the Wash Buffer Concentrate (10X) with 180 mL ddH₂O. Mix gently to avoid foaming. Wash Buffer (1X) can be stored at 2–8°C for up to 6 months.

Note: Additional Wash Buffer Concentrate (200 mL, Cat. No. EPX-66666-001) can be purchased separately for automated plate washers.

Prepare Standard Mixes

This kit is supplied with three lyophilized Standard Mixes for generation of standard curves. Two vials of each Standard Mix are provided to permit the user to run the assay twice if running a partial plate. For experiments measuring serum, plasma, or CSF samples, use 1X UAB as the diluent to reconstitute and dilute the standard.

Note: Change pipette tips after each dilution step and avoid air bubbles.


1. Centrifuge the three standard mix stock vials at 2,000 × *g* for 10 seconds.
2. Add 50 µL of diluent to each stock vial.
3. Vortex the vials at high speed for 30 seconds and centrifuge at 2,000 × *g* for 10 seconds to collect contents at the bottom of the vials.
4. Incubate on ice for 10 minutes to ensure complete reconstitution.
5. Combine the contents of each vial into one of the vials. Add 100 µL of diluent to the combined contents for a total volume of 250 µL.
6. Vortex the vial at high speed for 30 seconds and centrifuge at 2,000 × *g* for 10 seconds to collect contents at the bottom of the vial.

Prepare 4-fold serial dilution

1. Label the tubes in the 8-Tube Strip: Std1, Std2, Std3, Std4, Std5, Std6 and Std7.
2. Add 200 µL of the reconstituted standard mix into Std1 tube.
3. Add 150 µL of diluent into Std2–Std7 tubes.
4. Transfer 50 µL from Std1 tube into Std2 tube.
5. Mix by pipetting up and down 10 times.
6. Transfer 50 µL of the mixed standards from Std2 tube into Std3 tube using new pipette tip.
7. Mix by pipetting up and down 10 times.
8. Repeat steps 4–7 for tubes Std4–Std7, changing pipette tips between dilution steps.
9. Add 150 µL of diluent to the last tube of the 8-Tube Strip to serve as a background.

10. Keep tubes on ice until ready to use.

Note: Use the reconstituted standard immediately. The reconstituted standard cannot be stored. Discard unopened standard vials if the entire plate was used in a single experiment.

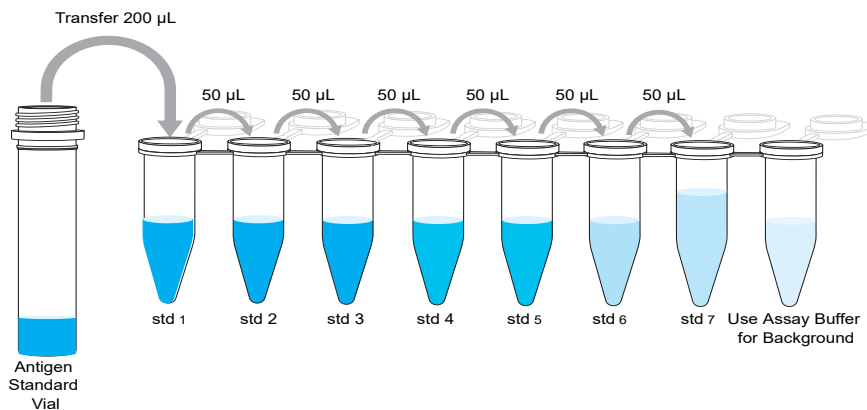

Assay protocol

1. Add Capture Bead Mix to the plate.
 - a. Vortex the 1X Capture Bead Mix vial for 30 seconds at high speed.
 - b. Using a multichannel pipette, add 50 μ L of the Capture Bead Mix to each well of the plate.
2. Wash beads using a Hand-Held Magnetic Plate Washer.

Note: To avoid loss of beads, secure the plate using the clamps on both sides of the Hand-Held Magnetic Plate Washer during this procedure.

Note: This protocol was developed using the Hand-Held Magnetic Plate Washer (Cat. No. EPX-55555-000). Other washers should be validated by the end user.

- a. Place the plate on the Hand-Held Magnetic Plate Washer and wait 2 minutes to allow the beads to settle on the bottom of each well.
 - b. Remove the liquid by quickly inverting the washer/plate assembly over a sink or waste container.
 - c. Gently blot the inverted washer/plate assembly onto several layers of paper towels or absorbent surface to remove any residual liquid.
 - d. Add 150 μ L of 1X Wash Buffer into each well and wait 30 seconds.
 - e. Remove the liquid by quickly inverting the washer/plate assembly over a sink or waste container.
 - f. Gently blot the inverted washer/plate assembly onto several layers of paper towels or absorbent surface to remove any residual liquid.
 - g. Remove the plate from the magnet and proceed to step 3.
3. Add samples and standards to the plate.
 - a. **Serum, plasma, and CSF:** Add 25 μ L of 1X UAB to each well followed by 25 μ L of prepared standards or samples as defined on the plate layout. Add an additional 25 μ L of 1X UAB to the wells designated as backgrounds.
 - b. Seal the plate using one of the provided Plate Seals and cover with the provided Microplate Lid. Shake at 600 rpm for 2 hours at room temperature.

Note: For those wishing to perform the assay over two days, the 96-well plate can be incubated overnight. Shake the 96-well plate for 30 minutes at room temperature at 600 rpm, then transfer the plate to 4°C and store on a level surface. After overnight incubation, shake the plate for an additional 30 minutes at room temperature at 600 rpm.

4. Remove and discard the Plate Seal. Wash the plate following the steps below.
 - a. Place the plate on the Hand-Held Magnetic Plate Washer and wait 2 minutes to allow particles to settle on the bottom of each well.

- b. Remove the liquid by quickly inverting the washer/plate assembly over a sink or waste container.
 - c. Gently blot the inverted washer/plate assembly onto several layers of paper towels or absorbent surface to remove any residual liquid.
 - d. Add 150 μL of 1X Wash Buffer into each well and wait 30 seconds.
 - e. Remove the liquid by quickly inverting the washer/plate assembly over a sink or waste container.
 - f. Gently blot the inverted washer/plate assembly onto several layers of paper towels or absorbent surface to remove any residual liquid.
 - g. Repeat steps 4d-4f once for a total of two washes.
 - h. Remove the plate from the magnet and proceed to the next step.
5. Add Biotinylated detection Antibody Mix to the plate.
 - a. Using a multichannel pipette, add 25 μL of the detection antibody solution to each well of the plate. Gently tap the plate to evenly distribute the solution in the wells.

Note: A narrow trough reservoir for multichannel pipetting is recommended to be used to prevent volume loss.

- b. Seal the plate using a new Plate Seal and cover with the provided Microplate Lid. Shake at 600 rpm for 30 minutes at room temperature.
6. Wash the plate following step 4.
7. Add Streptavidin-PE (SA-PE) to the plate.
 - a. Add 50 μL of SA-PE solution to each well.
 - b. Seal the plate using new Plate Seal and cover with the provided Microplate Lid. Shake at 600 rpm for 30 minutes at room temperature.
8. Wash the plate following step 4.
9. Prepare the plate for analysis on a xMAP™ instrument.
 - a. Add 120 μL of reading buffer into each well.
 - b. Seal the plate using new Plate Seal and cover with the provided Microplate Lid. Shake at 600 rpm for 5 minutes at room temperature.
10. Remove the Plate Seal and run the plate on a xMAP™ instrument.

Instrument settings

Follow the recommended guidelines and procedures for calibration and verification of the instrument. Laser-based systems require 30 minutes to warm up prior to use.

Instrument	Acquisition volume	Timeout (optional)	Bead type	DD gate	Reporter gain	Min. bead count
MAGPIX™	50 µL ^[1]	N/A	N/A	N/A	Standard PMT	50
INTELLIFLEX™	30 µL	40 sec	MagPlex™	4,000–13,000	Standard PMT	50
FLEXMAP 3D™ Luminex™ 100/200™	50 µL	60 sec	MagPlex™	7,500–25,000	Standard PMT	50
Bio-Rad™ Bio-Plex™	50 µL	60 sec	MagPlex™	5,000–25,000	Standard PMT	50

^[1] MAGPIX volume can be changed during the run to optimize bead count.

Note: To assure a good bead count, the probe height must be adjusted to the plate provided in the kit. We recommend using two 5.08 mm spacer disks to adjust the sample probe height for Mylar-bottom plates.

Analyze results

The concentration of the samples can be calculated by plotting the expected concentration of the standards against the NET MFI generated by each standard. For Bio-Plex™ Manager, plot standard concentrations against FI-Bkgd. A 4PL or 5PL algorithm is recommended for the best curve fit. Analyze the assayed samples according to the operation manual for the Luminex™ or Bio-Plex™ instrument.

We offer a free and robust analysis software package for data analysis. To analyze the data, follow the instructions below or contact our technical support.

1. Export the run data in .csv format and navigate to the ProcartaPlex™ Analysis App on Thermo Fisher Connect: <https://apps.thermofisher.com/apps/procartaplex>

Note: Before exporting .csv raw data from Bio-Plex™ Manager, please make sure to set 'Analytes Labels' under 'Document Export Properties' to 'Name (Region)'. The .csv raw data exported as Report Type 'xPONENT' from INTELLIFLEX™ instruments are supported.

2. Upload the .csv files to the ProcartaPlex™ Analysis App to analyze the run data. The intuitive software features 4PL/5PL curve fit optimization, group-wise statistical and heat map analysis. Users can export detailed reports including images for presentations and publications.

IMPORTANT! For ProcartaPlex™ getting started guides, technical literature, protocol support tools, and common troubleshooting questions visit thermofisher.com/procartaplexsupport. For more complete troubleshooting questions and answers, visit our FAQ database at thermofisher.com/procartaplexfaqs



Recommended plate layout

Standards		Samples									
1	1	1	1	9	9	17	17	25	25	33	33
2	2	2	2	10	10	18	18	26	26	34	34
3	3	3	3	11	11	19	19	27	27	35	35
4	4	4	4	12	12	20	20	28	28	36	36
5	5	5	5	13	13	21	21	29	29	37	37
6	6	6	6	14	14	22	22	30	30	38	38
7	7	7	7	15	15	23	23	31	31	39	39
Bkgd ^[1]	Bkgd	8	8	16	16	24	24	32	32	40	40

^[1] Background

	1	2	3	4	5	6	7	8	9	10	11	12
A												
B												
C												
D												
E												
F												
G												
H												

