

MHC Class II (I-A/I-E) Monoclonal Antibody (M5/114.15.2), APC, eBioscience™

Product Details	
Size	100 µg
Species Reactivity	Mouse
Published Species	Mouse, Human
Host/Isotype	Rat / IgG2b, kappa
Recommended Isotype Control	Rat IgG2b kappa Isotype Control (eB149/10H5), APC, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	M5/114.15.2
Conjugate	APC
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.1% gelatin
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_469455

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	-	1 Publication
Immunocytochemistry (ICC/IF)	-	1 Publication
Flow Cytometry (Flow)	0.03 µg/test	109 Publications

Product Specific Information

Description: The M5/114.15.2 monoclonal antibody reacts with the mouse major histocompatibility complex class II, both I-A and I-E subregion-encoded glycoproteins (I-A b, I-A d, I-A q, I-E d, I-E k, not I-A f, I-A k, or I-A s). It detects a polymorphic determinant present on B cells, monocytes, macrophages, dendritic cells, and activated T lymphocytes from mice carrying the H-2 b, H-2 d, H-2 q, H-2 p, H-2 r and H-2 u but not from mice carrying the H-2 s or H-2 f haplotypes. The M5/114 mAb is reported to inhibit I-A-restricted T cell responses of the H-2 b, H-2 d, H-2 q, H-2 u but not H-2 f, H-2 k, or H-2 s haplotypes.


Applications Reported: M5/114.15.2 has been reported for use in flow cytometric analysis.

Applications Tested: The M5/114.15.2 antibody has been tested by flow cytometric analysis of mouse splenocytes and can be used at less than or equal to 0.03 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Excitation: 633-647 nm; **Emission:** 660 nm; **Laser:** Red Laser.

Filtration: 0.2 µm post-manufacturing filtered.

Product Images For MHC Class II (I-A/I-E) Monoclonal Antibody (M5/114.15.2), APC, eBioscience™

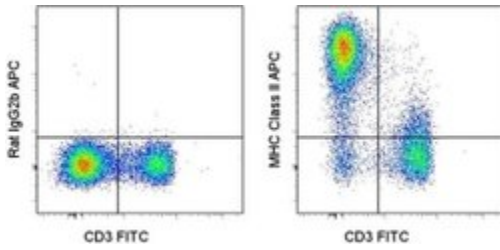

MHC Class II (I-A/I-E) Antibody (17-5321-82) in Flow

Staining of C57BL/6 splenocytes with Anti-Mouse CD3e FITC (Product # 11-0031-82) and 0.015 µg of Rat IgG2b K Isotype Control APC (Product # 17-4031-82) (left) or 0.015 µg of Anti-Mouse MHC Class II (I-A/I-E) APC (right). Cells in the lymphocyte gate were used for analysis.

[View more figures on thermofisher.com](http://thermofisher.com)

Immunohistochemistry (1)

Ophthalmology and eye diseases

Stratification of Antigen-presenting Cells within the Normal Cornea.

"17-5321 was used in Immunofluorescence on wholemounted tissues to investigate the presence of antigen-presenting cells within the non-inflamed cornea."

Authors: Knickelbein JE,Watkins SC,McMenamin PG,Hendricks RL

Species
Mouse

Dilution
Not Cited

Year
2009

Immunocytochemistry (1)

Nature

Sessile alveolar macrophages communicate with alveolar epithelium to modulate immunity.

"17-5321 was used in Immunofluorescence to show that a subset of alveolar macrophages (AMs) attached to the alveolar wall form connexin 43 (Cx43)-containing gap junction channels with the epithelium."

Authors: Westphalen K,Gusarova GA,Islam MN,Subramanian M,Cohen TS,Prince AS,Bhattacharya J

Species
Mouse

Dilution
Not Cited

Year
2014

Flow Cytometry (109)

International journal of biological sciences

MiR-103 protects from recurrent spontaneous abortion via inhibiting STAT1 mediated M1 macrophage polarization.

"Published figure using MHC Class II (I-A/I-E) monoclonal antibody (Product # 17-5321-82) in Flow Cytometry"

Authors: Zhu X,Liu H,Zhang Z,Wei R,Zhou X,Wang Z,Zhao L,Guo Q,Zhang Y,Chu C,Wang L,Li X

Species
Not Applicable

Dilution
Not Cited

Year
2021

Frontiers in immunology

Dietary Glucose Consumption Promotes RALDH Activity in Small Intestinal CD103⁺CD11b⁺ Dendritic Cells.

"Published figure using MHC Class II (I-A/I-E) monoclonal antibody (Product # 17-5321-82) in Flow Cytometry"

Authors: Ko HJ,Hong SW,Verma R,Jung J,Lee M,Kim N,Kim D,Surh CD,Kim KS,Rudra D,Im SH

Species
Not Applicable

Dilution
Not Cited

Year
2021

[View more Flow references on thermofisher.com](#)

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.