

# MHC Class II (I-A/I-E) Monoclonal Antibody (M5/114.15.2), APC-eFluor 780, eBioscience™

Product Details	
Size	100 µg
Species Reactivity	Mouse
Published Species	Mouse
Host/Isotype	Rat / IgG2b, kappa
Recommended Isotype Control	Rat IgG2b kappa Isotype Control (eB149/10H5), APC-eFluor 780, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	M5/114.15.2
Conjugate	APC-eFluor® 780
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_1548783

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	0.5 µg/test	73 Publications

## Product Specific Information

**Description:** The M5/114.15.2 monoclonal antibody reacts with the mouse major histocompatibility complex class II, both I-A and I-E subregion-encoded glycoproteins (I-A b, I-A d, I-A q, I-E d, I-E k, not I-A f, I-A k, or I-A s). It detects a polymorphic determinant present on B cells, monocytes, macrophages, dendritic cells, and activated T lymphocytes from mice carrying the H-2 b, H-2 d, H-2 q, H-2 p, H-2 r and H-2 u but not from mice carrying the H-2 s or H-2 f haplotypes. The M5/114 mAb is reported to inhibit I-A-restricted T cell responses of the H-2 b, H-2 d, H-2 q, H-2 u but not H-2 f, H-2 k, or H-2 s haplotypes.

**Applications Reported:** This M5/114.15.2 antibody has been reported for use in flow cytometric analysis.

**Applications Tested:** This M5/114.15.2 antibody has been tested y flow cytometric analysi sof mouse splenocytes. This can be used at less than or equal to 0.5 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10<sup>5</sup> to 10<sup>8</sup> cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

APC-eFluor 780 emits at 780 nm and is excited with the Red laser (633 nm). Please make sure that your instrument is capable of detecting this fluorocheme.


**Light sensitivity:** This tandem is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.

Fixation: Samples can be stored in IC Fixation Buffer (cat. 00-8222) (100  $\mu$ L cell sample + 100  $\mu$ L IC Fixation Buffer) or 1-step Fix /Lyse Solution (cat. 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency /compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Excitation: 633-647 nm; Emission: 780 nm; Laser: Red Laser.

Filtration: 0.2  $\mu$ m post-manufacturing filtered.

## Product Images For MHC Class II (I-A/I-E) Monoclonal Antibody (M5/114.15.2), APC-eFluor 780, eBioscience™

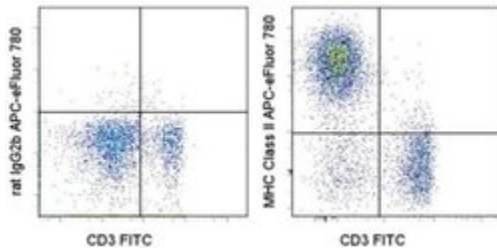

### MHC Class II (I-A/I-E) Antibody (47-5321-82) in Flow

Staining of C57BL/6 splenocytes with Anti-Mouse CD3e FITC (Product # 11-0031-82) and 0.25  $\mu$ g of Rat IgG2b K Isotype Control APC-eFluor® 780 (Product # 47-4031-82) (left) or 0.25  $\mu$ g of Anti-Mouse MHC Class II (I-A/I-E) APC-eFluor® 780 (right). Total viable cells were used for analysis.

[View more figures on thermofisher.com](#)

## 73 References

### Flow Cytometry (73)

International journal of biological sciences

#### MiR-103 protects from recurrent spontaneous abortion via inhibiting STAT1 mediated M1 macrophage polarization.

"Published figure using MHC Class II (I-A/I-E) monoclonal antibody (Product # 47-5321-82) in Flow Cytometry"

Authors: Zhu X,Liu H,Zhang Z,Wei R,Zhou X,Wang Z,Zhao L,Guo Q,Zhang Y,Chu C,Wang L,Li X

**Species**  
Not Applicable

**Dilution**  
Not Cited

**Year**  
2021

Theranostics

#### Pharmacological depletion of microglia and perivascular macrophages prevents Vascular Cognitive Impairment in Ang II-induced hypertension.

"47-5321 was used in Flow cytometry/Cell sorting to support the critical role of brain immune cells in the pathogenesis of hypertension-related cognitive impairment."

Authors: Kerkhofs D,van Hagen BT,Milanova IV,Schell KJ,van Essen H,Wijnands E,Goossens P,Blankesteyn WM, Unger T,Prickaerts J,Biessen EA,van Oostenbrugge RJ,Foulquier S

**Species**  
Mouse

**Dilution**  
Not Cited

**Year**  
2021

[View more Flow references on thermofisher.com](#)

## More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.