

CD11c Monoclonal Antibody (N418), Brilliant Ultra Violet 737, eBioscience™

| Product Details | |
|-----------------------------|---|
| Size | 100 µg |
| Species Reactivity | Mouse |
| Host/Isotype | Armenian hamster / IgG |
| Recommended Isotype Control | Armenian Hamster IgG Isotype Control (eBio299Arm), Brilliant Ultra Violet 737, eBioscience™ |
| Class | Monoclonal |
| Type | Antibody |
| Clone | N418 |
| Conjugate | Brilliant Ultra Violet™ 737 |
| Form | Liquid |
| Concentration | 0.2 mg/mL |
| Purification | Affinity chromatography |
| Storage buffer | PBS, pH 7.2, with BSA |
| Contains | 0.09% sodium azide |
| Storage conditions | 4° C, store in dark, DO NOT FREEZE! |

| Applications | Tested Dilution | Publications |
|-----------------------|-----------------|--------------|
| Flow Cytometry (Flow) | 0.5 µg/test | - |

Product Specific Information

Description: The N418 monoclonal antibody reacts with mouse CD11c, the integrin alphaX. CD11c non-covalently associates with beta2 integrin to form the CD11c/CD18 heterodimer. CD11c is expressed by dendritic cells, a subset of Intestinal Intraepithelial Lymphocytes (IEL) and some activated T cells. CD11c/CD18 binds to CD54, iC3b and fibrinogen and plays a role in leukocyte adhesive interactions. N418 binds to CD11c on splenic dendritic cells in the T-dependent areas of mouse spleen and precipitates a 150, 90 kDa heterodimer.

Applications Reported: This N418 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This N418 antibody has been tested by flow cytometric analysis of mouse splenocytes. This may be used at less than or equal to 0.5 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Brilliant Ultra Violet™ 737 (BUV737) is a tandem dye that emits at 732 nm and is intended for use on cytometers equipped with an ultraviolet (355 nm) laser. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright, Brilliant Violet™, Brilliant Ultra Violet™, or other polymer dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401) or Brilliant Stain Buffer (Product # 00-4409-75) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer or Brilliant Stain Buffer for more information.

Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.


Fixation: Samples can be stored in IC Fixation Buffer (Product # 00-8222) (100 μ L of cell sample + 100 μ L of IC Fixation Buffer) or 1-step Fix/Lyse Solution (Product # 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Our internal testing suggests that Brilliant Ultra Violet™ 737 (BUV737) is compatible with short-term methanol-based fixation, but should not be stored in buffers containing methanol for longer than one hour.

Excitation: 355 nm; Emission: 732 nm; Laser: Ultraviolet Laser.

BRILLIANT ULTRA VIOLET™ is a trademark or registered trademark of Becton, Dickinson and Company or its affiliates, and is used under license. Powered by Sirigen™

Product Images For CD11c Monoclonal Antibody (N418), Brilliant Ultra Violet 737, eBioscience™

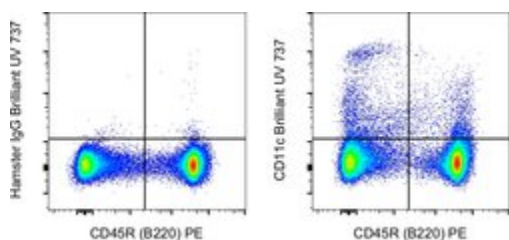

CD11c Antibody (367-0114-82) in Flow

C57BL/6 mouse splenocytes were stained with CD45R (B220) Monoclonal Antibody, PE (Product # 12-0452-82) and 0.25 μ g of Armenian Hamster IgG Isotype Control, Brilliant Ultra Violet 737 (BUV737) (Product # 367-4888-81) (left) or 0.25 μ g of CD11c Monoclonal Antibody, Brilliant Ultra Violet 737 (BUV737) (right). Total viable cells were used for analysis, as determined by 7-AAD (Product # 00-6993-50).

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.