

IL-13 Monoclonal Antibody (eBio13A), PE-Cyanine7, eBioscience™

Product Details	
Size	100 µg
Species Reactivity	Mouse
Published Species	Mouse, Human
Host/Isotype	Rat / IgG1, kappa
Recommended Isotype Control	Rat IgG1 kappa Isotype Control (eBRG1), PE-Cyanine7, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	eBio13A
Conjugate	PE-Cyanine7
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.1% gelatin
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2573530

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	0.5 µg/test	23 Publications

Product Specific Information

Description: The eBio13A antibody reacts with mouse IL-13. IL-13 is a cytokine produced mainly by Th2 cells, but also by antigen-primed CD8 T cells. IL-13 has a strong involvement in allergic inflammation and parasitic clearing and in cancer models has been shown to have either inhibitory or stimulatory activity depending on the tumor. In humans, IL-13 is found to play a role in isotype switching in B cells. IL-13 is implicating in down modulating macrophage activity, through the reduction of pro-inflammatory cytokines (IL-1, IL-6, IL-8, IL-10, IL-12)

Applications Reported: This eBio13A antibody has been reported for use in intracellular staining followed by flow cytometric analysis.

Applications Tested: This eBio13A antibody has been tested by intracellular staining of stimulated mouse splenocytes. This can be used at less than or equal to 0.5 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.

Fixation: Samples can be stored in IC Fixation Buffer (cat. 00-8222) (100 µL cell sample + 100 µL IC Fixation Buffer) or 1-step Fix /Lyse Solution (cat. 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency /compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific

performance should be determined empirically.


Excitation: 488-561 nm; Emission: 775 nm; Laser: Blue Laser, Green Laser, Yellow-Green Laser.

Filtration: 0.2 µm post-manufacturing filtered.

Advanced Verification Data

IL-13 Antibody (25-7133-82)


Figure 6 OX40L on ILC2s Orchestrates Adaptive Type 2 Immunity to Allergens (A-D) Mice of the specified genotypes were treated with PBS or papain (Pap) (i.n., day 0 and 1), followed by quantification of lung Th2, GATA3 + Treg, and GATA3 - Treg cells on day 5 (A-C). Whole lung cell suspensions were re-stimulated with PMA and ionomycin, followed by quantification of IL-13 + Th2 cells by intracellular staining (D). (E-G) Mice were treated with papain on days 0, 1, 14, and 21, followed by quantification on day 24 of lung eosinophils (E), and detection (F) and quantification (G) of RELMalpha + M2 alveolar (CD45 + SiglecF + CD11c + CD11b - F4/80 +) macrophages (MPhi). (H and I) Mice were treated with *A. alternata* (*A.alt*) (i.n., days 0 and 1), followed by quantification on day 9 of lung Th2, GATA3 + Treg, and GATA3 - Treg cells (H) and serum IgE concentration (I). Bar graphs indicate mean (+SEM). (A)-(D), ANOVA, two repeat experiments; (E), ANOVA, three repeat experiments; (F) and (G), ANOVA, two repeat experiments, representative gate shown in (F); (H) and (I), ANOVA, two repeat experiments. ** p <= 0.01, *** p <= 0.001, **** p <= 0.0001. See also Figure S7 . Cell treatment validation info.

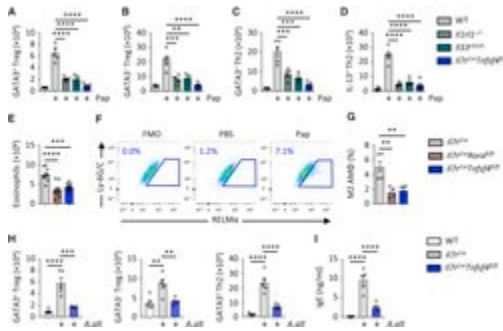

Product Images For IL-13 Monoclonal Antibody (eBio13A), PE-Cyanine7, eBioscience™

IL-13 Antibody (25-7133-82) in Flow

CD4+ BALB/c splenocytes were stimulated for 10-days under Th2 polarizing conditions. Cells were then restimulated with Cell Stimulation Cocktail (plus protein transport inhibitors) (Product # 00-4975-03) for 5 hours. Following restimulation, cells were fixed and permeabilized using the Intracellular Fixation & Permeabilization Buffer Set (Product # 88-8824-00) and then stained with Anti-Mouse CD4 eFluor® 450 (Product # 48-0042-82) and 0.25 µg of Rat IgG1 K Isotype Control PE-Cyanine7 (Product # 25-4301-82) (left) or 0.25 µg of Anti-Mouse IL-13 PE-Cyanine7 (right). Viable cells, as determined by Fixable Viability Dye eFluor® 660 (Product # 65-0864-14), were used for analysis.

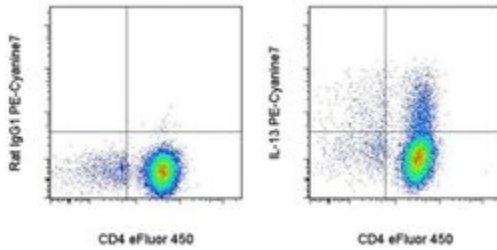

[View more figures on thermofisher.com](https://www.thermofisher.com)

Flow Cytometry (23)

British journal of pharmacology

Major histocompatibility complexes are up-regulated in glomerular endothelial cells via activation of c-Jun N-terminal kinase in 5/6 nephrectomy mice.

"25-7133 was used in Flow cytometry/Cell sorting to explore the mechanism underlying the up-regulation of major histocompatibility complex (MHC) proteins in glomerular endothelial cells in 5/6 nephrectomy mice."

Authors: Zhu D,Tang Q,Yu B,Meng M,Liu W,Li J,Zhu T,Vanhoutte PM,Leung SWS,Zhang Y,Shi Y

Species
Mouse

Dilution
1:100

Year
2020

Cell reports

T Cell-Intrinsic IRF5 Regulates T Cell Signaling, Migration, and Differentiation and Promotes Intestinal Inflammation.

"25-7133 was used in Flow cytometry/Cell sorting to identify a previously undefined key role for T cell-intrinsic IRF5. In mice, IRF5 in CD4+ T cells promotes Th1- and Th17-associated cytokines and decreases Th2-associated cytokines."

Authors: Yan J,Pandey SP,Barnes BJ,Turner JR,Abraham C

Species
Mouse

Dilution
Not Cited

Year
2020

[View more Flow references on thermofisher.com](#)

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.