

# CD4 Monoclonal Antibody (RPA-T4), eFluor™ 450, eBioscience™

## Product Details

Size	100 Tests
Species Reactivity	Human
Published Species	Human
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), eFluor™ 450, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	RPA-T4
Conjugate	eFluor™ 450
Form	Liquid
Concentration	5 µL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.1% gelatin, 0.2% BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_1272057

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	-	1 Publication
Immunocytochemistry (ICC/IF)	-	1 Publication
Flow Cytometry (Flow)	5 µL (0.25 µg)/test	37 Publications

## Product Specific Information

**Description:** The RPA-T4 monoclonal antibody reacts with human CD4, a 59 kDa cell surface receptor expressed by a majority of thymocytes, subpopulation of mature T cells (T-helper cells) and in low levels on monocytes. CD4 is a receptor for the human immunodeficiency virus (HIV). RPA-T4 blocks HIV binding and mixed lymphocyte reaction. The RPA-T4 antibody recognizes a different epitope than the OKT4 monoclonal antibody, and these antibodies do not cross-block binding to each other's respective epitopes.

**Applications Reported:** This RPA-T4 antibody has been reported for use in flow cytometric analysis.


**Applications Tested:** This RPA-T4 antibody has been pre-titrated and tested by flow cytometric analysis of human peripheral mononuclear cells. This can be used at 5 µL (0.25 µg) per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10<sup>5</sup> to 10<sup>8</sup> cells/test.

eFluor® 450 is an alternative to Pacific Blue®. eFluor® 450 emits at 445 nm and is excited with the Violet laser (405 nm). Please make sure that your instrument is capable of detecting this fluorochoime.

**Excitation:** 405 nm; **Emission:** 445 nm; **Laser:** Violet Laser.

Filtration: 0.2 µm post-manufacturing filtered.

## Product Images For CD4 Monoclonal Antibody (RPA-T4), eFluor™ 450, eBioscience™

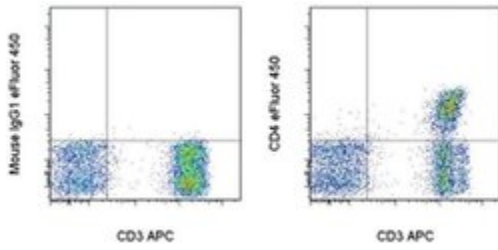

### CD4 Antibody (48-0049-42) in Flow

Staining of normal human peripheral blood cells with Anti-Human CD3 APC (Product # 17-0036-42) and Mouse IgG1 K Isotype Control eFluor® 450 (Product # 48-4714-82) (left) or Anti-Human CD4 eFluor® 450 (right). Cells in the lymphocyte gate were used for analysis.

[View more figures on thermofisher.com](http://thermofisher.com)

## Immunohistochemistry (1)

JCI insight

### CXCL13-producing TFH cells link immune suppression and adaptive memory in human breast cancer.

"Published figure using CD4 monoclonal antibody (Product # 48-0049-42) in Flow Cytometry"

Authors: Gu-Trantien C, Migliori E, Buisseret L, de Wind A, Brohée S, Garaud S, Noël G, Dang Chi VL, Lodewyckx JN, Naveaux C, Duvallier H, Goriely S, Larsimont D, Willard-Gallo K

**Species**  
Not Applicable

**Dilution**  
Not Cited

**Year**  
2017

## Immunocytochemistry (1)

Proceedings of the National Academy of Sciences of the United States of America

### Remarkably low affinity of CD4/peptide-major histocompatibility complex class II protein interactions.

"Published figure using CD4 monoclonal antibody (Product # 48-0049-42) in Immunofluorescence"

Authors: Jönsson P, Southcombe JH, Santos AM, Huo J, Fernandes RA, McColl J, Lever M, Evans EJ, Hudson A, Chang VT, Hanke T, Godkin A, Dunne PD, Horrocks MH, Palayret M, Screation GR, Petersen J, Rossjohn J, Fugger L, Dushek O, Xu XN, Davis SJ, Klenerman D

**Species**  
Not Applicable

**Dilution**  
Not Cited

**Year**  
2016

## Flow Cytometry (37)

Cancers

### Murlentamab, a Low Fucosylated Anti-Müllerian Hormone Type II Receptor (AMHRII) Antibody, Exhibits Anti-Tumor Activity through Tumor-Associated Macrophage Reprogramming and T Cell Activation.

"Published figure using CD4 monoclonal antibody (Product # 48-0049-42) in Flow Cytometry"

Authors: Prat M, Salon M, Allain T, Dubreuil O, Noël G, Preisser L, Jean B, Cassard L, Lemée F, Tabah-Fish I, Pipy B, Jeannin P, Prost JF, Barret JM, Coste A

**Species**  
Not Applicable

**Dilution**  
Not Cited

**Year**  
2021

Journal of Cancer

### HDAC Inhibitor, CG-745, Enhances the Anti-Cancer Effect of Anti-PD-1 Immune Checkpoint Inhibitor by Modulation of the Immune Microenvironment.

"Published figure using CD4 monoclonal antibody (Product # 48-0049-42) in Flow Cytometry"

Authors: Kim YD, Park SM, Ha HC, Lee AR, Won H, Cha H, Cho S, Cho JM

**Species**  
Not Applicable

**Dilution**  
Not Cited

**Year**  
2020

[View more Flow references on thermofisher.com](#)

## More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.