


## Washer Disinfector (Medical Automatic Glassware Washer)



Standard cleaning rack



Toughened glass door

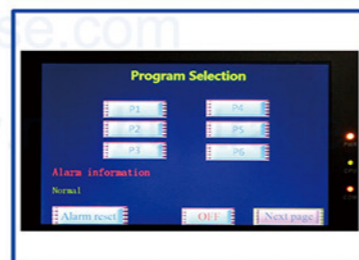

LCD display

### Introduction:

Automatic medical washing disinfector provides standard and automatic cleaning and disinfection for most of medical equipments, so as to eliminate the possibility of infection in cleaning and disinfecting.

### Features:

- \* Two toughened glass doors, prevent cross infection.
- \* Color LCD touch screen, advanced programmable control system with six standard programs.
- \* High efficient cleaning system, with rotary nozzles in top, bottom and among the shelves to make sure that the spray liquid can contact the surface of devices completely.
- \* High efficient drying system with HEPA filter, can clean and dry the inner of tubular things.
- \* Thermal disinfection, water heating system can make sure the water temperature keep above 93°C to reach the international disinfecting efficiency level.

### Standard Washing Procedure (can be self-designed by the user):

Pre-cleaning→Clean with enzymes →Rinse→Thermal disinfection→Glazing clean→Drying

### Technical Parameters:

Model	BK-LW220D	BK-LW320D	BK-LW420D
Capacity	220L	320L	420L
Water Consumption/Cycle	15 L	18L	25L
Inlet Water Pressure	0.03~0.8Mpa		
Noise	<60dB		
Material	Stainless steel		
Overall Power Consumption	18KW		
Power Consumption-Water Washing	2.8KW	25KW	
Power Consumption-Water Heating	15KW		
Power Consumption-Drying	3.6KW		
Power Supply	380V, 50Hz		
Standard Cleaning Racks	4 Layers	5 Layers	
	BK01-220*1	BK01-320*1	BK01-420*1
Standard Accessories	One bucket RBS A 150 leaner (1 L), One bucket RBS A 375 neutralizer (1 L)		
Optional Accessories	Trolley (BKO1)		
Washer Chamber Size(W*D*H)mm	600*629*658	600*629*820	688*700*703
External Size(W*D*H)mm	930*700*1840	930*700*1990	1100*900*2165
Package Size(W*D*H)mm	1100*1000*2200	1100*1000*2350	1270*1025*2340
Gross Weight(kg)	360	370	400