




ETL Certified Vertical Laminar Flow Cabinet



Email: support@biobase.com

Website: www.biobase.cc / www.meihuatrade.com

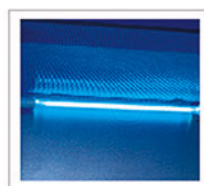

UV Lamp
Emission of 253,7nanometers
For most efficient decontamination



Waterproof Socket
America Standard



Footmaster Caster
Universal caster with brake
and leveling feet



Wind Speed Sensor
Direct measurement of wind speed

ADVANTAGE

1. ETL certified, Satisfy North American market.
2. With memory function in case of power-failure.
3. Audio and visual alarm(Abnormal airflow velocity).
4. Interlock function: UV Lamp and blower, LED Lamp.

TECHNICAL PARAMETERS

Model	BBS-V800	BBS-V1100	BBS-V1300	BBS-V1500	BBS-V1800
External Size(W*D*H)	800*700* (2090~2340)mm	1100*700* (2090~2340)mm	1300*700* (2090~2340)mm	1500*700* (2090~2340)mm	1800*700* (2090~2340)mm
Internal Size(W*D*H)	700*575*625mm	1000*575*625mm	1200*575*625mm	1400*575*625mm	1700*575*625mm
Work Surface Height	750~1010mm				
Max Opening	480mm				
Display	LCD Display				
Airflow Velocity	80~99fpm				
Material	Working area: 304stainless Steel Main Body: 1.2mm Cold-rolled steel with anti-bacteria powder coating. Front Window: Motorized,5mm toughened glass, anti-UV				
Pre-Filter	Polyester fiber				
HEPA Filter	One,99.995% efficiency at 0.3um				
Noise	≤67dB(A)				
Illumination	≥300lx				
Waterproof Sockets	One socket (Double type), America standard, Max.Power: 500W				
Power Supply	110V±10% 60HZ				
Consumption	750W	950W	950W	950W	1150W
UV Lamp	18W*1	30W*1	30W*1	40W*1	40W*1
LED Lamp	10W*1	14W*1	14W*1	18W*1	22W*1
Standard Accessory	Manually height adjustable base stand,LED lamp, UV lamp, wind speed sensor, waterproof socket				
Package Size (W*D*H)	940*970*1560mm	1240*970*1560mm	1440*970*1560mm	1640*970*1560mm	1940*970*1560mm
Gross Weight	190kg	230kg	265kg	280kg	320kg