
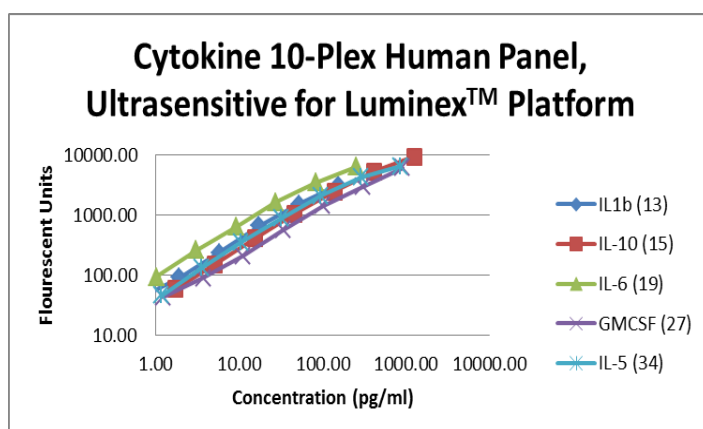


Cytokine 10-Plex Human Panel, Ultrasensitive Technical Data Sheet

Performance Characteristics

Typical Standard Curves (Note: Data presented here were collected as a 10-plex, but are depicted in two 5-plex groupings for clarity.)

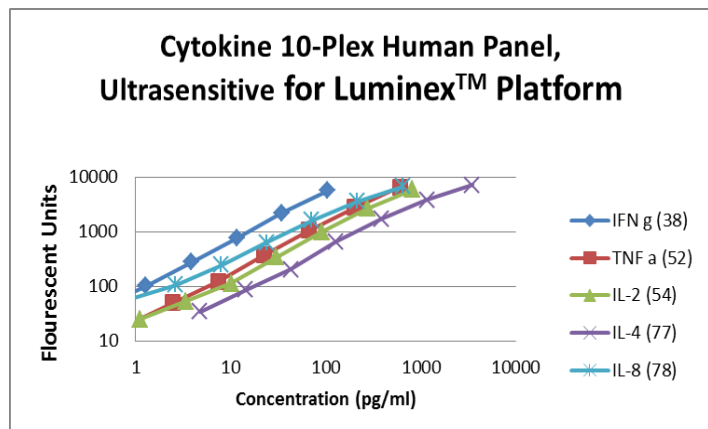

	IL-1β	IL-10	IL-6	GM CSF	IL-5	IFN γ	TNF α	IL-2	IL-4	IL-8
Magnetic Bead Region	13	15	19	27	34	38	52	54	77	78
Sensitivity	< 0.05	< 0.05	< 0.05	< 0.1	< 0.05	< 0.05	< 0.1	< 0.5	< 0.5	< 0.5
Recovery (Serum)	-	✓	✓	✓	✓	✓	-	✓	-	✓
Recovery (Heparin Plasma)	-	✓	-	✓	-	✓	-	-	-	✓
Recovery (EDTA Plasma)	✓	✓	✓	✓	✓	✓	✓	✓	-	✓
Recovery (Citrate Plasma)	-	✓	✓	✓	✓	✓	-	✓	-	✓
NIBSC calibration	1ng=25IU	1ng=13IU	1ng=30IU	1ng=6IU	1ng=0.2IU	1ng=0.7IU	1ng=0.2IU	1ng=7IU	1ng=65IU	1ng=40IU
NIBSC Preparation Number	86/680	93/722	89/548	88/646	90/586	82/587	88/786	86/564	88/656	89/520
Inter assay variation	8.8%	9.8%	7%	9.1%	7.5%	9%	8.3%	9.5%	8.7%	9.8%

Notes: 70-130% recovery (✓), 50-69% recovery (-), 131-150% recovery (+), <50% or >150% recovery (–not recommended)

Limited product warranty

Life Technologies Corporation and/or its affiliate(s) warrant their products as set forth in the Life Technologies' General Terms and Conditions of Sale found on Life Technologies' website at www.thermofisher.com/us/en/home/global/terms-and-conditions.html. If you have any questions, please contact Life Technologies at www.thermofisher.com/support.

Product label explanation of symbols and warnings

REF	Catalog Number	LOT	Batch code	Temperature limitation	Use by	Manufacturer	Consult instructions for use	Caution, consult accompanying documents
-----	----------------	-----	------------	------------------------	--------	--------------	------------------------------	---

Manufacturer's address: Bender MedSystems GmbH | Campus Vienna Biocenter 2 | 1030 Vienna, Austria

The information in this guide is subject to change without notice.

DISCLAIMER: TO THE EXTENT ALLOWED BY LAW, LIFE TECHNOLOGIES AND/OR ITS AFFILIATE(S) WILL NOT BE LIABLE FOR SPECIAL, INCIDENTAL, INDIRECT, PUNITIVE, MULTIPLE, OR CONSEQUENTIAL DAMAGES IN CONNECTION WITH OR ARISING FROM THIS DOCUMENT, INCLUDING YOUR USE OF IT.

Important Licensing Information: These products may be covered by one or more Limited Use Label Licenses. By use of these products, you accept the terms and conditions of all applicable Limited Use Label Licenses.

©2017 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific and its subsidiaries unless otherwise specified.

For Research Use Only. Not for use in diagnostic procedures.

MAN0017261 1.0 [00] Effective Date: 23 July 2017

Manufacturing site: Bender MedSystems GmbH | Campus Vienna Biocenter 2 | 1030 Vienna, Austria