
### invitrogen by Thermo Fisher Scientific

# eBioscience<sup>™</sup> Cell Stimulation Cocktail (plus protein transport inhibitors) (500X)

Catalog Number: 00-4975

For Research Use Only. Not for use in diagnostic procedures.

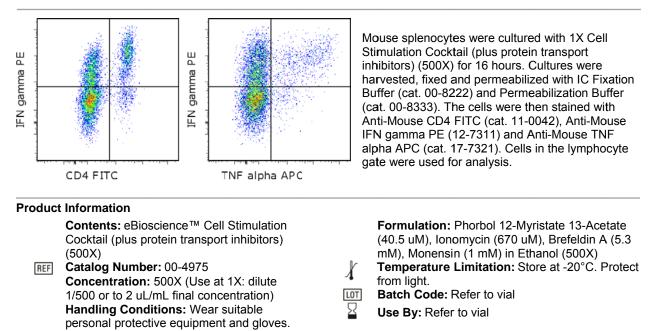

#### Description

The Cell Stimulation Cocktail (plus protein transport inhibitors) (500X) is a cocktail of phorbol 12-myristate 13-acetate (PMA), ionomycin, brefeldin A and monensin. Treatment with PMA and ionomycin is sufficient to induce activation of many cell types to produce cytokines. Brefeldin A and Monensin result in the accumulation of secreted proteins in the endoplasmic reticulum and Golgi apparatus. Thus, the Cell Stimulation Cocktail (plus protein transport inhibitors) (500X) can be used for induction and subsequent intracellular detection of cytokines and other secreted proteins in both *in vitro*-cultured and *ex vivo* cells.

For an unstimulated and vehicle control, see the Protein Transport Inhibitor Cocktail (500X) (cat. 00-4980).

#### **Applications Reported**

The Cell Stimulation Cocktail (plus protein transport inhibitors) (500X) has been reported for use in flow cytometric analysis.

#### **Applications Tested**

The Cell stimulation Cocktail (plus protein transport inhibitors) (500X) has been tested by treatment of mouse splenocytes for 5 to 18 hours followed by intracellular staining and flow cytometric analysis of IFN gamma. This cocktail is supplied as a 500X concentrate and can be directly added to culture medium at 1X (2  $\mu$ L/mL).

This product is formulated such that it will not freeze at the recommended storage temperatures. It is recommended that this product be kept tightly sealed and cold at all times to avoid excess evaporation.

#### **Related Products**

00-4970 eBioscience<sup>™</sup> Cell Stimulation Cocktail (500X) 00-4976 eBioscience<sup>™</sup> Lipopolysaccharide (LPS) Solution (500X) 00-4977 eBioscience<sup>™</sup> Phytohemagglutinin-L (PHA-L) Solution (500X)

Not for further distribution without written consent.

Copyright © 2016 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are property of Thermo Fisher Scientific and its subsidiaries unless otherwise specified.

Tel: 888.999.1371 or 858.642.2058 • Fax: 858.642.2046 • thermofisher.com/ebioscience • info@ebioscience.com



## eBioscience<sup>™</sup> Cell Stimulation Cocktail (plus protein transport inhibitors) (500X)

Catalog Number: 00-4975

For Research Use Only. Not for use in diagnostic procedures.

00-4978 eBioscience<sup>™</sup> Concanavalin A (Con A) Solution (500X) 00-4980 eBioscience<sup>™</sup> Protein Transport Inhibitor Cocktail (500X) 00-8222 eBioscience<sup>™</sup> IC Fixation Buffer 00-8333 eBioscience<sup>™</sup> Permeabilization Buffer (10X)

Not for further distribution without written consent. Copyright © 2016 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are property of Thermo Fisher Scientific and its subsidiaries unless otherwise specified. Tel: 888.999.1371 or 858.642.2058 • Fax: 858.642.2046 • thermofisher.com/ebioscience • info@ebioscience.com