

TLR4/MD-2 Complex Monoclonal Antibody (MTS510), Functional Grade, eBioscience™

Product Details	
Size	100 µg
Species Reactivity	Mouse
Published Species	Human, Mouse
Host/Isotype	Rat / IgG2a, kappa
Recommended Isotype Control	Rat IgG2a kappa Isotype Control (eBR2a), Functional Grade, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	MTS510
Conjugate	Functional Grade
Form	Liquid
Concentration	1 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2
Contains	no preservative
Storage conditions	4° C
RRID	AB_469285

Applications	Tested Dilution	Publications
Western Blot (WB)	-	2 Publications
Immunohistochemistry (IHC)	-	3 Publications
Immunohistochemistry (Frozen) (IHC (F))	-	1 Publication
Flow Cytometry (Flow)	1 µg/test	6 Publications
ELISA (ELISA)	-	1 Publication
Neutralization (Neu)	Assay-Dependent	5 Publications
Functional Assay (FN)	Assay-Dependent	-

Product Specific Information

Description: The MTS510 monoclonal antibody reacts with the mouse Toll-like receptor 4 (TLR4)/MD-2 complex. At least ten members of the Toll family have been identified. This family of type I transmembrane proteins is characterized by an extracellular domain with leucine-rich repeats and a cytoplasmic domain with homology to the type I IL-1 receptor. Two of these receptors, TLR2 and TLR4, are pattern recognition receptors and signaling molecules in response to bacterial lipoproteins and have been implicated in innate immunity and inflammation. TLR4 physically associates with MD-2, and together with CD14, this complex is responsible for LPS recognition and signaling. In the mouse, TLR4 is expressed by thioglycolate-elicited peritoneal macrophages. Incubation of peritoneal macrophages with LPS results in down regulation of surface TLR4/MD-2. The TLR4 gene is defective in C3H/HeJ and C57BL/10ScCr mice, both of which have been well characterized as hyporesponders to LPS.

The MTS510 monoclonal antibody co-immunoprecipitates MD-2 (~30 kDa) and TLR4 (~100 kDa), and preferentially reacts with TLR4 that is associated with MD-2. In comparison, binding of the UT41 monoclonal antibody occurs with and without formation of the TLR4/MD-2 complex. Please contact eBioscience Technical Support for further information.

Applications Reported: MTS510 has been reported for use in flow cytometric analysis and inhibition of LPS-induced cytokine production.

Applications Tested: The MTS510 antibody has been tested by flow cytometric analysis of mouse peritoneal exudates cells. This can be used at less than or equal to 1 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Storage and handling: Use in a sterile environment.


Filtration: 0.2 µm post-manufacturing filtered.

Purity: Greater than 90%, as determined by SDS-PAGE.

Endotoxin Level: Less than 0.001 ng/µg antibody, as determined by LAL assay.

Aggregation: Less than 10%, as determined by HPLC.

Product Images For TLR4/MD-2 Complex Monoclonal Antibody (MTS510), Functional Grade, eBioscience™

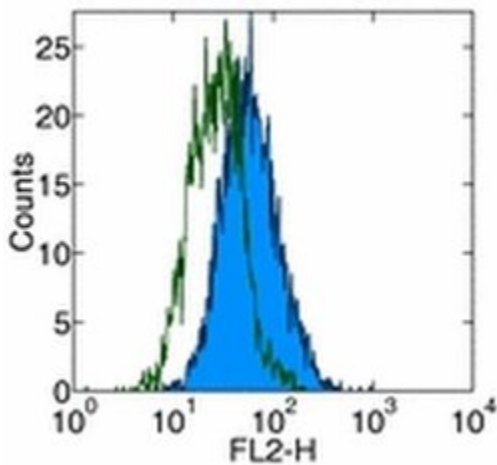

TLR4/MD-2 Complex Antibody (16-9924-82) in Flow

Staining of thioglycolate induced PECs with Anti-Mouse TLR4/MD-2 Complex PE. Appropriate isotype controls were used (open histogram). Cells in the large scatter population were used for analysis.

[View more figures on thermofisher.com](https://www.thermofisher.com)

18 References

Western Blot (2)

Physiology & behavior

Chronic exposure to low dose bacterial lipopolysaccharide inhibits leptin signaling in vagal afferent neurons.

"Published figure using TLR4/MD-2 Complex monoclonal antibody (Product # 16-9924-82) in Western Blot"

Authors: de La Serre CB,de Lartigue G,Raybould HE

Species
Not Applicable

Dilution
Not Cited

Year
2015

PloS one

Toll-like receptor 4 mediates inflammatory cytokine secretion in smooth muscle cells induced by oxidized low-density lipoprotein.

"Published figure using TLR4/MD-2 Complex monoclonal antibody (Product # 16-9924-82) in Western Blot"

Authors: Yang K,Zhang XJ,Cao LJ,Liu XH,Liu ZH,Wang XQ,Chen QJ,Lu L,Shen WF,Liu Y

Species
Mouse

Dilution
Not Cited

Year
2015

Immunohistochemistry (3)

EMBO molecular medicine

Identification of a novel mechanism of blood-brain communication during peripheral inflammation via choroid plexus-derived extracellular vesicles.

"Published figure using TLR4/MD-2 Complex monoclonal antibody (Product # 16-9924-82) in Immunofluorescence"

Authors: Balusu S, Van Wonterghem E, De Rycke R, Raemdonck K, Stremersch S, Gevaert K, Brkic M, Demeestere D, Vanhooren V, Hendrix A, Libert C, Vandembroucke RE

Species
Not Applicable

Dilution
Not Cited

Year
2016

Physiology & behavior

Chronic exposure to low dose bacterial lipopolysaccharide inhibits leptin signaling in vagal afferent neurons.

"Published figure using TLR4/MD-2 Complex monoclonal antibody (Product # 16-9924-82) in Western Blot"

Authors: de La Serre CB,de Lartigue G,Raybould HE

Species
Not Applicable

Dilution
Not Cited

Year
2015

[View more IHC references on thermofisher.com](#)

More applications with references on thermofisher.com

IHC (F) (1)

Flow (6)

ELISA (1)

Neu (5)

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.