

EOMES Monoclonal Antibody (Dan11mag), eFluor 660, eBioscience™

Product Details	
Size	100 µg
Species Reactivity	Mouse
Published Species	Mouse, Human
Host/Isotype	Rat / IgG2a, kappa
Recommended Isotype Control	Rat IgG2a kappa Isotype Control (eBR2a), eFluor 660, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	Dan11mag
Conjugate	eFluor® 660
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.1% gelatin
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2574227

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	-	9 Publications
Immunocytochemistry (ICC/IF)	-	5 Publications
Flow Cytometry (Flow)	0.25 µg/test	10 Publications

Product Specific Information

Description: This Dan11mag antibody recognizes Eomesodermin (EOMES), also known as T-box brain 2 (TBR2). Eomes is a T-box transcription factor that is highly homologous to T-bet, which is essential during trophoblast development and gastrulation in most vertebrates. In the immune system, Eomes controls the differentiation of effector and memory CD8+ T cells, as well as natural killer (NK) cells. Expression of Eomes in these cells correlates with high expression of CD122, the common beta-chain of the IL-2R and IL-15R.

Applications Reported: This Dan11mag antibody has been reported for use in intracellular staining followed by flow cytometric analysis.

Applications Tested: This Dan11mag antibody has been tested by intracellular staining and flow cytometric analysis of mouse splenocytes using the Foxp3/Transcription Factor Buffer Set (cat. 00-5523) and protocol. Please see Best Protocols Section (Staining intracellular Antigens for Flow Cytometry) for staining protocol (refer to Protocol B: One-step protocol for intracellular (nuclear) proteins). This can be used at less than or equal to 0.25 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸


cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

eFluor® 660 is a replacement for Alexa Fluor® 647. eFluor® 660 emits at 659 nm and is excited with the red laser (633 nm). Please make sure that your instrument is capable of detecting this fluorochrome.

Excitation: 633-647 nm; Emission: 668 nm; Laser: Red Laser.

Filtration: 0.2 µm post-manufacturing filtered.

Product Images For EOMES Monoclonal Antibody (Dan11mag), eFluor 660, eBioscience™

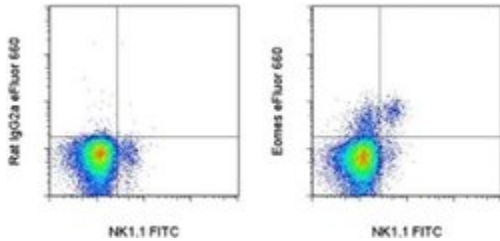

EOMES Antibody (50-4875-82) in Flow

Surface staining of C57Bl/6 splenocytes with Anti-Mouse NK1-1 FITC (Product # 11-5941-82), followed by fixation and permeabilization with the Foxp3/Transcription Factor Staining Buffer Set (Product # 00-5523-00) and intracellular staining with 0.125 µg of Rat IgG2a K Isotype Control eFluor® 660 (Product # 50-4321-82) (left) or 0.125 µg of Anti-Mouse EOMES eFluor® 660 (right). Cells in the lymphocyte gate were used for analysis.

[View more figures on thermofisher.com](https://www.thermofisher.com)

24 References

Immunohistochemistry (9)

Nature communications

LPA signaling acts as a cell-extrinsic mechanism to initiate cilia disassembly and promote neurogenesis.

"Published figure using EOMES monoclonal antibody (Product # 50-4875-82) in Immunohistochemistry"

Authors: Hu HB, Song ZQ, Song GP, Li S, Tu HQ, Wu M, Zhang YC, Yuan JF, Li TT, Li PY, Xu YL, Shen XL, Han QY, Li AL, Zhou T, Chun J, Zhang XM, Li HY

Species
Not Applicable

Dilution
Not Cited

Year
2021

Journal of neurodevelopmental disorders

Developmental exposure to diesel exhaust upregulates transcription factor expression, decreases hippocampal neurogenesis, and alters cortical lamina organization: relevance to neurodevelopmental disorders.

"Published figure using EOMES monoclonal antibody (Product # 50-4875-82) in Immunohistochemistry"

Authors: Cole TB, Chang YC, Dao K, Daza R, Hevner R, Costa LG

Species
Not Applicable

Dilution
Not Cited

Year
2020

[View more IHC references on thermofisher.com](#)

Immunocytochemistry (5)

Nature

Centrosome anchoring regulates progenitor properties and cortical formation.

"Published figure using EOMES monoclonal antibody (Product # 50-4875-82) in Immunocytochemistry"

Authors: Shao W, Yang J, He M, Yu XY, Lee CH, Yang Z, Joyner AL, Anderson KV, Zhang J, Tsou MB, Shi H, Shi SH

Species
Not Applicable

Dilution
Not Cited

Year
2020

Molecular autism

Loss of the neurodevelopmental disease-associated gene miR-146a impairs neural progenitor differentiation and causes learning and memory deficits.

"Published figure using EOMES monoclonal antibody (Product # 50-4875-82) in Immunocytochemistry"

Authors: Fregeac J, Moriceau S, Poli A, Nguyen LS, Oury F, Colleaux L

Species
Not Applicable

Dilution
Not Cited

Year
2020

[View more ICC/IF references on thermofisher.com](#)

More applications with references on thermofisher.com

Flow (10)

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.