

CD11c Monoclonal Antibody (N418), Alexa Fluor 532, eBioscience™

Product Details	
Size	100 µg
Species Reactivity	Mouse
Published Species	Human, Mouse
Host/Isotype	Armenian hamster / IgG
Recommended Isotype Control	Armenian Hamster IgG Isotype Control (eBio299Arm), Alexa Fluor 532, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	N418
Conjugate	Alexa Fluor® 532
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.1% gelatin
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_11217475

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	-	12 Publications
Immunohistochemistry (Frozen) (IHC (F))	-	3 Publications
Immunocytochemistry (ICC/IF)	-	5 Publications
Flow Cytometry (Flow)	0.25 µg/test	71 Publications
Miscellaneous PubMed (Misc)	-	1 Publication

Product Specific Information

Description: The N418 monoclonal antibody reacts with mouse CD11c, the integrin alpha X. CD11c non-covalently associates with beta 2 integrin to form the CD11c/CD18 heterodimer. CD11c is expressed by dendritic cells, a subset of Intestinal Intraepithelial Lymphocytes (IEL) and some activated T cells. CD11c/CD18 binds to CD54, iC3b and fibrinogen and plays a role in leukocyte adhesive interactions. N418 binds to CD11c on splenic dendritic cells in the T-dependent areas of mouse spleen and precipitates a 150, 90 kDa heterodimer.

Applications Reported: This N418 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This N418 antibody has been tested by flow cytometric analysis of mouse splenocytes. This can be used at less than or equal to 0.25 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.


Alexa Fluor® 532 is excited with the Green laser (532 nm) and emits at 561 nm. This cannot be used with the Yellow-Green laser (561 nm). We recommend using a 560/14 band pass filter. Please make sure that your instrument is capable of detecting this fluorochrome.

Excitation: 532 nm; Emission: 561 nm; Laser: Green Laser

Product Images For CD11c Monoclonal Antibody (N418), Alexa Fluor 532, eBioscience™

CD11c Antibody (58-0114-82) in Flow

Staining of BALB/c splenocytes with Anti-Human/Mouse CD45R (B220) APC (Product # 17-0452-82) and staining buffer (autofluorescence) (left) or 0.125 µg of Anti-Mouse CD11c Alexa Fluor® 532 (right). Total viable cells were used for analysis.

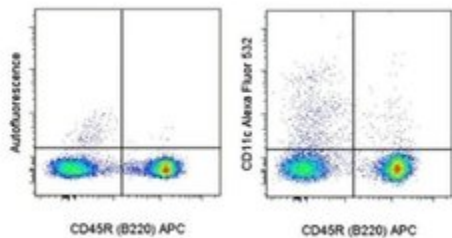

[View more figures on thermofisher.com](http://thermofisher.com)

Immunohistochemistry (12)

EMBO molecular medicine

Glufosinate constrains synchronous and metachronous metastasis by promoting anti-tumor macrophages.

"Published figure using CD11c monoclonal antibody (Product # 58-0114-82) in Flow Cytometry"

Authors: Menga A, Serra M, Todisco S, Riera-Domingo C, Ammarah U, Ehling M, Palmieri EM, Di Noia MA, Gissi R, Favia M, Pierri CL, Porporato PE, Castegna A, Mazzone M

Species
Not Applicable

Dilution
Not Cited

Year
2020

International journal of molecular medicine

Identification of different macrophage subpopulations with distinct activities in a mouse model of oxygen-induced retinopathy.

"Published figure using CD11c monoclonal antibody (Product # 58-0114-82) in Immunofluorescence"

Authors: Zhu Y, Zhang L, Lu Q, Gao Y, Cai Y, Sui A, Su T, Shen X, Xie B

Species
Not Applicable

Dilution
Not Cited

Year
2017

[View more IHC references on thermofisher.com](#)

Immunohistochemistry (Frozen) (3)

The Journal of experimental medicine

Essential roles of DC-derived IL-15 as a mediator of inflammatory responses in vivo.

"Published figure using CD11c monoclonal antibody (Product # 58-0114-82) in Immunohistochemistry"

Authors: Ohteki T, Tada H, Ishida K, Sato T, Maki C, Yamada T, Hamuro J, Koyasu S

Species
Not Applicable

Dilution
Not Cited

Year
2006

International immunology

The existence of CD11c+ sentinel and F4/80+ interstitial dendritic cells in dental pulp and their dynamics and functional properties.

Authors: Zhang J, Kawashima N, Suda H, Nakano Y, Takano Y, Azuma M

Species
Not Applicable

Dilution
Not Cited

Year
2006

[View more IHC \(F\) references on thermofisher.com](#)

More applications with references on thermofisher.com

ICC/IF (5)

Flow (71)

Misc (1)

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.