

Ki-67 Monoclonal Antibody (SolA15), PE-Cyanine7, eBioscience™

Product Details	
Size	100 µg
Species Reactivity	Dog, Cynomolgus monkey, Human, Mouse, Non-human primate, Rat
Published Species	Mouse, Human
Host/Isotype	Rat / IgG2a, kappa
Recommended Isotype Control	Rat IgG2a kappa Isotype Control (eBR2a), PE-Cyanine7, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	SolA15
Conjugate	PE-Cyanine7
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.1% gelatin
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_11220070

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	0.125 µg/test	19 Publications
Functional Assay (FN)	-	1 Publication

Product Specific Information

Description: The monoclonal antibody SolA15 recognizes mouse and rat Ki-67, a 300 kDa nuclear protein. Ki-67 is present during all active phases of the cell cycle (G1, S, G2, and mitosis), but is absent from resting cells (G0). Ki-67 is detected within the nucleus during interphase but redistributes to the chromosomes during mitosis. Ki-67 is used as a marker for determining the growth fraction of a given population of cells. In studies of tumor cells, the "Ki-67 labeling index" refers to the number of Ki-67 positive cells within the population and this is used to predict outcome of particular cancer types. Ki-67 has been shown to interact with the DNA-bound protein chromobox protein homolog 3 (CBX3) (heterochromatin).

The SolA15 antibody also recognizes human, non-human primate and canine Ki-67.

Applications Reported: This SolA15 antibody has been reported for use in intracellular staining followed by flow cytometric analysis.

Applications Tested: This SolA15 antibody has been tested by intracellular staining and flow cytometric analysis using the Foxp3 /Transcription Factor Staining Buffer Set (cat. 00-5523) and protocol. Please see Best Protocols Section (Staining intracellular Antigens for Flow Cytometry) for staining protocol (refer to Protocol B: One-step protocol for intracellular (nuclear) proteins). This can be used at less than or equal to 0.125 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.


Light sensitivity: This tandem dye is sensitive photo-induced oxidation. Please protect this vial and stained samples from light.

Fixation: Samples can be stored in IC Fixation Buffer (cat. 00-8222) (100 μ L cell sample + 100 μ L IC Fixation Buffer) or 1-step Fix/Lyse Solution (cat. 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Excitation: 488-561 nm; Emission: 775 nm; Laser: Blue Laser, Green Laser, Yellow-Green Laser.

Filtration: 0.2 μ m post-manufacturing filtered.

Product Images For Ki-67 Monoclonal Antibody (SoIA15), PE-Cyanine7, eBioscience™

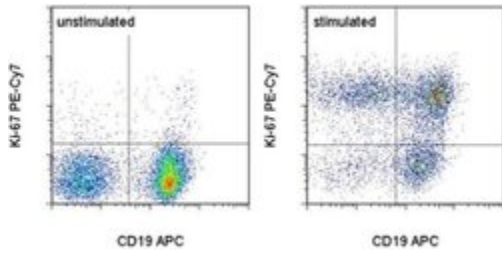

Ki-67 Antibody (25-5698-82) in Flow

C57Bl/6 splenocytes were unstimulated (left) or stimulated for 2 days with Anti-Mouse CD3 Functional Grade Purified (Product # 16-0031-82) (right). Cells were stained with Anti-Mouse CD19 APC (Product # 17-0193-82) followed by fixation and permeabilization using the Foxp3 Staining Buffer Set (Product # 00-5523-00) and subsequently stained with 0.06 μ g of Anti-Mouse/Rat Ki-67 PE-Cyanine7. Cells in the lymphocyte gate were used for analysis.

20 References

Flow Cytometry (19)

Nature communications

Thymic iNKT single cell analyses unmask the common developmental program of mouse innate T cells.

"25-5698 was used in Flow cytometry/Cell sorting to highlight the common requirements for the development of innate-like T cells with a focus on how Hivep3 impacts the maturation of these lymphocytes."

Authors: Harsha Krovi S,Zhang J,Michaels-Foster MJ,Brunetti T,Loh L,Scott-Browne J,Gapin L

Species
Mouse

Dilution
Not Cited

Year
2020

Nature communications

Luminal Galectin-9-Lamp2 interaction regulates lysosome and autophagy to prevent pathogenesis in the intestine and pancreas.

"25-5698 was used in Flow cytometry/Cell sorting to show in gut epithelial cells, galectin-9 is enriched in lysosomes and predominantly binds to lysosome-associated membrane protein 2 (Lamp2) in a Asn(N)-glycan dependent manner."

Authors: Sudhakar JN,Lu HH,Chiang HY,Suen CS,Hwang MJ,Wu SY,Shen CN,Chang YM,Li FA,Liu FT,Shui JW

Species
Mouse

Dilution
1:200

Year
2020

[View more Flow references on thermofisher.com](#)

Functional Assay (1)

PloS one

Myeloid cells expressing VEGF and arginase-1 following uptake of damaged retinal pigment epithelium suggests potential mechanism that drives the onset of choroidal angiogenesis in mice.

Authors: Liu J,Copland DA,Horie S,Wu WK,Chen M,Xu Y,Paul Morgan B,Mack M,Xu H,Nicholson LB,Dick AD

Species
Not Applicable

Dilution
Not Cited

Year
2015

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.