

CD107a (LAMP-1) Monoclonal Antibody (eBioH4A3), eFluor 450, eBioscience™

Product Details	
Size	100 Tests
Species Reactivity	Human
Published Species	Human
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), eFluor 450, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	eBioH4A3
Conjugate	eFluor® 450
Form	Liquid
Concentration	5 µL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.1% gelatin, 0.2% BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2574035

Applications	Tested Dilution	Publications
Immunocytochemistry (ICC/IF)	-	3 Publications
Flow Cytometry (Flow)	5 µL (0.25 µg)/test	5 Publications

Product Specific Information

Description: The eBioH4A3 monoclonal antibody reacts with human CD107a, also known as lysosomal-associated membrane protein-1 (LAMP-1). CD107a is a highly glycosylated protein of approximately 110kDa. It is predominantly expressed intracellularly in the lysosomal/endosomal membrane in nearly all cells. CD107a is transiently expressed on the cell surface of degranulating cytolytic T cells, and is also upregulated on the surface of activated platelets and some cancer cells.

Applications Reported: This eBioH4A3 antibody has been reported for use in intracellular staining followed by flow cytometric analysis.


Applications Tested: This eBioH4A3 antibody has been pre-titrated and tested by intracellular staining and flow cytometric analysis of Jurkat cells using the Intracellular Fixation & Permeabilization Buffer Set (cat. 88-8824) and protocol. Please refer to Best Protocols: Protocol A: Two step protocol for (cytoplasmic) intracellular proteins located under the Resources Tab online. This can be used at 5 µL (0.25 µg) per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test.

eFluor® 450 is an alternative to Pacific Blue®. eFluor® 450 emits at 445 nm and is excited with the Violet laser (405 nm). Please make sure that your instrument is capable of detecting this fluorochoime.

Excitation: 405 nm; Emission: 445 nm; Laser: Violet Laser.

Filtration: 0.2 µm post-manufacturing filtered.

Product Images For CD107a (LAMP-1) Monoclonal Antibody (eBioH4A3), eFluor 450, eBioscience™

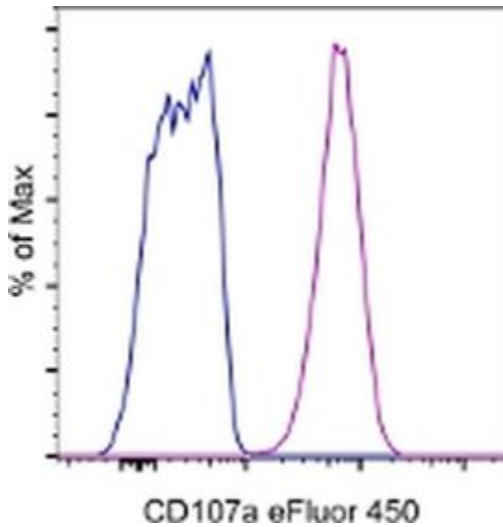

CD107a (LAMP-1) Antibody (48-1079-42) in Flow

Intracellular staining of Jurkat cells with Mouse IgG1 K Isotype Control eFluor® 450 (Product # 48-4714-82) (blue histogram) or Anti-Human CD107a (LAMP-1) eFluor® 450 (purple histogram) using the Intracellular Fixation & Permeabilization Buffer Set (Product # 88-8824-00) and protocol. Total viable cells, as determined by Fixable Viability Dye eFluor® 660 (Product # 65-0864-14), were used for analysis.

[View more figures on thermofisher.com](https://www.thermofisher.com)

8 References

Immunocytochemistry (3)

The Journal of clinical investigation

Targeting cancer metabolism by simultaneously disrupting parallel nutrient access pathways.

"Published figure using CD107a (LAMP-1) monoclonal antibody (Product # 48-1079-42) in Immunofluorescence"

Authors: Kim SM,Roy SG,Chen B,Nguyen TM,McMonigle RJ,McCracken AN,Zhang Y,Kofuji S,Hou J,Selwan E,Finicle BT,Nguyen TT,Ravi A,Ramirez MU,Wiher T,Guenther GG,Kono M,Sasaki AT,Weisman LS,Potma EO,Tromberg BJ, Edwards RA,Hanessian S,Edinger AL

Species
Human

Dilution
Not Cited

Year
2016

Cell death & disease

Fucosylation of LAMP-1 and LAMP-2 by FUT1 correlates with lysosomal positioning and autophagic flux of breast cancer cells.

"Published figure using CD107a (LAMP-1) monoclonal antibody (Product # 48-1079-42) in Immunofluorescence"

Authors: Tan KP,Ho MY,Cho HC,Yu J,Hung JT,Yu AL

Species
Human

Dilution
Not Cited

Year
2016

[View more ICC/IF references on thermofisher.com](#)

Flow Cytometry (5)

Oncoimmunology

T cells targeting NY-ESO-1 demonstrate efficacy against disseminated neuroblastoma.

"Published figure using CD107a (LAMP-1) monoclonal antibody (Product # 48-1079-42) in Flow Cytometry"

Authors: Singh N,Kulikovskaya I,Barrett DM,Binder-Scholl G,Jakobsen B,Martinez D,Pawel B,June CH,Kalos MD, Grupp SA

Species
Not Applicable

Dilution
Not Cited

Year
2021

PLoS pathogens

HLA-B*27:05 alters immunodominance hierarchy of universal influenza-specific CD8+ T cells.

"Published figure using CD107a (LAMP-1) monoclonal antibody (Product # 48-1079-42) in Flow Cytometry"

Authors: Sant S,Quiñones-Parra SM,Koutsakos M,Grant EJ,Loudovaris T,Mannerling SI,Crowe J,van de Sandt CE, Rimmelzwaan GF,Rossjohn J,Gras S,Loh L,Nguyen THO,Kedzierska K

Species
Not Applicable

Dilution
Not Cited

Year
2020

[View more Flow references on thermofisher.com](#)

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON-INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.