


c-Myc Monoclonal Antibody (9E10), Alexa Fluor 647

Product Details	
Size	50 µL
Species Reactivity	Human
Published Species	Human
Host/Isotype	Mouse / IgG1
Class	Monoclonal
Type	Antibody
Clone	9E10
Conjugate	Alexa Fluor® 647
Immunogen	Synthetic peptide A(408) E E Q K L I S E E D L L R K R R E Q L K H K L E Q L R N S C A(438) of human c-Myc
Form	Liquid
Concentration	1 mg/mL
Purification	Protein A
Storage buffer	PBS with 1mg/mL BSA
Contains	0.05% sodium azide
Storage conditions	4° C, do not freeze
RRID	AB_2609821

Applications	Tested Dilution	Publications
Western Blot (WB)	1:500-1:1,000	-
Immunocytochemistry (ICC/IF)	1:25-1:100	-
Flow Cytometry (Flow)	-	1 Publication

Product Images For c-Myc Monoclonal Antibody (9E10), Alexa Fluor 647

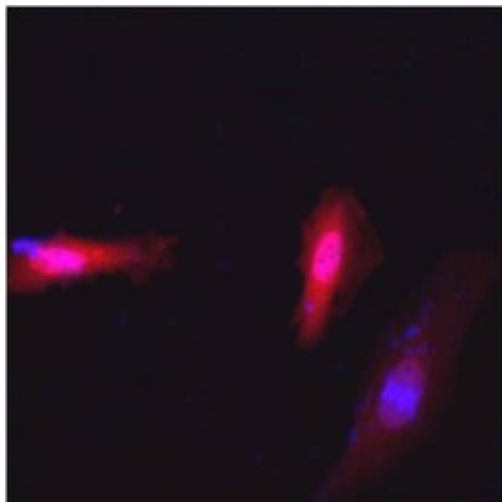


c-Myc Antibody (MA1-980-A647) in ICC/IF

Immunofluorescent analysis of HeLa cells transfected with a construct containing a MYC Epitope Tag. Cells were fixed with 4% paraformaldehyde, permeabilized with 0.1% Triton X-100 in PBS for 15 minutes at room temperature and blocked with 3% Blocker BSA (Product # 37525) for 15 minutes at room temperature. Cells were stained with AlexaFluor 647 conjugated MYC-tag monoclonal antibody (Product # MA1-980-A647) at a dilution of 1:25 for 1 hour at room temperature. Nuclei (blue) were stained with Hoechst 33342 dye (Product # 62249). Images were taken on a Thermo Scientific ToxInsight at 20X magnification.


c-Myc Antibody (MA1-980-A647) in WB

Western blot was performed using Anti-Myc Tag Polyclonal Antibody (Product # MA1-980-A647) by loading whole cell extracts of untransfected and transiently transfected HEK-293E lysates: untransfected, 40 μ g (Lane 1), empty vector control, 60 μ g (Lane 2), Myc-p65-V5, 60 μ g (Lane 3), Myc-p65-V5, 40 μ g (Lane 4), Myc-p65-V5, 20 μ g (Lane 5), Myc-p65-V5, 10 μ g (Lane 6), p65-HA, 60 μ g (Lane 7) and 25 ng of Positope (Lane 8) were electrophoresed using NuPAGE™ 4-12% Bis-Tris Protein Gel (Product # NP0321BOX). Resolved proteins were then transferred onto a nitrocellulose membrane (Product # IB23001) by iBlot® 2 Dry Blotting System (Product # IB21001). A ~65 kDa band corresponding to Myc-p65-V5 was observed in HEK293E transfected lysates on probing with the primary antibody (1:1000 dilution) and detected by chemiluminescence with Goat Mouse IgG (H+L) Superclonal™ Secondary Antibody, HRP (Product # A28177, 1:4000 dilution). Positope (Product # R90050) is a 53 kDa recombinant protein consisting of multiple epitope tags but lacks HA tag, which has been used as a positive control for c-Myc tag detection. No cross-reactivity was seen with HA-tagged p65 expressing lysate. Figure 1a is a non-overlay chemiluminescent image where construct Myc-p65-V5 has picked up 65 kDa band at different lysate concentrations and figure 1b is a fluorescent image of the same.

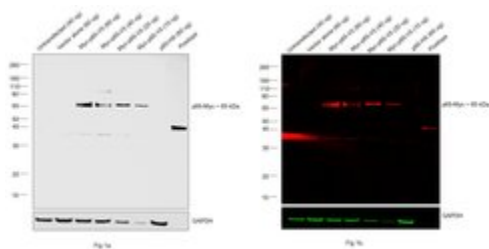

c-Myc Antibody (MA1-980-A647) in WB

Western blot analysis of Myc Epitope Tag was performed by loading various amounts of E. coli lysate containing a multi-epitope tagged protein per well onto a 4-20% Tris-HCl polyacrylamide gel. Proteins were transferred to a low fluorescence PVDF membrane and blocked with Sea Block blocking buffer for at least 1 hour. The membrane was probed with a AlexaFluor 647-conjugated Myc Epitope Tag monoclonal antibody (Product # MA1-980-A647) at a dilution of 1:500 for 1 hour at room temperature on a rocking platform and washed in TBS-0.1% Tween-20. Detection was performed using a fluorescence imaging system.

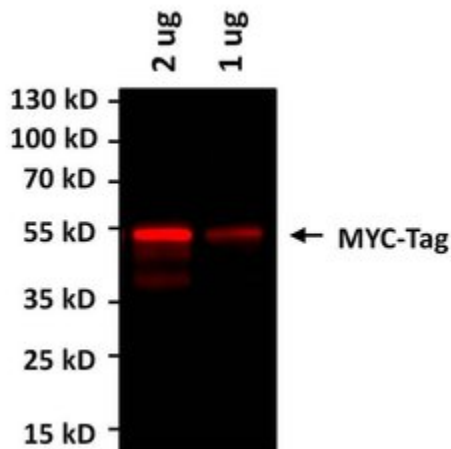

View more figures on thermofisher.com

Flow Cytometry (1)

Bioconjugate chemistry

Engineering Reversible Cell-Cell Interactions with Lipid Anchored Prosthetic Receptors.

"MA1-980-A647 was used in Flow cytometry/Cell sorting to develop a nongenetic method to rapidly, stably, and reversibly modify any cell membrane with a chemically self-assembled nanoring that can function as a prosthetic receptor."

Authors: Csizmar CM,Petersburg JR,Hendricks A,Stern LA,Hackel BJ,Wagner CR

Species
Human

Dilution
Not Cited

Year
2018

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.