



Stay in control with Thermal Cycler Fleet Control Software

Manage multiple users, methods, and instruments with a single intuitive interface


Applied Biosystems™ Thermal Cycler Fleet Control Software is as intuitive and easy to use as our thermal cyclers. With this tool, you will be able to view and control all Applied Biosystems™ thermal cyclers.

Secure—Thermal Cycler Fleet Control Software comes preinstalled on a server that sits securely behind your firewall

Flexible—create custom permissions and rules for users, methods, and instruments

Powerful—quickly view your thermal cycler inventory and availability right from the dashboard

Convenient—control your thermal cyclers, reports, methods, and firmware upgrades from your office, another lab computer, or any other networked device

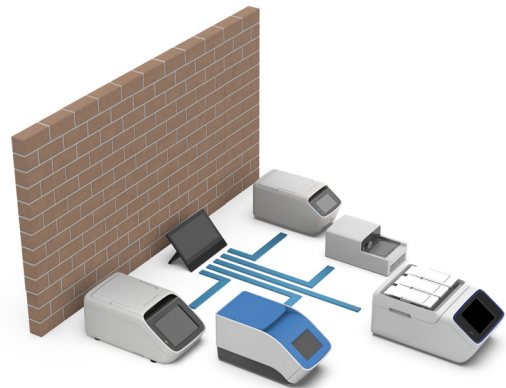


Enjoy peace of mind knowing that Thermal Cycler Fleet Control Software operates completely inside of your organization's firewall.



Superior service and support

Help reduce downtime with rapid exchange support

Thermal Cycler Fleet Control Software arrives preinstalled on a server and can be operated from the server or from any computer or device on the same network using the Google™ Chrome™ internet browser.

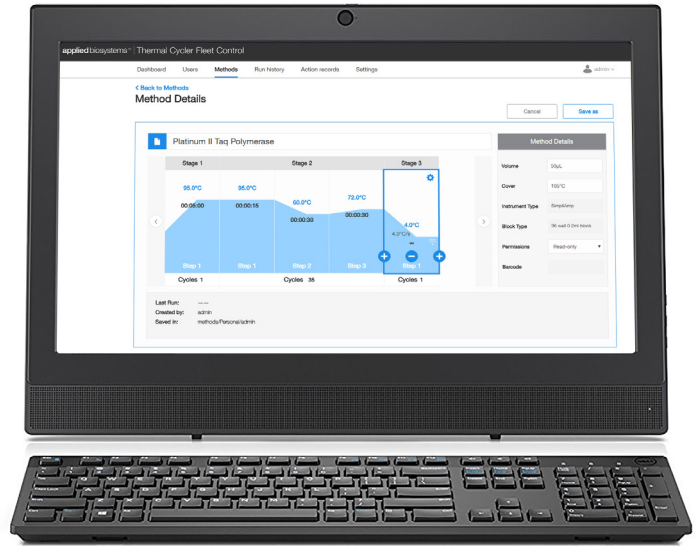

Specifications

Compatible instruments	ProFlex PCR System, Veriti Thermal Cycler, SimpliAmp Thermal Cycler, Automated Thermal Cycler, and MiniAmp Thermal Cyclers
Key capabilities	<ul style="list-style-type: none"> View all instruments together on a dashboard Automatically store and easily access run reports Manage multiple user accounts with password protection Track system changes with an audit log Identify methods or instruments via barcode data entry
Hardware	Dell™ Optiplex™ 3050 All-in-One Computer server with touchscreen, keyboard, and mouse
Hardware memory	128 GB solid-state hard drive
Hardware display interface	19.5-inch wide-viewing angle, HD, 1600 x 900 resolution
Hardware power	External power supply unit; 130 W adapter
Network connections	Ethernet LAN 10/100/1000 or Wi-Fi
Compatible browser	Chrome
Standard warranty	3 years

Ordering information

Product	Cat. No.
Thermal Cycler Fleet Control Software	A40070
ProFlex 96-Well PCR System*	4484075
SimpliAmp Thermal Cycler	A24811
Veriti 96-Well Thermal Cycler*	4375786
Automated Thermal Cycler, 96-well*	A31486
MiniAmp Plus Thermal Cycler*	A37835

* Thermal Cycler Fleet Control Software is compatible with all versions of Applied Biosystems thermal cyclers.

Find out more at thermofisher.com/fleetcontrol

ThermoFisher
SCIENTIFIC

For Research Use Only. Not for use in diagnostic procedures. © 2019 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific and its subsidiaries unless otherwise specified. Dell and Optiplex are trademarks of Dell Inc. Google and Chrome are trademarks of Google Inc. COL18674 0319