


Multi Sealer MS 500



Reliable sealing process

- RF Hand-held sealer
- Compact and handy
- Power cord or battery operation

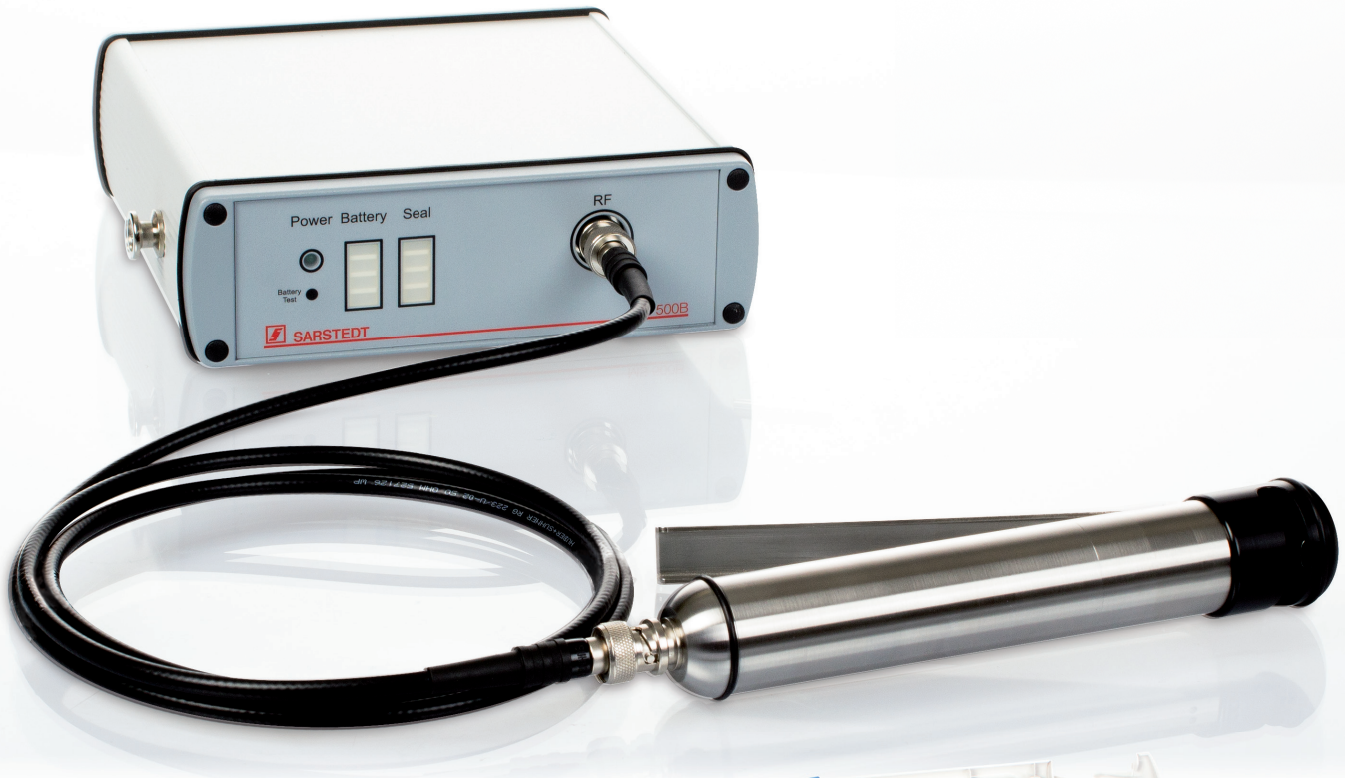

Multi Sealer MS 500

The new hand-held tube sealers, MS 500 M (mains) and MS 500 B (battery operated) are specially designed for use during blood donation or for the production of blood components.

The hand-held sealer features a specialised mechanical system that compensates for differences in user handling in order to ensure secure and consistent sealing.

Product features

- Simple segment separation
- Indicator light and audible signal after sealing
- Easy electrode replacement
- Easy to clean
- Replaceable battery (MS 500 B)
- Robust RF hand-held sealer
- Solid case

Technical data

	MS 500 M	MS 500 B
Power supply	115 - 230 V, 50/60 Hz	-
NiMH battery	-	24 V / 1700 mAh (replaceable)
Seals/charge	-	> 500
Full charge	-	ca. 3 hours
Charger	-	Supply voltage: 100 - 240 V, 50/60 Hz Output voltage: 41 V DC, 900 mA
Instrument dimensions (WxHxD)	17.0 x 5.5 x 23.5 cm	18.0 x 5.5 x 23.0 cm
Instrument weight	2.0 kg	2.0 kg
Hand sealer dimensions (Diameter x L)	3.1 x 18.5 cm	3.1 x 18.5 cm
Hand sealer weight	0.275 kg	0.275 kg
Case dimensions (WxHxD)	39.0 x 87.0 x 32.5 cm	39.0 x 87.0 x 32.5 cm
RF cable length	1.8 m	1.8 m

Ordering information

Order No.	Description
87.7708.850	MS 500 M (mains), includes hand-held sealer, RF cable, case
87.7708.860	MS 500 B (battery-operated), includes hand-held sealer, RF-cable, 1700 mAh battery, charger, strap, case

SARSTEDT AG & Co. KG
P.O. Box 12 20
D-51582 Nümbrecht
Phone: +49 2293 305 0
Fax: +49 2293 305 3992
export@sarstedt.com
www.sarstedt.com